



**Multi Connection System from V1.6.x**

## **Output Format for Host Connection**

Ref: RAA088BEN



## Output Format for Host Connection

---



CE

IVD



**HORIBA ABX SAS**  
Parc Euromédecine  
Rue du Caducée  
BP 7290  
34184 Montpellier Cedex 4  
FRANCE



# Contents

---

<b>1. Foreword</b> .....	<b>1</b>
1.1. Document Update.....	1
<b>2. HL7 Format</b> .....	<b>2</b>
2.1. Protocol Description.....	2
2.2. OUL and OML Messages Segments.....	6
2.3. Laboratory Testing Workflow Examples.....	28
2.4. Patient Detailed Structure.....	42
2.5. Minimal Lower Layer Protocol.....	46
2.6. References.....	47
<b>3. QC to LIS</b> .....	<b>48</b>
3.1. Overview.....	48
3.2. LIS to Yumizen P8000 ACK / NAK.....	48
3.3. MSH - Message Header Segment.....	49
3.4. OBR - Observation Request Segment.....	49
3.5. OBX - Observation Result Segment.....	50
3.6. ORU^R01.....	52
3.7. QC to LIS - Detailed Structure.....	54
<b>4. QC Export</b> .....	<b>56</b>
4.1. QC Results Availability.....	56



# 1. Foreword

## 1.1. Document Update

---

### 1.1.1. Revisions

Internal Reference	Software Version used for Documentation	Document Date Issued
RAA088AEN	1.6.x	03/2023
RAA088BEN	1.6.x	04/2024

This document applies to the latest software version listed and higher versions.

When a subsequent software version changes the information in this document, a new electronic edition is released and supplied by HORIBA Medical.

This document is only available online at [www.horiba-abx.com/documentation](http://www.horiba-abx.com/documentation).

### 1.1.2. What's New?

Update	Chapter
Parameter added: ESR	<ul style="list-style-type: none"><li>Parameters</li><li>OBR - Observation Request Segment</li></ul>
Control added: ESRtrol	<ul style="list-style-type: none"><li>OBR - Observation Request Segment</li><li>OBX - Observation Result Segment</li></ul>
New comments list updated	Suspected Pathologies

## 2. HL7 Format

### 2.1. Protocol Description

#### 2.1.1. Overview

This document is intended as a guide for software developers responsible for creating the interface between a hospital (or laboratory) information system (LIS/LIMS/EMR) application and Yumizen P8000. It is assumed that these software developers are familiar with the HL7 standard and have the HL7 specification documents available for reference.

#### Message structure

The following table briefly describes concepts used when describing the HL7 high level protocol. For further details refer to the original HL7 standard specification.

Concept	Definition
Message	A complete, self-contained entity of data. An example of a message is a complete patient test result including patient identification, order information, parameter values and error messages.
Segment	A message is composed with segments each containing related elements of data (attributes). Examples of segments are the patient information segment keeping all the patient data that is common to all tests and the order segment keeping data that is common for the individual test.

Concept	Definition
Field	Each segment has a number of fields each holding one or more data elements (attributes). For instance, the patient information segment has a field containing the patients name and a field holding the patients birth date.
Component field	A field may be divided into several component fields. The name field of the patient information segment has the components last name, first name and middle initials.

Messages consist of segments of various types that are listed in the table below:

Segment Type	Name
MSH	Message Header Segment
MSA	Message Acknowledgement Segment
PID	Patient Identification Segment
PV1	Patient Visit Segment
ORC	Common Request Segment
TQ1	Timing/Quantity Segment
OBR	Observation Request Segment
SPM	Specimen Segment
OBX	Observation Result Segment
NTE	Notes and Comments Segment
SAC	Specimen Container Segment

To report results Yumizen P8000 sends messages to the LIMS as a sequence of segments.

#### Delimiters

Delimiters are used to separate the segment into fields and field components. Delimiters may vary from implementation to implementation, and are defined as part of the header segment.

The following delimiters are used in Yumizen P8000:

Delimiter	Name	Dec. Code	Hex. Code
	Field Delimiter	124	7C
~	Repeat field delimiter	126	7E
^	Component field delimiter	94	5E
&	Sub component delimiter	38	26
	Escape character	91	5C

## Fields

A field is a string of characters.

HL7 does not take into account how systems store data within an application.

When fields are transmitted, they are sent as character strings.

Whether the field is required, optional, or conditional in a segment, this information is provided in the column labeled OPT.

The designations are:

Designation	Description
R	Required
RE	Required but may be Empty: The field or data type component description must stipulate when the field or data type component may be empty.
O	Optional
C	Conditional on the triggered event or on some other field(s). The field definitions following the segment attribute table should specify the algorithm that defines the conditionality for this field.
X	Not used with this triggered event
B	Left in for backward compatibility with previous versions of HL7. The field definitions following the segment attribute table should denote the optionality of the field for prior versions.
W	Withdrawn
CE	Conditional but it may be empty

## Date Format

Dates are always sent as: YYYYMMDD.

Times are always sent as: HHMMSS.

Dates and Times together are sent as: YYYYMMDDHHMMSS. When Yumizen P8000 receives a Date or Date/Time, the following formats are allowed: YYYYMMDD, YYYYMMDDHHMM and YYYYMMDDHHMMSS.

Decimal values are transmitted with a “.” character as the decimal separator.

### 2.1.2. Interface Description

Yumizen P8000 needs to populate its own database with data regarding samples worked by laboratory instruments.


This information is usually received from an external system.

The Yumizen P8000 HL7 interface provides a correct data flow from that system to Yumizen P8000 database.

- The external system must encapsulate requests information in an HL7 message and send it to Yumizen P8000 HL7 interface. This process is called *requests downloading* in the Yumizen P8000 worklist.
- Yumizen P8000 is capable to send results regarding samples to an external system when result information is ready. The Yumizen P8000 HL7 interface encapsulates sample results data in HL7 messages and sends the messages to the external system. This process is called *results uploading* in the Yumizen P8000 workflow.
- The charset settings available are:
  - UTF-8 (used by default)
  - windows-1252
  - ISO-8859-15
- LIS must communicate with two lines (the connection remains after each result sent to LIS):
  - Port 10001 as client for requests sent from LIS to Yumizen P8000. The communication closes after the end of the message.
  - Port 10002 as server for results sent from Yumizen P8000 to LIS.

### 2.1.3. Automation Manager Scenario

This integration scenario applies when Yumizen P8000 plays the role of the Automation Manager. The connected external system (usually a LIS) is an Order Filler, as shown in the following diagram.

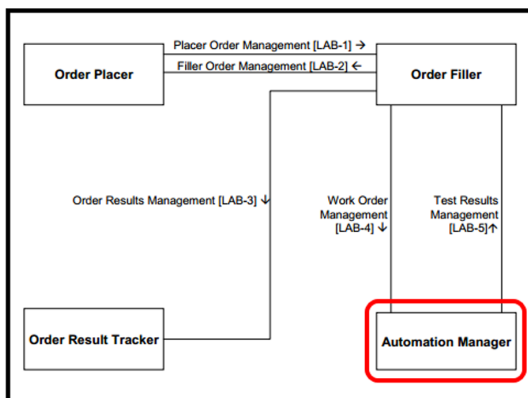

The transactions supported for this scenario are:

- The external system sends orders to Yumizen P8000 (requests downloading): the interaction is similar to one described in the transaction “*Work Order Management (LAB-4)*” of IHE Laboratory Technical Framework. As stated in the IHE document “*This transaction is used if the Order Filler issues a new order to the Automation Manager*”.
- Yumizen P8000 sends results to the external system (results uploading): the interaction is similar to one described in the transaction “*Test Results Management (LAB-5)*” of IHE Laboratory Technical Framework. As stated in the IHE document “*This transaction is used when Automation Manager transmits test results to Order Filler*”.

### 2.1.4. Message Structure

The following table lists the message types supported by Yumizen P8000 for the laboratory testing workflow transactions:

Message	Description
OML^O33 Laboratory order for multiple order related to a single specimen	Specimen oriented order information used to send request from an external system to Yumizen P8000 (requests downloading).
ORL^O34 Laboratory order response	Application Acknowledgement/Rejection of a request message send from Yumizen P8000 to an external system.
OUL^R22 Unsolicited Specimen Observation Message	Specimen oriented result values used to send results from Yumizen P8000 to an external system (results uploading).
ORU^R01 Unsolicited Observation Message	Specimen oriented QC result values used to send QC results from Yumizen P8000 to an external system (QC results uploading).
ACK Generic Acknowledgement Message	Commit Acknowledgement/Rejection of a message used for both request and result message.

#### 2.1.4.1. OML^O33^OML\_O33

The following table lists the detailed structure for message OML^O33 used to send request from an external system to Yumizen P8000:

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
[	--- PATIENT begin	R	[1..1]
PID	Patient Identification	R	[1..1]
PV1	Patient Visit	RE	[0..1]
]	--- PATIENT end		
{	--- SPECIMEN begin	R	[1..*]

Segment	Meaning	Usage	Card.
SPM	Specimen	R	[1..1]
{	--- ORDER begin		
ORC	Common Order	R	[1..1]
[[TQ1]]	Timing Quantity	R	[0..1]
[	--- OBSERVATION_REQUEST begin	O	[0..*]
OBR	Observation Request	R	[1..1]
[[OBX]]	Observation Request Result (Clinical Info)	O	[0..*]
{{	--- PRIOR_RESULT begin	O	[0..*]
PV1	Patient Visit - previous result	R	[1..1]
[ ORC ]	Common Order - previous result	R	[1..1]
OBR	Order Detail - previous result	R	[1..1]
[ OBX ]	Observation Result - previous result	R	[1..*]
[[ NTE ]]	Notes and Comments - previous result	C	[0..*]
}}	--- PRIOR_RESULT end		
]	--- OBSERVATION_REQUEST end		
}	--- ORDER end		
}	--- SPECIMEN end		

#### 2.1.4.2. OUL^R22^OUL\_R22

The following table lists the detailed structure for message OUL^R22 used to send results from Yumizen P8000 to an external system.

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
[PID]	Patient Identification	R	[1..1]

Segment	Meaning	Usage	Card.
PV1	Patient Visit	R	[0..1]
{	--- SPECIMEN begin	R	[1..*]
SPM	Specimen information	R	[1..1]
[[OBX]]	Observation Result (for Specimen)	O	[0..*]
{	--- ORDER begin	R	[1..*]
OBR	Observation order	R	[1..1]
ORC	Common order	R	[1..1]
{{	--- RESULT begin	R	[0..*]
OBX	Observation Result	R	[1..1]
[[NTE]]	Notes and Comments	O	[0..*]
}}	--- RESULT end		
}	--- ORDER end		
}	--- SPECIMEN end		

#### 2.1.4.3. ORU^R01

The following table lists the detailed structure for message ORU^R01 used to send QC results from Yumizen P8000 to an external system.

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
{	--- ORDER begin	R	[1..*]
OBR	Observation order	R	[1..1]
{{	--- RESULT begin	R	[0..*]
OBX	Observation Result	R	[1..1]
}}	--- RESULT end		
}	--- ORDER end		

#### 2.1.4.4. ACK / ORL^O34^ORL\_O34

The following table lists the detailed structure for message ACK used to acknowledge message reception by both Yumizen P8000 and an external system:

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
MSA	Message Acknowledgement	R	[1..1]
[ERR]	Error	O	[0..1]

## 2.2. OUL and OML Messages Segments

### 2.2.1. MSH - Message Header Segment

The MSH segment defines the intent, source, destination, and some specifics of the syntax of a message.

The Message Header segment contains general information and identifies the sender. The Message Header segment is always the first record in a transmission.

#### 2.2.1.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>MSH</b>
1	1	R	[1..1]	Field Separator	(Pipe)

SEQ	LEN	OPT	Card.	Element name	Example
2	4	R	[1..1]	Encoding characters	^~\& ^: sub field delimiter ~: repeat sub field delimiter \: ESCAPE sequence &: sub filed component delimiter
3	20	R	[1..1]	Sending Application	<b>LIS</b> NameSpace ID: LIS
4		R	[1..1]	Sending Facility	<b>LIS</b> NameSpace ID: LIS
5		R	[1..1]	Receiving Application	<b>YP8K</b> NameSpace ID: YP8K
6		R	[1..1]	Receiving Facility	<b>YP8K</b> NameSpace ID: YP8K
7	14	R	[1..1]	Date/Time of message	<b>20160416090430</b>
8		X		Security	
9		R	[1..1]	Message Type	<b>OML^O33^OML_O33</b> : request Message code ID: OML Trigger event ID: 033 Message structure ID: OML_033
10		R	[1..1]	Message Control ID	<b>18698910009</b>
11	1	R	[1..1]	Processing ID	<b>P</b> P (Production) D (Debugging) Default setting depending on user profile: ■ Tech: D ■ Others: P
12		R	[1..1]	Version ID	<b>2.5</b>
13		X		Sequence number	

SEQ	LEN	OPT	Card.	Element name	Example
14		X		Continuation Pointer	
15		X		Accept Acknowledgment Type	
16		X		Application Acknowledgment Type	
17		X		Country Code	
18		X		Character Set	UNICODE UTF-8
19		X		Principal language of message	
20		X		Alternate Character Set Handling Scheme	
21		X		Message Profile Identifier	

**Example of a message sent by LIS:**

```
MSH|^~\&|LIS|LIS|YP8K|YP8K|20160416090430||OML^O33^OML_O33|18698910009|P|2.5|||||
```

**2.2.1.2. From Yumizen P8000 to LIS**

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>MSH</b>
1	1	R	[1..1]	Field Separator	(Pipe)
2	4	R	[1..1]	Encoding characters	^~\& ^: sub field delimiter ~: repeat sub field delimiter \: ESCAPE sequence &: sub filed component delimiter
3	20	R	[1..1]	Sending Application	<b>YP8K</b> NameSpace ID: YP8K
4		R	[1..1]	Sending Facility	NameSpace ID : Empty

SEQ	LEN	OPT	Card.	Element name	Example
5		R	[1..1]	Receiving Application	NameSpace ID: Empty
6		R	[1..1]	Receiving Facility	NameSpace ID: Empty
7	14	R	[1..1]	Date/Time of message	<b>20160705100955</b>
8		X		Security	
9		R	[1..1]	Message Type	<b>OUL^R22^OUL_R22</b> : results Message code ID: OUL Trigger event ID: R22 Message structure ID: OUL_R22
10		R	[1..1]	Message Control ID	<b>YP8K20160705100955</b>
11	1	R	[1..1]	Processing ID	<b>P</b> P (Production) D (Debugging) Default setting depending on user profile: ■ Tech: D ■ Others: P
12		R	[1..1]	Version ID	<b>2.5</b>
13		X		Sequence number	
14		X		Continuation Pointer	
15		X		Accept Acknowledgment Type	
16		X		Application Acknowledgement Type	
17		X		Country Code	
18		X		Character Set	UNICODE UTF-8
19		X		Principal language of message	
20		X		Alternate Character Set Handling Scheme	
21		X		Message Profile Identifier	

**Example of a message sent by Yumizen P8000:**

```
MSH|^~\&|YP8K|^|^|^|20160705100955||OUL^R22^OUL_R22|
YP8K20160705100955|P|2.5|||||
```

**2.2.2. PID - Patient Identification Segment**

The PID segment is used by all applications as the primary means of communicating patient identification information. This segment contains permanent patient identifying and demographic information.

**2.2.2.1. From LIS to Yumizen P8000**

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>PID</b>
1		X		Set ID - PID	
2		X		Patient ID	
3	34	R	[1..*]	Patient Identifier List	<b>P0002^^^LIS^PI</b>
	25				ID number: P0002
					Check Digit: empty
					Check Digit Scheme: empty
	3				Assigning authority: LIS
	2				Identifier Type Code: PI
4		X		Alternate Patient ID - PID	

SEQ	LEN	OPT	Card.	Element name	Example
5		O	[1..1]	Patient Name	<b>DOE^JOHN^^</b>
	30				Family Name: DOE
	30				Given Name: JOHN
					Second and Further Given Names or Initials Thereof: empty
				Suffix: empty	
6		X		Mother's maiden name	
7	8	R	[0..1]	Date/Time of Birth	<b>19601206</b> Date of Birth Format = YYYYMMDD
8	1	R	[1..1]	Administrative Sex	<b>M</b> M : Male F: Female U: Unknown
9		X		Patient Alias	
10		X		Race	
11		O	[0..*]	Patient Address	<b>Main Street^^Springfield^NY^65466^USA^ATC1</b>
	120				Street address: Main Street
					Other designation: empty
	100				City: Springfield
	50				State or province: NY
	20				Zip or postal code: 65466
	15				Country: USA
	5				Address type: ATC1
12		X		Country Code	
13	50	O	[0..1]	Phone Number - Home	<b>0033412364567</b>
14		X		Phone Number - Business	
15		X		Primary Language	
16		X		Marital Status	

SEQ	LEN	OPT	Card.	Element name	Example
17		X		Religion	
18		X		Patient Account Number	
19		X		SSN Number Patient	
20		X		Driver License Number Patient	
21		X		Mother's Identifier	
22		X		Ethnic Group	
23		X		Birth Place	
24		X		Multiple Birth Indicator	
25		X		Birth Order	
26		X		Citizenship	
27		X		Veterans Military Status	
28		X		Nationality	
29		X		Patient Death Date and Time	
30		X		Patient Death Indicator	
31	1	R	[0..1]	Identity Unknown Indicator	<b>N</b>
32	2	R	[0..*]	Identity Reliability Code	<b>AL</b>
33		X		Last Update Date/Time	
34		X		Last Update facility	
35		X		Species Code	
36		X		Breed Code	
37		X		Strain	
38		X		Production Class Code	
39		X		Tribal Citizenship	

**Example of a patient identification segment sent by LIS:**

PID||P0002^^^LIS^PI||DOE^JOHN^^|19601206|M||Main  
Street^^Springfield^NY^65466^USA^ATC1||0033412364567||||||||||||||||N|AL

**2.2.2.2. From Yumizen P8000 to LIS**

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>PID</b>
1		X		Set ID - PID	
2		X		Patient ID	
3	34	R	[1..*]	Patient Identifier List	<b>P0002^^^LIS^PI</b>
	25				ID number: P0002
					Check Digit: empty
					Check Digit Scheme: empty
	3				Assigning authority: LIS
	2				Identifier Type Code: PI
4		X		Alternate Patient ID - PID	
5		O	[1..1]	Patient Name	<b>DOE^JOHN^^</b>
	30				Family Name: DOE
	30				Given Name: JOHN
					Second and Further Given Names or Initials Thereof: empty
					Suffix: empty
6		X		Mother's maiden name	
7	8	R	[0..1]	Date/Time of Birth	<b>19601206</b> Date of Birth Format = YYYYMMDD
8	1	R	[1..1]	Administrative Sex	<b>M</b> M: Male F: Female U: Unknown
9		X		Patient Alias	
10		X		Race	

SEQ	LEN	OPT	Card.	Element name	Example
11		O	[0..*]	Patient Address	<b>Main Street^^Springfield^NY^65466^USA^ATC1</b>
	120				Street address: Main Street
					Other designation: empty
	100				City: Springfield
	50				State or province: NY
	20				Zip or postal code: 65466
	15				Country: USA
	5			Address type: ATC1	
12		X		Country Code	
13	50	O	[0..1]	Phone Number - Home	
14		X		Phone Number - Business	
15		X		Primary Language	
16		X		Marital Status	
17		X		Religion	
18		X		Patient Account Number	
19		X		SSN Number Patient	
20		X		Driver License Number Patient	
21		X		Mother's Identifier	
22		X		Ethnic Group	
23		X		Birth Place	
24		X		Multiple Birth Indicator	
25		X		Birth Order	
26		X		Citizenship	
27		X		Veterans Military Status	
28		X		Nationality	
29		X		Patient Death Date and Time	
30		X		Patient Death Indicator	

SEQ	LEN	OPT	Card.	Element name	Example
31	1	RE	[0..1]	Identity Unknown Indicator	<b>Y</b>
32	2	CE	[0..*]	Identity Reliability Code	
33		X		Last Update Date/Time	
34		X		Last Update facility	
35		X		Species Code	
36		X		Breed Code	
37		X		Strain	
38		X		Production Class Code	
39		X		Tribal Citizenship	

**Example of a patient identification segment sent by Yumizen P8000:**

PID|||P0002^^^LIS^PI||DOE^JOHN^^^^^|19601206|M|||Main Street^^Springfield^NY^65466^USA^ATC1|||||ABC123|||||||Y

### 2.2.3. PV1 - Patient Visit Segment

This segment is used to communicate the location and requestor.

#### 2.2.3.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	O		Segment ID	<b>PV1</b>
1		X		Set ID - PV1	
2	1	RE	[1..1]	Patient Class (to categorize patients by site)	<b>N</b> E: Emergency I: Inpatient N: Not applicable U: Unknown

SEQ	LEN	OPT	Card.	Element name	Example
3		X		Assigned Patient Location	
4		X		Admission Type	
5		X		Preadmit Number	
6		X		Prior Patient Location	
7		X		Attending Doctor	
8		X		Referring Doctor	
9		X		Consulting Doctor	
10		X		Hospital Service	
11		X		Temporary Location	
12		X		Preadmit Test Indicator	
13		X		Remission Indicator	
14		X		Admit Source	
15		X		Ambulatory Status	
16		X		VIP Indicator	
17		X		Admitting Doctor	
18		X		Patient Type	
19		X		Visit Number	
20		X		Financial Class	
21		X		Charge Price Indicator	
22		X		Courtesy Code	
23		X		Credit Rating	

SEQ	LEN	OPT	Card.	Element name	Example
24		X		Contract Code	
25		X		Contract Effective	
26		X		Contract Effective	
27		X		Contract Period	
28		X		Interest Code	
29		X		Transfer to Bad Debt Code	
30		X		Transfer to Bad Debt Date	
31		X		Bad Debt Agency Code	
32		X		Bad Debt Transfer Amount	
33		X		Bad Debt Recovery Amount	
34		X		Delete Account Indicator	
35		X		Delete Account Date	
36		X		Discharge Disposition	
37		X		Discharged to Location	
38		X		Diet Type	
39		X		Servicing Facility	
40		X		Bed Status	
41		X		Account Status	
42		X		Pending Location	
43		X		Prior Temporary Location	
44	14	RE	[0..*]	Admit Date/Time	<b>20160416090430</b>
45	14	RE	[0..*]	Discharge Date/Time	<b>20160416090430</b>
46		X		Current Patient Balance	
47		X		Total Charges	
48		X		Total Adjustments	
49		X		Total Payments	
50		X		Alternate Visit ID	
51		X		Visit Indicator	
52		X		Other Healthcare Provider	

**Example of a patient visit segment sent by LIS:**

PV1||N|||||||||||||||||||||||||||||||||||||||||20160416090430|20160416090430

**2.2.3.2. From Yumizen P8000 to LIS**

This segment is used to communicate the location and requestor.

SEQ	LEN	OPT	Card.	Element name	Example
0	3	O		Segment ID	<b>PV1</b>
1		X		Set ID - PV1	
2	1	R	[1..1]	Patient Class (to categorize patients by site)	<b>N</b> E: Emergency I: Inpatient N: Not applicable U: Unknown
3		X		Assigned Patient Location	
4		X		Admission Type	
5		X		Preadmit Number	
6		X		Prior Patient Location	
7		X		Attending Doctor	
8		X		Referring Doctor	
9		X		Consulting Doctor	
10		X		Hospital Service	
11		X		Temporary Location	
12		X		Preadmit Test Indicator	
13		X		Remission Indicator	
14		X		Admit Source	

SEQ	LEN	OPT	Card.	Element name	Example
15		X		Ambulatory Status	
16		X		VIP Indicator	
17		X		Admitting Doctor	
18		X		Patient Type	
19		X		Visit Number	
20		X		Financial Class	
21		X		Charge Price Indicator	
22		X		Courtesy Code	
23		X		Credit Rating	
24		X		Contract Code	
25		X		Contract Effective	
26		X		Contract Effective	
27		X		Contract Period	
28		X		Interest Code	
29		X		Transfer to Bad Debt Code	
30		X		Transfer to Bad Debt Date	
31		X		Bad Debt Agency Code	
32		X		Bad Debt Transfer Amount	
33		X		Bad Debt Recovery Amount	
34		X		Delete Account Indicator	
35		X		Delete Account Date	
36		X		Discharge Disposition	
37		X		Discharged to Location	
38		X		Diet Type	
39		X		Servicing Facility	
40		X		Bed Status	
41		X		Account Status	
42		X		Pending Location	
43		X		Prior Temporary Location	

SEQ	LEN	OPT	Card.	Element name	Example
44	14	RE	[0..*]	Admit Date/Time	20160416090430
45	14	RE	[0..*]	Discharge Date/Time	20160416090430
46		X		Current Patient Balance	
47		X		Total Charges	
48		X		Total Adjustments	
49		X		Total Payments	
50		X		Alternate Visit ID	
51		X		Visit Indicator	
52		X		Other Healthcare Provider	

**Example of a patient visit segment sent by Yumizen P8000:**  
PV1||N|||||||||||||||||||||||||||||||||||||||||20160416090430|20160416090430

## 2.2.4. SPM - Specimen Segment

The intent of this segment is to describe the characteristics of a specimen.  
It contains information about the sample.

### 2.2.4.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	SPM
1	1	R	[1..1]	Set ID - SPM	1
2	16	R	[0..1]	Sample ID	201604163002
3		X		Specimen Parent IDs	
4	15	R	[1..1]	Specimen Type	BLOOD or BODY FLUID (editable)
5		X		Specimen Type Modifier	

SEQ	LEN	OPT	Card.	Element name	Example
6		X		Specimen Additives	
7		X		Specimen Collection Method	
8	20	RE	[0..1]	Specimen Source Site	MAIN LAB
9		X		Specimen Source Site Modifier	
10		X		Specimen Collection Site	
11		X		Specimen Role	
12		X		Specimen Collection Amount	
13		X		Grouped Specimen Count	
14		X		Specimen Description	
15		X		Specimen Handling Code	
16		X		Specimen Risk Code	
17	14	R	[0..1]	Specimen Collection Date/Time	201604160904
18	14	R	[0..1]	Specimen Received Date/Time	201604160904
19		X		Specimen Expiration Date/Time	
20		X		Specimen availability	
21		X		Specimen Reject Reason	
22		X		Specimen Quality	
23		X		Specimen Appropriateness	
24		X		Specimen Condition	
25		X		Specimen Current Quantity	
26		X		Number of Specimen Containers	
27		X		Container Type	
28		X		Container Condition	
29		X		Specimen Child Role	

**Example of a specimen segment sent by LIS:**  
SPM|1|201604163002||BLOOD||||MAIN LAB|||||||201604160904|201604160904||||

### 2.2.4.2. From Yumizen P8000 to LIS

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>SPM</b>
1	1	R	[1..1]	Set ID - SPM	<b>1</b>
2	16	R	[0..1]	Sample ID	<b>201604163002</b>
3		X		Specimen Parent IDs	
4	15	R	[1..1]	Specimen Type	<b>BLOOD or BODY FLUID</b> (editable)
5		X		Specimen Type Modifier	
6		X		Specimen Additives	
7		X		Specimen Collection Method	
8	20	RE	[0..1]	Specimen Source Site	<b>MAIN LAB</b>
9		X		Specimen Source Site Modifier	
10		X		Specimen Collection Site	
11		X		Specimen Role	
12		X		Specimen Collection Amount	
13		X		Grouped Specimen Count	
14		X		Specimen Description	
15		X		Specimen Handling Code	
16		X		Specimen Risk Code	
17	14	X	[0..1]	Specimen Collection Date/Time	
18	14	X	[0..1]	Specimen Received Date/Time	
19		X		Specimen Expiration Date/Time	
20		X		Specimen availability	
21		X		Specimen Reject Reason	
22		X		Specimen Quality	
23		X		Specimen Appropriateness	
24		X		Specimen Condition	
25		X		Specimen Current Quantity	
26		X		Number of Specimen Containers	

SEQ	LEN	OPT	Card.	Element name	Example
27		X		Container Type	
28		X		Container Condition	
29		X		Specimen Child Role	

#### Example of a specimen segment sent by Yumizen P8000:

SPM|1|201604163002||BLOOD|||MAIN LAB

### 2.2.5. OBR - Observation Request Segment

In the reporting of clinical data, the OBR serves as the report header. It identifies the observation set represented by the following atomic observations. It includes the relevant ordering information when that applies. It contains many of the attributes that usually apply to all of the included observations.

The Universal Service Identifier field corresponds to any parameters or compatible panels:

CBC
DIF
RET
SLIDE
CBF
RBC_PLTO
ESR

- To request a CBR, send CBC in one group of ORC, TQ1 and OBR frames and RET in another group of ORC, TQ1 and OBR frames.
- To request a DIR, send CBC in one group of ORC, TQ1 and OBR frames, a DIF in another group of ORC, TQ1 and OBR frames and RET in another group of ORC, TQ1 and OBR frames.

2.2.5.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>OBR</b>
1	1	R	[0..1]	Set ID - OBR	<b>1</b>
2	20	R	[0..1]	Placer Order Number (RequestID)	<b>L604163002</b>
3	20	R	[0..1]	Filler Order Number (RequestID)	<b>L604163002</b>
4	31	R	[1..1]	Universal Service Identifier	<b>CBC^CBC profile^YP8K</b>
					Code: CBC
					Name: CBC profile
					Sending application: YP8K
5		X		Priority	
6		X		Requested Date/Time	
7		X		Observation Date/Time #	
8		X		Observation End Date/Time #	
9		X		Collection Volume	
10		X		Collector Identifier	
11		R		Specimen Action Code	<b>BLOOD or BODY FLUID</b>
12		X		Danger Code	
13		X		Relevant Clinical Information	
14		X		Specimen Received Date/Time	
15		X		Specimen Source	
16	20	O	[0..1]	Ordering Provider	<b>DR HOUSE</b>
17		X		Order Callback Phone Number	
18		X		Placer Field 1	
19		X		Placer Field 2	
20		X		Filler Field 1	
21		X		Filler Field 2	

SEQ	LEN	OPT	Card.	Element name	Example
22		X		Results Rpt/Status Chng Date/Time	
23		X		Charge to practice	
24		X		Diagnostic Serv Sect ID	
25	1	R	[1..1]	Results Status	<b>P</b>
26		X		Parent Result	
27		X		Quantity/Timing	
28		X		Result Copies To	
29		X		Parent	
30		X		Transportation Mode	
31		X		Reason for Study	
32		X		Principal Results Interpreter	
33		X		Assistant Results Interpreter	
34		X		Technician	
35		X		Transcriptionist	
36		X		Scheduled Date / Time	
37		X		Number of Sample Container	
38		X		Transport Logistics of Collected Sample	
39		X		Collector Comment	
40		X		Transport Arrangement Responsibility	
41		X		Transport Arranged	
42		X		Escort Required	
43		X		Planned Patient Transport Comment	
44		X		Procedure Code	
45		X		Procedure Code Modifier	
46		X		Placer Supplemental Service Information	

SEQ	LEN	OPT	Card.	Element name	Example
47		X		Filler Supplemental Service Information	
48		X		Medically Necessary Duplicate Procedure Reason	
49		X		Result Handling	
50		X		Parent Universal Service Identifier	



SEQ 2 & 3 fields cannot be empty and must be unique.

**Example of an observation result segment sent by LIS:**

OBR|1|L604163002|L604163002|CBC^CBC profile^YP8K|||||BLOOD|||||DR HOUSE|||||P

**2.2.5.2. From Yumizen P8000 to LIS**

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>OBR</b>
1	1	O	[0..1]	Set ID - OBR	<b>1</b>
2	20	RE	[0..1]	Placer Order Number	<b>L604163002</b>
3	20	RE	[0..1]	Filler Order Number	<b>L604163002</b>
4	31	R	[1..1]	Universal Service Identifier	<b>RDW-SD^RDW-SD^YP8K</b> Code: RDW-SD Name: RDW-SD profile Sending application: YP8K
5		X		Priority	
6		X		Requested Date/Time	
7		R		Observation Date/Time #	<b>20160416090400</b>

SEQ	LEN	OPT	Card.	Element name	Example
8		X		Observation End Date/Time #	
9		X		Collection Volume	
10		X		Collector Identifier	
11		X		Specimen Action Code	
12		X		Danger Code	
13		X		Relevant Clinical Information	
14		X		Specimen Received Date/Time	
15		X		Specimen Source	
16	20	RE	[0..1]	Ordering Provider	<b>DR HOUSE</b>
17		X		Order Callback Phone Number	
18		X		Placer Field 1	
19		X		Placer Field 2	
20		X		Filler Field 1	
21		X		Filler Field 2	
22		X		Results Rpt/Status Chng Date/Time	
23		X		Charge to practice	
24		X		Diagnostic Serv Sect ID	
25	1	RE	[1..1]	Results Status	<b>F</b>
26		X		Parent Result	
27		X		Quantity/Timing	
28		X		Result Copies To	
29		X		Parent	
30		X		Transportation Mode	
31		X		Reason for Study	
32		RE		Principal Results Interpreter	<b>ruleResult</b>
33		X		Assistant Results Interpreter	
34		X		Technician	

SEQ	LEN	OPT	Card.	Element name	Example
35		X		Transcriptionist	
36		X		Scheduled Date / Time	
37		X		Number of Sample Container	
38		X		Transport Logistics of Collected Sample	
39		X		Collector Comment	
40		X		Transport Arrangement Responsibility	
41		X		Transport Arranged	
42		X		Escort Required	
43		X		Planned Patient Transport Comment	
44		X		Procedure Code	
45		X		Procedure Code Modifier	
46		X		Placer Supplemental Service Information	
47		X		Filler Supplemental Service Information	
48		X		Medically Necessary Duplicate Procedure Reason	
49		X		Result Handling	
50		X		Parent Universal Service Identifier	

**Example of an observation result segment from:**

```
OBR|1|L604163002|L604163002|RDW-SD^RDW-SD^YP8K|||20160416090400|||||||DR HOUSE|||||||F||||||ruleResult
```

**Example of a comment linked to the sample:**

```
OBR|1|000004|000004|ORDER_COMMENT^Request comment^HALIA|||20191023104806|||||||Dr Dree|||||||F||||||
```

### 2.2.5.3. OBR - Segment Specificities for Order Comments

OBR segment is different depending on type of information sent:

- order comments
- test code results

The **OBR 4** field contains different data according to the type of information sent.

Order comments are always sent first, then the test codes values are sent.

First order comment	OBR 1 0007954944 0007954944 ^Commentaire de la demande^HALIA   20221008085830       INC       F<CR>
	ORC SC 0007954944 0007954944 0007954944 A   20221008101514<CR>
	TQ1       20221008085833  S<CR>
	OBX 1 ST ORDER_COMMENT^Commentaire de la demande   Slide done     F   20221008101514<CR>
Second order comment	OBR 2 0007954944 0007954944 ^Commentaire de la demande^HALIA   20221008085830       INC       F<CR>
	ORC SC 0007954944 0007954944 0007954944 A   20221008101514<CR>
	TQ1       20221008085833  S<CR>
	OBX 1 ST ORDER_COMMENT^Commentaire de la demande  Manual count     F   20221008101514<CR>
First result value	OBR 3 0007954944 0007954944 MCV^VGM^HALIA   20221008085830       INC       F      SARAH<CR>
	ORC SC 0007954944 0007954944 0007954944 A   20221008101547<CR>
	TQ1       20221008085833  S<CR>
	OBX 1 NM MCV^VGM  87.4 fL 82.0 - 98.0     F   20221008091605   Yumizen H1500-1<CR>

## 2.2.6. ORC - Common Order Segment

The Common Order segment (ORC) is used to transmit fields that are common to all orders.

### 2.2.6.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>ORC</b>
1	2	R	[1..1]	Order Control	<b>NW:</b> New order. Event request in OML message sent by the order placer in transaction LAB-1 or in OML message sent by the order filler in transaction LAB-4. <b>PR:</b> Previous values. <b>RF:</b> Rerun
2	20	R	[0..1]	Placer Order Number (RequestID)	<b>L604163002</b>
3	20	R	[0..1]	Filler Order Number (RequestID)	<b>L604163002</b>
4	20	R	[0..1]	Placer Group Number	<b>L604163002</b>
5		X		Order Status	
6		X		Response Flag	
7		X		Quantity/Timing	
8		X		Parent	
9	14	R	[0..1]	Date/Time of Transaction	<b>20160416090430</b>
10		X		Entered By	
11		X		Verified By	
12		X		Ordering Provider	
13		RE		Enterer's Location (WardID)	<b>ward01</b>
14		X		Call Back Phone Number	
15		X		Order Effective Date/Time	
16		X		Order Control Code Reason	

SEQ	LEN	OPT	Card.	Element name	Example
17		X		Entering Organization	
18		X		Entering Device	
19		X		Action By	
20		X		Advanced Beneficiary Notice Code	
21	129	RE		Ordering Facility Name (Ward ID and Wardname)	<b>hematology^^^^^^^^^ward01</b>
	100				Ward Name: hematology
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	20				Ward ID: ward01
22		X		Ordering Facility Address	
23		X		Ordering Facility Phone Number	
24		X		Ordering Provider Address	
25		X		Order Status Modifier	
26		X		Advanced Beneficiary Notice Override Reason	
27		X		Filler's Expected Availability Date/Time	
28		X		Confidentiality Code	
29		X		Order Type	
30		X		Enterer Authorization Mode	
31		X		Parent Universal Service Identifier	



SEQ 2, 3 & 4 fields cannot be empty and must be unique.  
Ward can be automatically learned from Yumizen P8000. ORC.13 and ORC.21 are required in that case.

**Example of a common order segment:**

ORC|NW|L604163002|L604163002|L604163002||||20160416090430||||ward01|||||||  
hematology^^^^^^^^^^^ward01

**2.2.6.2. From Yumizen P8000 to LIS**

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>ORC</b>
1	2	R	[1..1]	Order Control	From Yumizen to LIS: only <b>SC</b> : Status changed
2	20	R	[0..1]	Placer Order Number (RequestID)	<b>L604163002</b>
3	20	R	[0..1]	Filler Order Number (RequestID)	<b>L604163002</b>
4	20	R	[0..1]	Placer Group Number (RequestID)	<b>L604163002</b>
5		R		Order Status	<b>A</b>
6		X		Response Flag	
7		X		Quantity/Timing	
8		X		Parent	
9	14	R	[0..1]	Date/Time of Transaction	<b>20160705100647</b>
10		X		Entered By	
11		X		Verified By	
12		X		Ordering Provider	
13		X		Enterer's Location	
14		X		Call Back Phone Number	
15		X		Order Effective Date/Time	

SEQ	LEN	OPT	Card.	Element name	Example
16		X		Order Control Code Reason	
17		X		Entering Organization	
18		X		Entering Device	
19		X		Action By	
20		X		Advanced Beneficiary Notice Code	
21	129	RE		Ordering Facility Name (Ward ID and Wardname)	<b>hematology^^^^^^^^^^^ward01</b>
	100				Ward Name: hematology
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	1				not used
	20				Ward ID: ward01
22		X		Ordering Facility Address	
23		X		Ordering Facility Phone Number	
24		X		Ordering Provider Address	
25		X		Order Status Modifier	
26		X		Advanced Beneficiary Notice Override Reason	
27		X		Filler's Expected Availability Date/Time	
28		X		Confidentiality Code	
29		X		Order Type	

SEQ	LEN	OPT	Card.	Element name	Example
30		X		Enterer Authorization Mode	
31		X		Parent Universal Service Identifier	

**Example of a common order segment from Yumizen P8000:**

ORC|SC|L604163002|L604163002|L604163002|A||||20160705100647|||||||||^^^^^^^^^^

SEQ	LEN	OPT	Card.	Element name	Example
12		X		Conjunction	
13		X		Occurrence duration	
14		X		Total occurrence's	

**Example of a timing/quantity segment from LIS:**

TQ1|||||20160416090430||S

## 2.2.7. TQ1 - Time/Quantity Segment

The timing/quantity segment contains information about the priority and timing of an order.

### 2.2.7.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>TQ1</b>
1		X		Set ID - TQ1	
2		X		Quantity	
3		X		Repeat Pattern	
4		X		Explicit Time	
5		X		Relative Time and Units	
6		X		Service Duration	
7	14	R	[0..1]	Start date/time	<b>20160416090430</b>
8		X		End date/time	
9	1	R	[0..*]	Priority	<b>S</b> S: urgent R: normal
10		X		Condition text	
11		X		Text instruction	

### 2.2.7.2. From Yumizen P8000 to LIS

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>TQ1</b>
1		X		Set ID - TQ1	
2		X		Quantity	
3		X		Repeat Pattern	
4		X		Explicit Time	
5		X		Relative Time and Units	
6		X		Service Duration	
7	14	R	[0..1]	Start date/time	<b>20160416090430</b>
8		X		End date/time	
9	1	R	[0..*]	Priority	<b>S</b> S: urgent R: normal
10		X		Condition text	
11		X		Text instruction	
12		X		Conjunction	
13		X		Occurrence duration	
14		X		Total occurrence's	

**Example of a timing/quantity segment from Yumizen P8000:**

TQ1|||||20160416090430||S

## 2.2.8. NTE - Notes and Comments Segment

The NTE segment is defined here for inclusion in messages defined in other chapters. It is commonly used for sending notes and comments.

NTE segment is always and exclusively linked to a parameter numerical value.

### 2.2.8.1. From Yumizen P8000 to LIS



The instrument flag linked to a parameter is sent with the possible following values:

- \*: if suspected
- X: if rejected
- Characters string

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>NTE</b>
1		R	[1..1]	Set ID - NTE	<b>1</b>
2		X		Source of Comment	
3		RE	[0..1]	Comment	<b>Lymphoproliferative disorder or viral infection suspicion</b>
4		X		Comment Type	

#### Example of a comment linked to a parameter:

```
OBX|1|NM|HGB^HGB||12.4|g/L|130.0 - 170.0|L~LL||F|||20180827160639|||Yumizen 1
```

```
NTE|1||*<CR>
```

```
NTE|2||Anemia<CR>
```

```
NTE|3||Panic value<CR>
```

## 2.2.9. OBX - Observation Result Segment

The OBX segment is used to transmit a single observation or observation fragment. It represents the smallest indivisible unit of a report.



The LIS can transmit clinical information and/or previous values (according on the ORC.1) using the OBX segment.

### 2.2.9.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>OBX</b>
1		R	[1..1]	Set ID - OBX	<b>1</b>
2	2	C	[0..1]	Value Type	<b>CE</b> NM = Numeric Results CE = Coded Entry SN = Structured Numeric ST = String Data
3	31	RE	[1..1]	Observation Identifier	<b>CLL</b> Code: CLL Name: empty
4		X		Observation Sub-ID	
5	50	C	[0..1]	Observation Value	<b>Clinical comment.</b>
6	20	C	[0..1]	Units	
7		X		References Range	
8		X		Abnormal Flags	
9		X		Probability	
10		X		Nature of Abnormal Test	
11	1	R	[1..1]	Observation Result Status	<b>F</b>

SEQ	LEN	OPT	Card.	Element name	Example
12		X		Effective Date of Reference Range	
13		X		User Defined Access Checks	
14	14	R	[0..1]	Date/Time of the Observation	<b>20160728150751</b>
15		X		Producer's ID	
16		X		Responsible Observer	
17		X		Observation Method	
18		X		Equipment Instance Identifier	
19		X		Date/Time of the Analysis	
20		X			
21		X			
22		X			
23		X		Performing Organization Name	
24		X		Performing Organization Address	
25		X		Performing Organization Medical Director	

For field 3 and 5, the LIS must only send one type of element at the same time: comment code or free text.



- To apply rules on the clinical comments, send a comment code in field 3 and let field 5 empty.
- To have free comment, send free text in field 5 and let field 3 empty.

#### Example for clinical information:

##### Used for rules (with code only)

```
OBX|1|CE|CLL|||||F|||20180827160639|||
```

##### Used without rules (freetext only)

```
OBX|1|CE||Clinical information|||||F|||20180827160639|||
```

#### Example for previous value:

```
ORC|PR|L503111235|L503111235|L503111235|||||20150311110927||||ward01|||||
hematology^^^^^^^^^^^ward01<CR>
OBR|1|L503111235|L503111235|WBC^^P8000|||||BLOOD|||||DR HOUSE|
0033412364566|||||P<CR>
OBX|1|ST|WBC||35|s|||||F|||20150311121323|||||<CR>
```

#### 2.2.9.2. From Yumizen P8000 to LIS



The OBX segment can be used to transmit a result and/or a comment linked to the sample.

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Sement ID	<b>OBX</b>
1		R	[1..1]	Set ID - OBX	<b>1</b>
2	2	C	[0..1]	Value Type	<b>NM</b> NM = Numeric Results CE = Coded Entry SN = Structured Numeric ST = String Data ED = Encapsulated Data
3	31	R	[1..1]	Observation Identifier	<b>RDW-SD^RDW-SD</b> Code: RDW-SD Name: RDW-SD ORDER_COMMENT^Request comment
4		X		Observation Sub-ID	
5	50	C	[0..1]	Observation Value	<b>45.0</b> Delta-Run ^IM^PNG^Base64^Hexadecimal_Image
6	20	C	[0..1]	Units	<b>fl</b>

SEQ	LEN	OPT	Card.	Element name	Example
7		X		References Range	
8		X		Abnormal Flags	
9		X		Probability	
10		X		Nature of Abnormal Test	
11	1	RE	[1..1]	Observation Result Status	<b>F</b>
12		X		Effective Date of Reference Range	
13		X		User Defined Access Checks	
14	14	R	[0..1]	Date/Time of the Observation	<b>20160705100630</b>
15		X		Producer's ID	
16		X		Responsible Observer	
17		X		Observation Method	
18		C	[0..1]	Equipment Instance Identifier	<b>Yumizen H2500-SPS</b>
19		X		Date/Time of the Analysis	
20		X			
21		X			
22		X			
23		X		Performing Organization Name	
24		X		Performing Organization Address	
25		X		Performing Organization Medical Director	

**Example of a generic observation/result segment from Yumizen P8000:**  
 OBX|1|NM|RDW-SD^RDW-SD||45.0|f|||||F|||20160705100630|||Yumizen H2500-SPS

**Example of a comment linked to the sample:**

```
OBR|1|34534534589|34534534589|ORDER_COMMENT^Request comment^P8000|||
20191113110805|||||||||||||F|||||<CR>
ORC|SC|34534534589|34534534589|34534534589|A|||||20191114114332|||
WardCode|||||Wardname^^^^^^^^^WardCode<CR>
TQ1|||||20191114110805||R<CR>
OBX|1|ST|ORDER_COMMENT^Request comment||Delta-Run MPV|||||F|||
20191023144806|||<CR>
```



Images are also sent using PNG format with base64 encoding.

Code (see OBX.3)	Description
RBC	Red blood cells histogram
PLT	Platelets histogram
BAS%	White blood cells histogram
NEU%	LMNE matrix
PLT-Ox	Optical platelets matrix
BFMN%	Body fluid matrix
BFRBC	Body fluid red blood cells histogram

**Example of a BASO curve from Yumizen P8000:**

```
OBX|1|ED|BAS%||
^IM^PNG^Base64^iVBORw0KGgoAAAANSUheUgAAAQQAAB2CAIAAACh93ojAAAH
mkIEQVR42u2d208UVxyA980H46sPPHjm0/1T8CYiNwUixclLgi7Ui/
V1jatynJZUcplpVZl1WJLja3bsq7YqnRBI1CsXKUgKJRSbWLUB2LbTTRFJbo9YXDcrrhsl
Zk5u35fJpthmLjMmd/n7/xmzpxw
+QFgHBNNAIAMAMgAgAwAyACADADIAIAMAMgAgAwAyAAgkwWlpWVfHneFucTExIS/
M0uELgaeZfHVzc3NhsLw0Hnk/KWuMBeTyRT+zhaLxRpAQkKC9b+85pbE5SvC/
2NYNDrL0/7VRspQWlauRTM1d/QWFhZqmhMLi4q/
Pf0jsYsM0yaD3W7XopkudfVpLcPwjZub3/+12EUG2TND59VBrWUQxCUuJ3aRYdpk2Lf/
oBbN1PbLtT179mgtw9ZtH567eJnwRQapM0NL2xWn06m1DG5PrWP/
YcIXGaSuGerON7pcLq1l+POvv7NyNhG+yCB1Zjj1/Tmv16vDfZb4ZcmELzJIXTMc
```

```
+9rV3t6ugwzWdzYTvsaggdWaoPHRkaGhIBxk2bN5K
+CKD1DVDSWn5yMiIDjLk23efPHueCEYGeTODLS9/
bGxMBxmOVh87VPONEyWm8tYMO3Nz9Rmo2HChcVdJBRGMDBJnBluePjIM3/
hJywfbIWBkklDmKCgoOEeGf0YfmrNziGBkkDczWCwWfWR4+vSpOWs9EYwM8tYMIuHo
9nzTUobrIYPMmUFPGegmIYPUNQMylAOZARmQgZoBGZCBzGC4DMwMgAzUDBMw
MwAykBNUm9DMDIAM1AzPYGYAZCAzTMDMAMhAzTABMwMgA5lhAmYgQAZqhufEJ
zEzADKQGcaxbmBmAGSgZhgNr7DlfaaBUEYG6TKD1WrVWQZxjcdx4AihjAwh9zCZ9K8
ZdJh1Olih4RvvbvUyUEaGV5dBi8ygw8sZXuTRo8dr11kJZWR4Hvfi35o/f/
7MmTORq6vVjXV1dWLjjBkzpkKS7t27p3XNoMPLGV5yQYkRe8gQIENKSrt27eHhoZmz5
6tbty4cePdu3fFutrzcJ0TozdPQO6F9AC9ZzQQkZAmW4f/
9+UO8oqJsOZ84crWuGlo4eh8Ohvzw59t1cUEKGSqDI8kgelCBP1bs2z/tzXSxpbWqqkp/
Gwo8tVxQQoYpZOjv71e3z507V+ua4Yz3Qk1Njf4y/Pb7TS4olcMUMsybN294eFis9/
T0iKJC65rBXftDfX29/
jI8fjzGBSVkCCXDrFmz3G73ggULRAcpLi7u1q1bWtcMXx0/0dbW5jeCRN5dggxS3YE
+9PnRgYEBQ2TwnDr9ScUBAhoZZBmbVOaouHPnjt8gkpYlx8YIEtPIIEVmsO8qevDggVE
y+Hy
+txYuZIoAZJDieYbtO3b6DYUpApBBIsyQa7P5jYyApBBippBtzeVhIapApBBisyg25tKQs
AUAcggRc2g/
5M9L8lUAchgfGYw5GGGSUngBhwyGFsztHZfk0QG5mBFBomZQ0fvQFFRkQwyclEVGQ
yuGX5q7Xl6nX45IE
+kp4QMxmWGs94LbrdbEhnoKSGDkTWDq8bzOsdDTwkZoiczHK764vr1635piB+/
Fd3V92tmZibv+EEGXWuGktLykZEREWRQBnWLz/r6egZ4I4OumWFnu3Jkyd
+mUhanpyWlqasx8YnEe7loFPNIMMovSDMZrPP52PMEjLonRny8/
P9EsOYJWtQqWZo7ui1WCwyy8CYJWtQKTNO9A7IMGR1iutLvNYEGXSoGvVrarIRWvko
uw6Yt7zW0dBL0yKbtZnDVeJqamiSX4bODzuoTHoleGbStGcrKHVLdZJiUzq7uHQW7CX
rZZGhq74mqzJCba/
NLz8OHj8xZ6w162WTovDoYPTXDpa6+7OxsfySwNIHpAqST4ef27ujJDN19g8XFxREhgzk
7RxmtZLVaGbAkiQzfnayNnprB7altbGyMFBmU0Up
+ZqSURoZyx97oyQyFdvvo6GhEyCASgkgLzwcVLUtuae00m81i+4q3V9Z6G4MOrad/
8M2ZstloGWy2vCipGUSvlyMjwx+Z+Hy+RysWKeOXxGdsXII4InVC78tX+tPT03fY8t
+Q54QMkUE0clZWWpRkBrXXEQUoHafVa9Yo3afWju7Kyso35zkhQ2RQGjkaagbRi1DHS
EchSm0tuk+iy7S34INIUnFi+Lc4WNGhUnpTynpg58pisbyr4UMIRalkSM
+M4gO0uLFi9Ux0tGEOCgR6+r1Yus4yoBw5VdLlixRO1eKfcoEasqPQdepluVpu3BkelWX
wojDj42NnfTwx9Ec/mfvbs5gmuGaOogadHdCmoo+S9eibM85V
+o9iHDX0lcvhpcBSuQnW3JXJe11mwW/RzxmZqapq4r21NSUtT1mJiYf/
dJT18rjBIZRvS4xCfrgesZGZmrVq9W2iojM1PZvnLVqsA2F+VqPup2sSX0PoHnJZxzF/
4+4iylcxp6H+W8/6/vUuNEtEnQPmpqUq7ia9+nTelm/
gPGAAZAJABABkAkAEAGQCQAQAZAJABABkAkAEAGQCQAQAZAJABABkAkAEAGQ
```

```
CQAQAZAJABABkAkAEAGQCQAQAZAJABABkAkAEgevgX4Tpp/u
+X2ikAAAAASUVORK5CYII=|||N||F|||20170328101027|||702M2SH00019
```

## 2.2.10. MSA - Message Acknowledgment Segment

The MSA segment contains information sent while acknowledging another message.

### Acknowledgment codes explanation

Value	Description
AA	Application Accept - Enhanced mode: Application acknowledgment: Accept
AE	Application Error - Enhanced mode: Application acknowledgment: Error
AR	Application Reject - Enhanced mode: Application acknowledgment: Reject

### 2.2.10.1. From LIS to Yumizen P8000

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>MSA</b>
1	2	R	[1..1]	Acknowledgment Code	<b>AA</b>
2		R	[1..1]	Message Control ID	<b>YP8K20160705100955</b>
3		X		Text Message	
4		X		Expected Sequence Number	
5		X		Delayed Acknowledgment Type	
6		X		Error Condition	

**Example of a message acknowledgment segment from LIS:**  
MSA|AA|YP8K20160705100955

### 2.2.10.2. From Yumizen P8000 to LIS

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>MSA</b>
1	2	R	[1..1]	Acknowledgment Code	<b>AA</b>
2		R	[1..1]	Message Control ID	<b>18698910009</b>
3		X		Text Message	
4		X		Expected Sequence Number	
5		X		Delayed Acknowledgment Type	
6		X		Error Condition	

**Example of a message acknowledgment segment from Yumizen P8000:**  
MSA|AA|18698910009

### 2.2.11. SAC - Specimen Container Segment

#### 2.2.11.1. From Yumizen P8000 to Host

The goal of this segment is to inform that the sample is currently in analysis inside the instrument.

To configure this segment, refer to the RAS1117 - Configuring the LIS Connection procedure.

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	SAC
7	14	R		Event date	20170725121422

<sup>a</sup>: This field is configurable by user.

SEQ	LEN	OPT	Card.	Element name	Example
8		R		Event code <sup>a</sup>	CHECK_IN
		R		Event description <sup>a</sup>	CheckIn / Manual CheckIn
		R		User or device that triggered the event	YH2500-YH2500
10	6	C		Rack number	748565
11	2	C		Position in rack	7

<sup>a</sup>: This field is configurable by user.

#### Example of an automatic check in:

SAC|||||20220516143527|CHECK\_IN^CheckIn^YH2500 - YH2500||748265|7

#### Example of a manual check in:

SAC|||||20220613103833|CHECK\_IN^Manual CheckIn^SUPPORT

### 2.2.12. Special Characteristics for HORIBA Medical Data

#### 2.2.12.1. Parameters

LOINC Code: Logical Observation Identifiers Names & Codes

CBC Codes	LOINC Code	Definition
<b>RBC</b>	789-8	Red Blood Cells
<b>HGB</b>	718-7	Hemoglobin Concentration
<b>HCT</b>	4544-3	Hematocrit
<b>MCV</b>	787-2	Mean Corpuscular Volume
<b>MCH</b>	785-6	Mean Corpuscular Hemoglobin
<b>MCHC</b>	786-4	Mean Corpuscular Hemoglobin Concentration
<b>RDW-SD</b>	21000-5	Red Distribution Width Standard Deviation

CBC Codes	LOINC Code	Definition
<b>RDW-CV</b>	788-0	Red Distribution Width
<b>MIC</b>	X-MIC	Microcytic Red Blood Cells percentage (versus RBC)
<b>MAC</b>	X-MAC	Macrocytic Red Blood Cells percentage (versus RBC)
<b>PLT</b>	777-3	Platelets
<b>PCT</b>	51637-7	Plateletcrit
<b>PDW</b>	51631-0	Platelets Distribution Width
<b>MPV</b>	32623-1	Mean Platelet Volume
<b>P-LCC</b>	96354-6	Platelets - Large Cell Count
<b>P-LCR</b>	48386-7	Platelets - Large Cell Ratio
<b>WBC</b>	6690-2	White Blood Cells
<b>NRBC#</b>	771-6	Nucleated Red Blood Cells absolute value
<b>NRBC%</b>	58413-6	Nucleated Red Blood Cells percentage
<b>TNC</b>	50774-9	Total Nucleated cells

RBC_PLTO Codes	LOINC Code	Definition
<b>PLT-Ox</b>	97995-5	Platelets from optical channel
<b>LPF</b>	97994-8	Large Platelet Fraction

DIF Codes	LOINC Code	Definition
<b>LYM#</b>	731-0	Lymphocytes absolute value
<b>LYM%</b>	736-9	Lymphocytes percentage
<b>MON#</b>	742-7	Monocytes absolute value
<b>MON%</b>	5905-5	Monocytes percentage
<b>NEU#</b>	751-8	Neutrophils absolute value
<b>NEU%</b>	770-8	Neutrophils percentage
<b>EOS#</b>	711-2	Eosinophils absolute value
<b>EOS%</b>	713-8	Eosinophils percentage
<b>BAS#</b>	704-7	Basophils absolute value
<b>BAS%</b>	706-2	Basophils percentage

DIF Codes	LOINC Code	Definition
<b>IMG#</b>	53115-2	Immature Granulocytic cells absolute value
<b>IMG%</b>	71695-1	Immature Granulocytic cells percentage
<b>IMM#</b>	X-IMM#	Immature Monocytic cells absolute value
<b>IMM%</b>	X-IMM%	Immature Monocytic cells percentage
<b>IML#</b>	X-IML#	Immature Lymphocytic cells absolute value
<b>IML%</b>	X-IML%	Immature Lymphocytic cells percentage
<b>ALY#</b>	43743-4	Atypical Lymphocytes absolute value
<b>ALY%</b>	42250-1	Atypical Lymphocytes percentage
<b>LIC#</b>	55432-9	Large Immature Cells absolute value
<b>LIC%</b>	55433-7	Large Immature Cells percentage

RET Codes	LOINC Code	Definition
<b>RET#</b>	14196-0	Reticulocytes absolute value
<b>RET%</b>	17849-1	Reticulocytes percentage
<b>RET-L</b>	X-RET-L	Reticulocytes with a low RNA content
<b>RET-M</b>	X-RET-M	Reticulocytes with a medium RNA content
<b>RET-H</b>	X-RET-H	Reticulocytes with a high RNA content
<b>CRC</b>	X-CRC	Corrected Reticulocyte Count
<b>MRV</b>	48706-6	Mean Reticulocyte Volume
<b>RHCC</b>	X-RHCC	Reticulocyte Hemoglobin Cellular Content
<b>IRF</b>	X-IRF	Immature Reticulocyte Fraction

Extended DIF Codes	LOINC Code	Definition
<b>Blast#</b>	X-BLAST#	Blasts absolute value
<b>Blast%</b>	X-BLAST%	Blasts percentage
<b>Other#</b>	X-OTHER#	Other cells absolute value
<b>Other%</b>	X-OTHER%	Other cells percentage
<b>Band Cel#</b>	X-BANDNEU#	Bands absolute value
<b>Band Cel%</b>	X-BANDNEU%	Bands percentage
<b>VarLym#</b>	X-VARLYM#	Variant lymphocytes absolute value

Extended DIF Codes	LOINC Code	Definition
<b>VarLym%</b>	X-VARLYM%	Variant lymphocytes percentage
<b>LGLym#</b>	X-LGLYM#	Large granular lymphocytes absolute value
<b>LGLym%</b>	X-LGLYM%	Large granular lymphocytes percentage
<b>Hairy cell#</b>	X-HAIRY#	Hairy cells absolute value
<b>Hairy cell%</b>	X-HAIRY%	Hairy cells percentage
<b>Sezary Cel#</b>	X-SEZARY#	Sezary cells absolute value
<b>Sezary Cel%</b>	X-SEZARY%	Sezary cells percentage
<b>Plasma cell#</b>	X-PLASMA#	Plasma cells absolute value
<b>Plasma cell%</b>	X-PLASMA%	Plasma cells percentage
<b>ProLympho#</b>	X-PROLYM#	Prolymphocytes absolute value
<b>ProLympho%</b>	X-PROLYM%	Prolymphocytes percentage
<b>Promonocyte#</b>	X-PROMON#	Promonocytes absolute value
<b>Promonocyte%</b>	X-PROMON%	Promonocytes percentage
<b>Promyelo#</b>	X-PROMYE#	Promyelocytes absolute value
<b>Promyelo%</b>	X-PROMYE%	Promyelocytes percentage
<b>Myelocyte#</b>	X-MYELO#	Myelocytes absolute value
<b>Myelocyte%</b>	X-MYELO%	Myelocytes percentage
<b>Metamyelo#</b>	X-METAMYE#	Metamyelocytes absolute value
<b>Metamyelo%</b>	X-METAMYE%	Metamyelocytes percentage
<b>IMMEos#</b>	X-IMEOS#	Immature eosinophils absolute value
<b>IMMEos%</b>	X-IMEOS%	Immature eosinophils percentage
<b>IMMBas#</b>	X-IMBAS#	Immature basophils absolute value
<b>IMMBas%</b>	X-IMBAS%	Immature basophils percentage
<b>TAG</b>	-	Thrombocyte aggregates
<b>GT</b>	-	Giant thrombocytes
<b>SMU</b>	-	Smudge cells
<b>Artefact</b>	-	Artefact
<b>Not Classed%</b>	-	Not classed

CBF Codes	LOINC Code	Definition
<b>BFRBC</b>	X-BFRBC	Red Blood Cells absolute value
<b>BFWBC</b>	X-BFWBC	White Blood Cells absolute value
<b>BFMN#</b>	X-BFMN#	Mononuclear absolute value
<b>BFMN%</b>	X-BFMN%	Mononuclear percentage
<b>BFPN#</b>	X-BFPN#	Polymorphonuclear absolute value
<b>BFPN%</b>	X-BFPN%	Polymorphonuclear percentage
<b>BFOther%</b>	X-BFOTHER#	Other BF percentage
<b>BFOther#</b>	X-BFOTHER%	Other BF absolute value

ESR Codes	LOINC Code	Definition
<b>ESR</b>	82477-1	Erythrocyte Sedimentation Rate

### 2.2.12.2. Suspected Pathologies

Suspected pathologies are transmitted through one Comment record located after the corresponding Result record. If several pathologies are suspected, they are transmitted through several NTE frames.

Refer to [NTE - Notes and Comments Segment](#).

Here is the full list of comments:

[Full comments list](#)

Here are the new comments added in this version:

[New comments list](#)

## 2.3. Laboratory Testing Workflow Examples

### 2.3.1. LIS to Yumizen P8000 standard test request

```
<VT>
MSH|^~\&|LIS|LIS|YP8K|YP8K|20160416090430||OML^O33^OML_O33|18698910009|P|2.5|||||
<CR>
PID||P0002^^^LIS^PI||DOE^JOHN^^|19601206|M||Main
Street^^Springfield^NY^65466^USA^ATC1||0033412364567||||ABC123^^LIS||||||||N|AL<CR>
PV1||N||||||||||||||||20160416090430|20160416090430<CR>
SPM|1|201604163002||BLOOD|||MAIN LAB|||||201604160904|201604160904||||<CR>
ORC|NW|L604163002|L604163002|L604163002|||||20160416090430|||ward01|||||
hematology^^^^^^^^^ward01<CR>
TQ1|||||20160416090430||S<CR>
OBR|1|L604163002|L604163002|CBC^CBC profile^YP8K|||||DR HOUSE|||||P<CR>
OBX|1|CE|CLL|Clinical comment.||||F||20160728150751||<CR>
ORC|NW|L604163002|L604163002|L604163002|||||20160416090430|||ward01|||||
hematology^^^^^^^^^ward01<CR>
TQ1|||||20160416090430||S<CR>
OBR|1|L604163002|L604163002|DIF^DIF profile^YP8K|||||DR HOUSE|||||P<CR>
<FS>
<CR>
```

### 2.3.2. LIS to Yumizen P8000 standard test request with previous results

```
<VT>
MSH|^~\&|LIS|LIS|YP8K|YP8K|20160416090430||OML^O33^OML_O33|18698910009|P|2.5|||||
<CR>
PID||P0002^^^LIS^PI||DOE^JOHN^^|19601206|M||Main
Street^^Springfield^NY^65466^USA^ATC1||0033412364567||||ABC123^^LIS||||||||N|AL<CR>
PV1||N||||||||||||||||20160416090430|20160416090430<CR>
SPM|1|201604163002||BLOOD|||MAIN LAB|||||201604160904|201604160904||||<CR>
ORC|NW|L604163002|L604163002|L604163002|||||20160416090430|||ward01|||||
hematology^^^^^^^^^ward01<CR>
TQ1|||||20160416090430||S<CR>
OBR|1|L604163002|L604163002|CBC^CBC profile^YP8K|||||BLOOD|||||DR HOUSE|||||P<CR>
OBX|1|CE|CLL|Clinical comment.||||F||20160728150751||<CR>
ORC|NW|L604163002|L604163002|L604163002|||||20160416090430|||ward01|||||
hematology^^^^^^^^^ward01<CR>
TQ1|||||20160416090430||S<CR>
OBR|1|L604163002|L604163002|DIF^DIF profile^YP8K|||||DR HOUSE|||||P<CR>
PV1||N||||||||||||||||20150131110927|20150131110927<CR>
ORC|PR|L503111235|L503111235|L503111235|||||20150311110927|||ward01|||||
hematology^^^^^^^^^ward01<CR>
OBR|1|L503111235|L503111235|WBC^^P8000|||||BLOOD|||||DR HOUSE|0033412364566|||||
P<CR>
OBX|1|ST|WBC||35|s||||F||20150311121323|||||<CR>
<FS>
<CR>
```

### 2.3.3. Yumizen P8000 to LIS acknowledgment

```
<VT>
MSH|^~\&|YP8K|^|^|^|20160705095243||ORL^O34^ORL_O34|YP8K20160705095243|P|2.5|||||
<CR>
MSA|AA|18698910009
```

### 2.3.4. Yumizen P8000 to LIS standard result

```
<VT>
MSH|^~\&|YP8K||LIS||20220330094150||OUL^R22^OUL_R22|18344563693096|P|
2.5<CR>
PID||0002^^^LIS^PI||PATIENT 2^TEST^^^^^|19260607|F||Main
Street^^Springfield^NY^65466^USA^ATC1|||||ABC123|||||Y<CR>
PV1||N|WARD00002^^^|ATD^|||||ADD^|ABC123|||||20220330114144|
20220330114144<CR>
SPM|1|202203300002||BLOOD|||MAIN LAB<CR>
OBR|1|2203300002|2203300002|^Commentaire de la demande^HALIA||
20220330113916|||||F<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330054620|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|ST|ORDER_COMMENT^Commentaire de la demande||Woman: >21|||||F||
20220330054620<CR>
OBR|2|2203300002|2203300002|PCT^PCT^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|NM|PCT^PCT||0.179|%|0.15 - 0.4|||||F|||20220330152729|||||H2500ID<CR>
```

```
NTE|1||*<CR>
OBR|3|2203300002|2203300002|IMG% ^IMG% ^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|NM|IMG% ^IMG% ||0.1|%|0.0 - 2.0|||||F|||20220330114559|||||H2500ID<CR>
OBR|4|2203300002|2203300002|MCV^MCV^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|NM|MCV^MCV||104.1|fL|78.0 - 100.0|H|||||F|||20220330152729|||||H2500ID<CR>
OBR|5|2203300002|2203300002|NEU#^NEU#^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|NM|NEU#^NEU#||2.34|1E09/L|1.5 - 7.0|||||F|||20220330114559|||||H2500ID<CR>
OBR|6|2203300002|2203300002|P-LCR^P-LCR^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|NM|P-LCR^P-LCR||49.4|%|18.0 - 50.0|||||F|||20220330152729|||||H2500ID<CR>
NTE|1||*<CR>
OBR|7|2203300002|2203300002|IMG# ^IMG# ^HALIA|||20220330113916|||||F|||||
SUPPORT<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
```

OBX|1|NM|IMG#^IMG#||0.00|1E09/L|0.0 - 999.9||||F|||20220330114559||||H2500ID<CR>  
OBR|8|2203300002|2203300002|LPF^LPF^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|LPF^LPF||7.3|%|0.0 - 99.9||||F|||20220330152729||||H2500ID<CR>  
OBR|9|2203300002|2203300002|NEU%^NEU%^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|NEU%^NEU%||50.9|%|40.0 - 75.0||||F|||20220330114559||||H2500ID<CR>  
OBR|10|2203300002|2203300002|PLT-Ox^PLT-Ox^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|PLT-Ox^PLT-Ox||131|1E09/L|150.0 - 400.0|L||F|||20220330152729||||  
H2500ID<CR>  
NTE|1||Heterogeneous run<CR>  
OBR|11|2203300002|2203300002|RDW-CV^RDW-CV^HALIA|||  
20220330113916|||||||||||||F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|RDW-CV^RDW-CV||15.8|%|12.0 - 18.0||||F|||20220330152729||||H2500ID<CR>  
NTE|1||\*<CR>  
OBR|12|2203300002|2203300002|TNC^TNC^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>

ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|TNC^TNC||4.59|1E09/L|3.5 - 10.0||||F|||20220330114559||||H2500ID<CR>  
OBR|13|2203300002|2203300002|ALY%^ALY%^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|ALY%^ALY%||0.5|%|0.0 - 2.5||||F|||20220330114559||||H2500ID<CR>  
OBR|14|2203300002|2203300002|ALY#^ALY#^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|ALY#^ALY#||0.02|1E09/L|0.0 - 0.2||||F|||20220330114559||||H2500ID<CR>  
OBR|15|2203300002|2203300002|NRBC#^NRBC#^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|NRBC#^NRBC#||0.00|1E09/L|0.0 - 0.15||||F|||20220330114559||||  
H2500ID<CR>  
OBR|16|2203300002|2203300002|MPV^MPV^HALIA|||20220330113916|||||||||||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MPV^MPV||11.8|fL|7.4 - 12.0||||F|||20220330152729||||H2500ID<CR>  
NTE|1||\*<CR>

**HL7 Format**

Laboratory Testing Workflow Examples



OBR|17|2203300002|2203300002|RBC^RBC^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|RBC^RBC||4.22|1E12/L|3.8 - 5.8|||F|||20220330152729|||H2500ID<CR>  
OBR|18|2203300002|2203300002|P-LCC^P-LCC^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|P-LCC^P-LCC||75|1E09/L|44.0 - 140.0|||F|||20220330152729|||H2500ID<CR>  
NTE|1|\*<CR>  
OBR|19|2203300002|2203300002|MON#^MON#^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MON#^MON#||0.41|1E09/L|0.2 - 0.8|||F|||20220330114559|||H2500ID<CR>  
OBR|20|2203300002|2203300002|PLT^PLT^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|PLT^PLT||152|1E09/L|150.0 - 400.0|||F|||20220330152729|||H2500ID<CR>  
NTE|1|\*<CR>  
NTE|2||SLIDE: Schistocytes suspicion optical PLT<CR>  
NTE|3||Large Platelets suspicion<CR>  
NTE|4||PLT abn. histogram - Large PLT?<CR>  
NTE|5||PLTO abn. matrix - Schistocyte?<CR>

OBR|21|2203300002|2203300002|IML#^IML#^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|IML#^IML#||0.00|1E09/L|0.0 - 0.05|||F|||20220330114559|||H2500ID<CR>  
OBR|22|2203300002|2203300002|WBC^WBC^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|WBC^WBC||4.59|1E09/L|3.5 - 10.0|||F|||20220330114559|||H2500ID<CR>  
OBR|23|2203300002|2203300002|LIC%^LIC%^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|LIC%^LIC%||0.1|0.0 - 3.0|||F|||20220330114559|||H2500ID<CR>  
OBR|24|2203300002|2203300002|MON%^MON%^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MON%^MON%||8.8|4.0 - 13.0|||F|||20220330114559|||H2500ID<CR>  
OBR|25|2203300002|2203300002|LIC#^LIC#^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130|||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|LIC#^LIC#||0.00|1E09/L|0.0 - 0.2|||F|||20220330114559|||H2500ID<CR>

OBR|26|2203300002|2203300002|IML%<sup>^</sup>IML%<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|IML%<sup>^</sup>IML%||0.0%|0.0 - 0.2||||F|||20220330114559||||H2500ID<CR>  
OBR|27|2203300002|2203300002|LYM#<sup>^</sup>LYM#<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|LYM#<sup>^</sup>LYM#||1.73|1E09/L|1.0 - 4.0||||F|||20220330114559||||H2500ID<CR>  
OBR|28|2203300002|2203300002|HGB<sup>^</sup>HGB<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|HGB<sup>^</sup>HGB||156|g/L|120.0 - 160.0||||F|||20220330152729||||H2500ID<CR>  
OBR|29|2203300002|2203300002|PDW<sup>^</sup>PDW<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|PDW<sup>^</sup>PDW||24.6|fL|11.0 - 20.0|H~HH||F|||20220330152729||||H2500ID<CR>  
NTE|1||\*<CR>  
OBR|30|2203300002|2203300002|SLIDE<sup>^</sup>SLIDE<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|ST|SLIDE<sup>^</sup>SLIDE||SLIDE\_EJECTED|-||||F|||20220330152740||||H2500ID<CR>

NTE|1||Slide Count Expected<CR>  
OBR|31|2203300002|2203300002|MIC<sup>^</sup>MIC<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MIC<sup>^</sup>MIC||0.7%|0.0 - 20.0||||F|||20220330152729||||H2500ID<CR>  
OBR|32|2203300002|2203300002|LYM%<sup>^</sup>LYM%<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|LYM%<sup>^</sup>LYM%||37.8%|15.0 - 45.0||||F|||20220330114559||||H2500ID<CR>  
OBR|33|2203300002|2203300002|RDW-SD<sup>^</sup>RDW-SD<sup>^</sup>HALIA|||  
20220330113916|||||F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|RDW-SD<sup>^</sup>RDW-SD||49.0|fL|37.0 - 56.0||||F|||20220330152729||||H2500ID<CR>  
NTE|1||\*<CR>  
NTE|2||RBC abn. histogram - Abnormal distribution<CR>  
OBR|34|2203300002|2203300002|MAC<sup>^</sup>MAC<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MAC<sup>^</sup>MAC||6.8%|2.0 - 10.0||||F|||20220330152729||||H2500ID<CR>  
OBR|35|2203300002|2203300002|BAS%<sup>^</sup>BAS%<sup>^</sup>HALIA|||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^^W00002<CR>

TQ1|||||20220330113916||R<CR>  
OBX|1|NM|BAS% ^BAS%||1.2|0.0 - 2.0|||F||20220330114559||||H2500ID<CR>  
OBR|36|2203300002|2203300002|BAS# ^BAS# ^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|BAS# ^BAS#||0.05|1E09/L|0.0 - 0.2|||F||20220330114559||||H2500ID<CR>  
OBR|37|2203300002|2203300002|NRBC% ^NRBC% ^HALIA||20220330113916|||||F|||||  
F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|NRBC% ^NRBC%||0.0|0.0 - 2.0|||F||20220330114559||||H2500ID<CR>  
OBR|38|2203300002|2203300002|MCH ^MCH ^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MCH ^MCH||37.1|pg|25.0 - 32.0|H~HH|||F||20220330152729||||H2500ID<CR>  
OBR|39|2203300002|2203300002|MCHC ^MCHC ^HALIA||20220330113916|||||F|||||  
F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|MCHC ^MCHC||356|g/L|320.0 - 350.0|H|||F||20220330152729||||H2500ID<CR>  
OBR|40|2203300002|2203300002|HCT ^HCT ^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>

OBX|1|NM|HCT ^HCT||0.439|L/L|0.36 - 0.48|||F||20220330152729||||H2500ID<CR>  
OBR|41|2203300002|2203300002|IMM# ^IMM# ^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|IMM# ^IMM#||0.00|1E09/L|0.0 - 0.1|||F||20220330114559||||H2500ID<CR>  
OBR|42|2203300002|2203300002|EOS# ^EOS# ^HALIA||20220330113916|||||F|||||  
SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|EOS# ^EOS#||0.06|1E09/L|0.0 - 0.5|||F||20220330114559||||H2500ID<CR>  
OBR|43|2203300002|2203300002|EOS% ^EOS% ^HALIA||20220330113916|||||F|||||  
F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|EOS% ^EOS%||1.3|0.5 - 7.0|||F||20220330114559||||H2500ID<CR>  
OBR|44|2203300002|2203300002|IMM% ^IMM% ^HALIA||20220330113916|||||F|||||  
F|||||SUPPORT<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||20220330094130||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|NM|IMM% ^IMM%||0.0|0.0 - 0.5|||F||20220330114559||||H2500ID<CR>  
OBR|45|2203300002|2203300002|^BAS% ^HALIA||20220330113916|||||F<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||||W00002|||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|ED|BAS% ^BAS%||  
^IM ^PNG ^Base64 ^iVBORw0KGgoAAAANSUHEUgAAAgAAADsCAIAAACe1MmJAAAV

## HL7 Format

### Laboratory Testing Workflow Examples



YkIEQVR42u2d6VOU157HfZcXVt7Nzb1x7k2UjU1f9gqqbmTSpVc  
+May6mpeOPC0oCoiYJxgYbu1gYBG2RRFBEQFFGRBhFlaSlo  
+6KcJeB2oxUTHWO5xTJBNMg8pk2LDbl0D/2cfs7nW5/  
q6ss1hPCccc26ePud3zowhQgghZFhm8CsgbBCCGAgghCAGQgghilEQQghilQQqghIIYQ  
gBkIIiYiBEEIiYICEEiYCCGEIAZCCCF6FENCYILu/  
qLpY09uwSZjrPI6f8WtUjaZisrKzutx4SuCPeLRyAP/tU1wJdNQumwGoshNX2Ho7Fz  
+jhUVvIv//6x8jqt/xa12Ja8/datW7r8e4RhxSq/eATy4F9dA3zZLQXQ2LSNIq/G2X+dP/  
+fv2KYU3kN5V1LfQ9xACIATFMjs1brL/  
88osuxZCakZIXZKfvlQZADlhckRFGwchB3UphorK6vjKDPOeYgDEwBrD5Ni4KUqvexWu  
fXdjdcQG+h5iAMsG/YyvhObM2vXRyqtPjhoo1GvYhgYeLY0MIS  
+Jw7+1TXAI01CezFYLBaehxujMUBHG5wRA4BfoP8ZA2IQJ0Eh4TxIAMSg/  
RoDYhAnYSu/  
4hEDIAZmDJPDZDLpWAYxvFuvRCgdPGQAXsMYwCcxms47FUHCgKGPpp4yAGJgxoA  
YXqWlvTPKHMdTBkAMrDEghle59+BhUOgqnlAYqCOYcIPo71b32IYHBxcEmCg11HHA  
NQxaLzG4EflnY2dTn2LQcnsQvplQ  
+A5XPNJWEGF4nNHw13Q8xAGLgrCTE8DrxidsO2k/  
QAxEDIAZmDBOivbtP92KorKm1Jm2nBylGQAysMUyITudl3Yvh7r37QWFsTEIMgBiYMS  
CGP/LixYslgWxMQgyAGFhjmOhHSb1Wq3V17+EOpCQAiGGcPKSms1n7IRi6emw2m  
+7fShbdhopvm+ie1DEADQyarTH4EU0d5zMyMnQvhtz8/TuyC3jAFQ  
+a7bG4EfuNBVIZ2frXgWxenoJN8bwuAEQg2ZrDH5E9akz+/  
fv170Ynj17zsEYAliBgCOEqHTUFRcXD0mQuQs4GAMAMbDGMAGOn6w+fvY4DGJYE/  
HNIvPNPHEAxMCMYRyOlP7HA4ZxFSesyWvpsnDoAYWGMYh8JDxQ0NDTKI4dGjn4  
NCV/  
LEARCDNjMGP9qsXVBY1NHRMSRHAoLD6H7UMQB1DNqsMfhReefe3H0XLlyQRAxrlI  
moPIZqHzWaMbgR60/c/eeK1euSCKG4pJSIhkQAYAgbdYY/Kj1b0/L+P777yURw/  
0HDwNZZKAMgBiMYyNLtnlp59+Gplms+d9TidEDIAWGMYi/itCT//  
LM8YuCaT8QAiIEZwzjEmswDAwPylFrPhEDIAbWGMZ46aolZnCNZ  
+IARADdQyI4Y1wzSd1DEAdgzZrDH5EtDFGKjFwzScAlc/  
azBj8CKNkYhjmK8AxDBqdmRm8RikFQPXfAlgBmYMiOGNcM0nAGJgjWEszGazbGLgm  
k8AxMCMATG8Ea75BEAMGtQx  
+BEGg2Fivnw2dwGPHgAxUMfAJOF1AtmYRB0DUMfg4zUGfynvbOi4IKcYgijlcpnoPLZxz  
MGf2n9jZ1OOcUQwIF6iAEQg4/XGBCD4llcv  
+I4bQN9EJAYmDG4Enz2R45xWBN2FZYUkGfRAyAGFhJ8KTTeVIOMezem7crr5A  
+iRgAMTBjQAYvUlXanpiaSZ9EDIAWGPwpl2712q1SiiGhqZwo2UrfRlXAGKgjMHE4nN  
Ht81mk1AMV65991XkRvokdQxAHYpV1hj8hYa2czt37pRQDPcePAowrKABAFD57LsZg7  
9wqqE1JydHqjH82v90aWAlDQAAMfhujcFOOmolywsIFAMAwPPEAMAYmDGMAoVwB  
V2u31yiAGAMTAGsMolB0/  
qUROMSwPDqMBACCG10IJTecxvNwoVmyvq6uTUwwcsAqAGFhJGIX8A0UtlS1yisGwY  
hUNAAAxUMfgSV7+gfPnz8sphvDva6rPtNMbqWMA6hh8tMbgL+Wdu/

fsvXTpkpxiiDZtPlrhoFts+QxUPvtoXuAvrT8tY8eNGzfkFENy6o7cwqN0S8QAI MFHawz  
+0voTk2x37tyRUwwFhYdSd+XQLREDIAZmDG8QF5/  
w6NEjOcVQ6fh2c0IK3RIxAGJgjeENTGbL06dP5RRD5/  
nudVEmuiViAMTAjOENNm6KGp1P9z6v5CVa  
+iWiAEQA2sMiOFVHj1+wqkYiAEQA3UMnkQbjdKKgQNWqWMA6hh8usbglXiNMdKKYX  
BwEDEAUPnsuxkDYvCLzJm/  
kDYAgBh8tMaAGPwiHJcEgBiYMYzYrmoyySyGNZHfVNa10AwAEANrDK8xm80yiyE1IzO  
vyE4zAEAMzBgQw6tUVFbHJ2fQDAAQA2sMiOFVrn13Y3XEBpoBAGKjGjuGPx9DeLbkYB  
gaesWVOVogajoEjMV7T2OmUXAxKEAOVz0DIM2JADG8kiJufEQMgBt+sMSAGf0nYyq/  
omYgBEAMzhle0d/chhliLIQs+EQMgBI/UMfhF6+90XIZ  
+D5KLoeBAUcaefXROxAClgrNDs5rPOhMSEiXQ0t7Z5Q5js6JGAAsMbwkrrmjTUVm  
nFcO/Bw6DQVXROxAClgtqG38XQ1JaVISW5GAYHB5cEGOic1DEAdQyclfSSmvrG/  
Pz8lekzd/7nNAYADfkm2uzgrCRBOF5Ve+TIEcRQYi/  
bmsKJSQAafXTRet5VZ8pZSUJw9NiJ8vJyxKDKH0sDaA8AWhVUiSIGZgwKh0vKampqsA  
KTBgDEwBrDKwoPFTc0NGAFV/4+Zz5NAsD3NKkuhhkzZhw  
+fNjKx5vmDG8jdz8Ax0dHSjBIRCU  
+QTQgo4L19QXw4IFC48eDBZMBDGoJCId3T04MSXIncsKmitpFWAeBjOp2X1RdDb2+  
vwWAYUQwnT5786KOP3nnnnXnz5t27d486Bg8yMnddu3YNJbiStc3K7+iXkodA/  
j6o6SIOrv/4z/9Ss47BpYH169e7PysfLobw8PDbt28r76uqqkJDQ325xuAX5Z3J29N+  
+OEHI0BK3ZIGUxwfmFL5DL6m+ITDv/zpPTUm10a6Ov//TTTwcGBsb4KGNwRfM  
+nDH4ReuP25rgMZGSOXfv3TeCwQ3YgCf1RV1v7pvT  
+rLwYI5eXIMTEy4hh5syZvlxj8lvWb7JsefLkCurbAzEABpSUh7yvb  
+8Py1iULJw4cLe3t7hYhi+svrBBx8M/wdTUtnp/VHRxt9+  
+w0luDN73kJ6KWIAx4uh9Phf3p81XWL48ccfP/nkk+Fi  
+PDDD12Lq11dXYsWLWKNwYONm6KQATtWEQNoXID1/qx/nS4xvJwHpKs4v/Luu  
+8Wfxd//PHHM2fO/Oyzz27evMka2IYO7EWa/  
HxGjoqXQN8XFA1669/4z4GgT5KQgbDc6Sk1JaRRUela4Avyc7d99e/fcB9DKJgNMYgg  
+G5+s/  
rqyM20FHpguDTGUNewf9+uZz7GBCD0Hn27DkbkwB8LYZ9BbvzDjo4XVUQYmNNyMA  
jy4PDaBgAvmRXVnbewRIH9zElgtlsxgRsTALQlh2ZuwpLKpgxiMLwM6YIG5MANMGWnH  
K0wuHgpYR+Lb5HDMGNIYBal41Lr7ccYYZgxAOdFxADGxMATB  
+sdNkrj7TzhqDED2OhEDG5MARNge6XpDHYp2tJ7vRQYjZu78Uu5M6nRe/u/  
Z8+jD1DGA6qxfv8HVJPS/xiB+eacy0im/BDQwMus3RZdV13v8ukrsZbPncik0lc+gPI9/  
vcbVJPQ/YxC/  
9Td1diclJaGbkTnd0BRrTfSqaERERHzitoP2E3RjxAdqonQuUcTAWUmnmzsyMzPRwMg  
MDg5+GRDsMV2orq6urKm1Jm2nGyMGUJflyEhmDKJQU9dQUFCABkaNISzc/  
Yvqunhi8eLFQ79f8RYURvkbYgD9ioE1htKKrvdjgNGze7sHNfhlFvtXcnJyX19fcoXX7x4s  
SSQDUuIAVQuqGLGIBDF9mMOhwMHjJrbd+6ErPxa  
+S2daWxJS0tzf0oJjYejBhA3X3zAolR2aW5K0//0BRS0sLDnhbgoKCIgnlmsXLht9+unb

## HL7 Format

### Laboratory Testing Workflow Examples



```
dhopvm
+jMiAH0KQbqGHZn5zqdTgQwqeTm79+RXUBnpo4BVBRDVFSULHUM4pOWsfP69euM
9ZPKhZ7eyl0xdGYAtWjv7JObzbJUPotPFELi3bt3GesnfQ7MAFC90IYgMXBWUozJ3N/
fz1g/
2cwZ7cAMANCDGJgxbNwUxSjvRb6OWHfiVDP9GUAVzi28KpAYWGNADN6luKTUlR6b/
gygCs5L/7RarcwYRCEq2sgo70XuP3gYGLqS/gygCh3dvcnJyawxilLRGMMo711mz/uc/
gygznbVtnOZmZmzyBjE36yNGLxOylpV9GfQGEAVTje37dq1i/sYRIFberyOdWsS529T
+QyqUFvfoliBs5IQg9/nZJWD87cRA6hCRZUjLy+P+xiE4FTLOctgdTh/
GzGAWpSUHS8qKmlGLIMr5Jljb63D+NmlAtTh81F5aWsp9DlhBD+H8bcQAqrC/
sKiysplZgxh7hy9cUtT1
+O51OH8bMYAQZOfk1dXVscYgBE2d3UIJSYzvXofztxEDqEL6jsyWIhbuYxCCusa2zMxM
xnevW/
nb1DGAKmxPTbt69Sr3MQhBTV1DQUEB47vX4fxtAFXYmpj08OFDUSqfU1LTZX4YpRVV
drud8X0q4fxtgKijMluePn0qihgkPyup2H7M4XAwuE8lnL8NMHVcxzyLlgbJT1fNP1DU0tL
C4D6VcP42wNSJNhoFEoPkawxZe3OdTieD+1TC
+dsAU8d1miczBjG2iGXsvH79OoP7FLM8OJSODaAfMUj
+xhCfkHj371GdpYZABADdQzunQCb+/v7GdlZZqCOATQei0ymx48fcx
+DQDsBCMsMVD6DtpjN5lu3bnFWkhBsiopmWFclXPOJGGAqGAwGgcQg
+VIJ3OupVklLDV903EQN4x+n2bmYMiEGHiU/cxjWfiAG8w3X
+v0BikHyNgsY1EplTS3XfCIG8I727j5mDGJ9rseYrkq45hMxgNd0Oi8rf0dnjUGgz/
UY01UJ13wiBvB+xtDVY7PZJJoxiLxZu
+XcRcSgYrjmkoG8I66prbs7GyJ6hhEpcq2rNy4ujgFdrXDNJ4B3VNXWHZlyZlizkktGVGPrm
j17GNDVcTtd8AnhHybGK6urqlc5KEoGyiuqysjlGdLXCNZ8A3IF46EhrayszBiHIP1DU3NmZ
gK5WuOYTWDuysnP6+voQgxck79x17do1BnQVExwcbLFYDAaD
+5WLPwHGZxtq2q1btXCDEFi2xD1+/JjRfFrDidwA47J5i/XJkyesMQgBR6v6IjzIDTAu0X
+czSPLjEHkzduS1bJtIYTualjgHFxHdrgFqXCEBsby8Dtg8zhRG4qn2FMTCaT0IM4K0kIK
Hv2TUrSztTmuj8iAHGHos4K0I7GjoulAaf5R9LA+j8iAFGPw
+j9bxwYpB2xtB6vhcxMGIADKA5rSsYxBKDtGsMZ1rPpaSkMGT7LH
+fM48hADHASJrP9jBjEIWTjvrCwkLGayYNiAG0pdN5mTUGUThcUuZwOBivWWIADKC5
GCwWi1wzBmE3a+/Jyevu7mawZtJAHQNoS0f3RZvNjlcgdg7AkBUu+ffs2g7WPs3z5co/
DIFyvcxcsdK3CBQUFuf8ngBS7kn6/
pcfVQWSpfBa5Bv358+eM1EJNJs45L8XHzO3AKIw39IzxFljmhMVzXkYwq1AFJeU1tXV
sSABUuG+pYcZg/bExHAehncThiVLlvb397MgAVKx/+Ah1y09AolBzjWG
+rau4OBgxmLRsmzZsrEXJFh7AP3hVqWHGYpGnO254jq1irCvCUBbbMnbXbf0sMag9Ta
AhtacnBzGWT9dimAoAT3hVqVniPsYtCV//8GOjg4GWSYN1DGA5rhvDBOojkHOIzEsmzf/
+uuvDLj+Gn0cu0TIM7i3zrsaNmclaiUit09mYmBiGv/
+NeUu8DsZTxAaUjAKKQcKzkprP9nB3m1/
n2nfXV61djhAH7g3wjBj0PRkkq6LCQkJKDK9+ndlzFyAG0MNGmNbzBoNBODFIeB9D
+YmqiooKxlaWoBEDaE5bV5/7xjBmDFqyPS3j5s2bjk3sW0UMoDldvVetVitrDCLsAeCUJ
```

```
CYNiAHE4HTjvn37pJxsiLZZ23WwM6OqPPTWO52XxdzeSh0DKBQdKamvr3e1Z
+5j0Az3wc5EkkmD8mdiYk0cpAFikpq+48aNG+4mzVIJ2jD8YGei
+5UGZboQERExxEaIGwRQ0zs8lthOCTJG8wWi/tgZ6L7SYPy/
7pOuf0PRCT2Ng3TvNkxjD
+R8OuU5fzFftOX6tq7AwEAGU32vNLhbjsFgUN5MakECwJcoI5J7r6pYyHb2jcH9d72vl
9adONWsjmy6e1lg0P2kwd1ydLaLCfRHe3efoGIQc8bg/
mhYSXFJqS19tyrf9nCxnQUGfa80DG85k12QAPC9GNxFDKwxTGK6oOT
+g4eBoSvVOZPEbGaBQd+ThrdNF5g0gDU1L8uYpBrxjDGZu3hnhwV7vA7/
aFhJcHDwqH9sUhc9nm7v5IBVHcfVjxaztv
+mEdbGqMhjWylal0vSh2DPLha0ciWM7yIYUiqOoYxyjvH/
svd1HekjDw7jwUGMtmGNLKVqjXtoPJZHlytaGTLScvYObYlgbOSxvks2IsPI6lglKo3pLe1
UIXWKhCDPNMFdyvyaDkeRQyclaTCdGGyf32jgoFMtiFN6wYnxCDVdGHUIhMT88bFMA
KJYW1k5PKAQENIyOLFi5XXMd4vWbo0ICBAeV0RvnKM9x7/bFBw8KJF/6O8enz9yy
+XKE5S5isWi0V5Vd6npKa73w//+rjvly1b7vqe4/xsS5ZO/
HtmZGZN6meYvvcutP0Zpv6MVH+v/CSqf0+ITXq0GeUrb/vzX3zxhbs9K3/Yi/
4SGrZi4cKfW7+PKn1w2fKAcB/PhPrLmP16lt9T+Une9j0n+N87qZ9z7O/jfkYT/O
+d1M859vfxGOs8/
tnh7SoufntYwJjyqnxFYzEQQgjRTRADIYQQxEAIQQxEEIIQQyEEIIQyAeEMRACCEE
MRBCEEMhBBCEAMhBDEQAgBDEQqghBDIQQQhADIYQQxEAIQQxEEIIQQyEEII
QAYGCCCEAMRBCCEAMRBCCEEMhBBCCGlgBCCGAgghCAGQgghilEQQghilYQ
gBkllYiBEEKI7/L/yNqfkbQj9gAAAAASUVORK5CYII=||N||F|||20220330054605|||
H2500ID<CR>
OBR|46|2203300002|2203300002|^NEU%^HALIA|||20220330113916|||F<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||||W00002|||||
WARD00002^^^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|ED|NEU%^NEU%||
^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUgAAAgAAAEACAIAAABK8lkwAAAOe
0IEQVR42u3dOZlccxmG4b4GHNm6hhw4OgUdWbTHmogx4egC4B3kyaJLEzwRHQkM
RDSa6Ons3Nen4jUYyE8vIVxvFw/XnUpfL1tvL69uVi63X9r9/fefRv9+/5tFfb9/9S
+84vvh08crDkWBzbanTmw2W53t19+
+YE400K0AtBpQ9yh1ploSaDUQAEQJtBoIAHLk2AKOhgMBQI4QgPQHAUCUQKuBACB
HoOFAAJAj0GogAlgSaDUQAOTIyd3IGg4EADlyogC0GggAokSrAQQAOLXLS/
ABnEwgAciRHABsowdkeAoAcyQn01QXgVAIBQJTsXOqhABAA5Ag0HAqAcgRaDQQAUA
aLVAAKAHDmrP8B5BAKAHDIRAE4iEABEifE/AAHg7BwhAIAAIEc0HEAAs/Ly
+nZfjkCrgQAIQI4oAREACABYhAAIAAQAOQKtBgKAHIGGAwFAjmg4uwIEADmi4QACgB
DRcAABQI5oOIAAIEc0HEAAEClaDiAAbWD/
```

### Output Format for Host Connection

Ref: RAA088BEN

## HL7 Format

### Laboratory Testing Workflow Examples



aziAAOQINBxAAEIEGg4gADkCDQcQgBxZmm0WiXPWgACwW4i0DmgCAAgAhwpAww  
EEADmi4QACECLQdGABCBFoOIAA5Mi5TN4PoeFAAJAjBAAQADYNESOCCAAEEAAJwe  
Eh/EADkilbTcCAACBFtBxAAHliGAWgAB  
+TlsX0DzhQQAE4PEQIACADLh4hxPtlfjimpHNgE4AAgABAhouik7UAAECIEoO1AABAiG  
k7DgQCwXYjoCSAAEIDdQQAOWIAAUUCIEIb0BwFg4xDJi3ICAAgAWwIAIYgAQAA4NEQ  
IQPqDACBEoO1AABAi2k7bgQAgRJKKRdsUIJwdIABikBMF4NQAUAUCIaDuAACBETB1A  
AJAg2g4gAAiRDecZaDsQAPYPkT4zGqkAIADpf6gANB9AAL15eX27lkF2Ha0v/  
QECIIAeAjhBGwQAAoAQIQCAACBBtB1AABAi2g4gAAiR98cRLV0gckaAAHB6gsRH  
+U4CcDqAACBEDh0k6nQAUCIaDuAACBBNB9AADgkQY6dlexcAAGAAI4TWOXiXAB  
4NT0Pz0CQAEAAJoK4Dhqw+FQ9+5AAKAAaBpj33PmzHQrWcXx5s+AMgrAtgt/  
cuvl9sJYGxt52nrOBdAADhdAAr6AAHAIMJ5JZR9mQ+AAKR/  
w0BvJACHDxDKQKYp5gz6pt0TnznAggAuwlgrV6BgZf5zGUQAHYLjvkFMEIhx7kAAkBP  
cBy7ls7q1XznAggAuwlgrV6BgZf5zGUQAHYLjvkFMEIhx7kAAkBP  
fo0Pnz5ecRQBBAEAB9V2CAAgAALYNvGX2y3OBRAA1ouMnkWen379eGWzar5zAQQA  
Aghd5t8LQGcsCBEAAW13vB0J/18B1LoAAMHte9F9p51YAAhcgAHZpupt6QWAALB/  
fk01aMexBxAAAARipqUEBAhAW3RM/  
XN9vUf1v1KNAAAABHBoW2Zf55QJIDXQCAAiAADqF/n3g1o3g7HczDBQgAOnf9jK/  
tQAEIkaABLBMfy4BACCAGZNRyKCdIUpw7AEEIP17hH5eZ+  
+0HbaEBBDaqjFxebYNvMa/fWW8NjaQhHMBBIBWMBHkKz4WCcAgACK/yyX  
+Wt1CBsGChDAVgKYKvHDA0DjSz0CFyAAvJNZtS7z  
+yy6QAAACKAosFrUdvqvuhOvBwIACMD  
+WmxNzdQxPGEbD8k1LkAAkCsAGbrpC0f6Hnlom8EABBAfvrPNhzwPnDzijx508QIAC  
AAAhg2QDOmjFNy7V/4DQ0DBQhAgbjO/UF8QKdW/  
wUuQAB21giCSBUDAQAEYHcsHFwPmiiBAAAC0AEwtQBSS0apU8ZoHiCA4/  
bUPCP3U4d11i0QEQAAC7/BwtgoYIPAQAQABrVPkP7ANwLoAAMD4XUuO4/  
LqeAAACwBQhVWsn5KKvytugADW4+X17cqWIVW3gNO/J+D2nXvu2w  
+fPI5x8gMEMNH4n1E9Af2HgY7qZH4qAHckIADMJYCx43wMAwUIYI2rdeN/  
CIAAAAjYwWb9OoFbPAPAMFCAAjYXwNir3f4Tvob8UoELEAABjZ2NVAMBADoBJ5+7  
wzdPxlwLEHEIDxP61kUPJX89wPKUkBBLD55X/  
rjt9ybRzS6A5+ElBd02Pn1L1yJwACAahgmSwof2Rj68p+3owBw0ABAIcAqTt  
+6wb6BglACABfDhn5E/+JBNANK9aBAHbuAjiz7FMyW5gACAAEYADoYgKotVyEOz  
+AAFz+Vy71jBoAeuAMAIAAaCtAOYz6JIXGilAgACKwCbFn  
+FBbxho5IWlKxRtjy57ZGAHELL0myvupLgPb/YYah5d9IG3SjQ/EBT+VgTQJ  
+7v95LzF7UOMCWgAeN/1uoJIIAhcR+5I8gAqXGvD6CrAOKr7atc  
+1fxWgzhVf4dWhRwYv6EDyR+RqM7SiqM/  
0ntVl1r5E8jAVwKtnkEMGe9ngzEPQGMF8D8nb15Gmh6m3WptLUQwKJ9s3xwVMGQAE  
ZwAFav73eYA9ytzJxJFr3fP7ZRgJLbURf4BNC7Q3IPAbQu+0w4qnLmuG/  
36XxwVnWtWdTBOb47imACY+0Va7u+3wfmJgh8Qmgd0LNeU/QbaT/nNG/  
1rBUPcningBm7P5dXQB7zKsgADKQ+ASwgABmjv66j4c8fKLvsZEn7lf9OY7plQ9/n

+eqv64AHFdkIPEJgACGaeDRJ9ad2SD9sZMPTIszjQC6VipWWfE/  
WwBPD6dfv7PFX0AS3+rDWRg1VUCqNwBsN8jX  
+qWfdoJoMo7E8DePrAoMgG0LVnsN8s38pXDjyICIANxTWBTJP85asbvTocQASztA4IPAA  
aAduo01t+LsTIQ9wQw3eXVUSv8OM0wQz  
+tECOAWSLsEAE45TBcBnYIAUwhgBVXAJ2kz1y13Z0WCGD28T9119Rcd63/2901+Yh  
+AgAloIcA1tJAAtgB+2FcEQAAggOPG/6z4qPfUO5g+GTShQggABEAAPRX/  
+QUQL4I2R8goAWx57wlQwJgLqxULPknze3/Ys6sHKAEApwsg41fHFH/  
WXfMnvnkboCKY4AAfHPAxos9PK2MxUd2zCunEoA+ABCA8T87IHpSx/  
7fBvE94ciOeeUSV/oEAAALQ/ftlyzuDcJd1WADIV/EEQAAGgMUESGLdP1UAj/  
bGctfvBAAQQM0zaubBrnfHfp7yzNyyAkv8LgACWuaTaZom3ikM84/  
sMNpbB58vnK5IFBHDkbf6KAnjU7iUBXSKADcRAACCAVTsA8qZKzTalq2Rmb3kEbyAAf  
QAggBMHgJLYNTdQPz3D3d9lxRq4kcxQ8hFbUEAALi5+X17Ur1AaDLLfbw9ffmTeDK  
G+KZKgMCIAAQwAICWHEA6LffWx7leaG/xLU/YAA9ADvN84vn3LZ//  
miYAMFQUIYLERqfIH/  
sRf9aeO5Y8v6RAAQAAGgHYd6Bnz4MbUhr9iXr9o2QcggMtrV3angk/  
edXRqoPcRwHA96AMAAZwogHZr/pS8T8zf3v66dt28qdFdyZB678IABDAifO/  
ngpg7F1CpAAiGzG1L/rvqFP  
+xIAQACHn1qS8pNM8gpM78qOyNTIrisA8wAIAAQwWAC1F11oMaUrvqJVPgInflRP6I8tlQ  
CAAEE4c/xO3IMLlv2KU8O2VLZ7cm7eCf+uO33YPkwHQujC+9/b7DZEvu/+T8L  
+U8H27FCa/vnvK/  
ffteonv5IAEAHQREt1BfDnf9+Gfkaa2Xx5dvcK8WONphWAPgAoAU19mFb/  
IPvCTkxk13qf+9c/EsCjPVBBrIbZRV/E6gQkABND7UwJpXn5FH+yG+F/  
3w91A9uDOFtO7+giggylAASwvABisng44a/OKPpvBVA  
+7LJdCWjaqV4EABDAqgK43QO1qvAtCkF1C0c6jQECSPuUkqL8JHoIC6DufN3ye4L4Aa  
mpcU8AAAEciYB3034DP79ulMfMDW4dx+lElDnnxlelhNWwpzX  
+19/2iMB3P97XpmlvIATE9OLhrg+ABDAAGL49kF5A29m1sAPariN  
+9YCaHeXsJAMCAAEsFUn8Fp1/4AASnZy3rSsA58CRgAgglUFsGLoJ9Wybn9+6ko7dSe  
Lxb+m3Z1BdYUQAaiAAOYAbrQvgP5zgEtG+JTfE+guBvYUwB9//9uU  
+A8iglXvemqJIT6gt3k0/OfL5yvSBARwyY71p6+8/  
ZfbD9om9FMfKhdF3Pqx70t3AhMACGAKAYT/  
dIEBxN2yXAJDP1MFkLcuf7uZwEYB6QMAAYQECPuad1//9VP+KoMI5/o+  
+trh1UDjBVB32YNNHvdiAMACQBhAdz+/YCSB0F1K6Dd/  
vpYAQAAPj6MA3cGezRHxC/FMTMD3RcdLklAgABND9MY67r78P90Ws26/6tuxZQuXE  
+JVPpphUAQABdBRBm407gjFFAgWeBITwPoMwYorXPBHBKJ3CqBh7V/  
TcoAd3vw7yHQdat1PeUBICdBFBoSv/x3++E+/+  
da7YRwP2Xj38acMI84P5DNpdY5V8fAAigvQDUwz2C+y  
+9vQB9nD5SPzymb15RacJIIUAIIcWJaCs6A8rYzTf38L/N9w5XNJlWrdXoFb/  
hKglAGx0B1A1++xRNDnrfoZfKcLdTSN0IEAEMBcAgj/7bRTwLkV9y9x2228PnrNy  
+vblbx3sNlsp2zll4BS+wDixDNH7UEELOFBWCz2WzrlYC2EUDqPwN9JmjqvU

### Output Format for Host Connection

Ref: RAA088BEN

HL7 Format

Laboratory Testing Workflow Examples



```
+AJSARgmg0TD8AaG/qACIBCCA5+VFeRPQAB5Q0IjbyBi3jz1fdbt7iMAGADaCuD
+HaqnZNXjWlqu67cUAAACaFUC6paSfQRgwB8AE8FCQZ8kgC2DkgAAbCuAul
+aAACAAAgAAAIiAAABgbwFs3B6OSwAEQAAAQABnBCUBACAAAAGAAIAAACAAAgAA
AiAAACAAAgAAAIiAAACAAAGAAIAAACAAAgAAIAAAAAARAAAQA4GQB1NpW/
d42m8123vZ/xpGtUtyDYFEAAAAASUVORK5CYII=|||N||F||20220330054605|||
H2500ID<CR>
OBR|47|2203300002|2203300002|^PLT^HALIA||20220330113916|||||||||F<CR>
ORC|SC|2203300002|2203300002|2203300002|A|||||W00002|||||
WARD00002^^^^^^^W00002<CR>
TQ1|||||20220330113916||R<CR>
OBX|1|ED|PLT^PLT|
^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUGAAAgAAADsCAIAAAcE1MmJAAQ
y0IEQVR42u3dW2xV55nGcW6iXKDCzoxygaJlo5mMKo3mamY00xjb4PpEOBV82LZDKb
UMQSGDgkhidsY1U4fDONQgcGhCmKkY1sZGsApleUQU2xjGxLjCKIV1AuiNDdOkYgU
WelNZ6FNNrvmB6G977XW79FfS7tuRNrlvb6H93u+911z/
kxERJShOW4BERExBiliYgxERMQYiilMRAREWMgliLGQEREjGliBgDERExBiliYgxERBR
FY2jf9dMDh34BAMgh+w/8/IXmZHANPn/44Yc5NoZdHXs++M0IACCH/O8vf/W3f/
dEcA0+594Ytm3f4VcCAIzhplpaWvxKAIAXqBgAgDHIGAAgR7zSvidY99P/8e3/6/3nf/33f/
m37771zjEVAwDeJu4jJwqKF2YAq23D6sAqtu3aW7dqTeonR070PbfpeAqYwCA6FO
+aMnPt6VaQzf+cd/OtE39Kszg8GHW/95FQMARJn/2Nz8UstPgg+ZxvD3//CdSR/
yyxhkDAawcR+kEm6YnhfxYA0F0aPIZQtDb+AGDpE+nNi1Zptu/b
+ZGdHOMpIL2OQMcyYkzQ1VRz8n5fb2ne7G0Bsjsf1KikgfSpJxRBHhi/1d/
e0VVUVjI93/+IPvbX1Ve4JgNsiY4hFIZBIFDY0FPf2tgaWkClwCUUDgDw1BhXDLOWdpaq
ESZRWLG/APQxxKtKSCZr6+qKbusKigYAmeh8jqwZLFtWuLvz1btUCZQNADIR2PYS/
c1v5lsbhn19ramE+Z7omgAoGKIjHfH+laufPI+zWAS9fXfYWRtxVMVbiPAGLyPISKcOT
+0fXvj2FjFFwhza6Oza8f6nYzAcagYogCp/
sOtLc3TscVAq5efffpH65xMwHGol8h9JwbGaytLZyYODpNYwhYt37VsZPn3FKAMagYQ
m4Mg+90dDRN3xUCBs7v39RsWw/
QxyBjCDk7dq67cuVQVowhoLyy5P2+IXcViDkqhvCRmngUOk20atXCbLICwEcf7V+/
caM7DDAGGUPi+hWqqgoypx5Ii5bWtW+81eM+A4xBxRCyLrYZcoWAIYmjdxVFIU
+Vu9sAY5Ax5HX/
WrBY333wUXbREQ0wBhVDxPvXpoAxSgBjkDFEUx9N0QBAxRCdA0jV1QVZ6V
+bAmUVC5xeBfQxyBjji
+Cv7TN3AMnpVQBPdD6HplzYsGfJrlzB6VWAMcgyYIat3PL2aTNZWV5fs7nzVLwVgDCq
GnHF2eChYkXPrCpMIXKpyUfHh46f9dgdGwBhywMjydbWhrwyhDx8e5E3WLeADAGx
pADunp2fvBBW74ZQ8obSsuLnVYCGIOMYbZjvtzw1VeH89AYnFYCGIOKIQdcHotLJlry
0xWcVgLOMehjmG1yNQPjQU8rVSwqFDYAkUTfkhfkgagaGsAFavhiDjCGTcyODNTXzc
```

```
UDQ9gAQMwQf5x
+8+DBjaFwhRStW5uEDQBjkDHMFcf7R5qbE9euHQ6RMQTFTWm5WXsAY1AxzFhTWz
JZGyJXsKEEMAYZw8yyd9+m0dF9oTMGp1cBxqBimBGL/
UnEoVhdIXUhtKy5YwM7AH6GGQM2SQfZqlOf8pefX1xMIm7aHGp1AEIFzqf87FcyPmrF
6QOAGMwK0m5IHUAoGLISy6MHquqKoiSKzjGCjAGGcO0XKGPqWJ8vDtixmBDCWAM
U1f7ro7Y/iYGLg5WV8+PpCvYUAIYg4xhSoRtAMYUNpS
+V1bqkQMYg4zhvjg1MBS6ARhT4OUfP5v6qgHQxyBjuAchHYDxoPz2d28/89wmTx2g81
nFcG/
2dT4f0gEYD8qKIQs0RQOMQcZw78NINTXz4+AK6abosgqnVwHGGoGKI3xVp1cBxiBjml
orrFjx3bi5gtOrAGNQMDyei2N9K1c+GU9X0A4NMAYZw2TOnB/
avr1xbKwznq5gQwlgDCqYzUO9De3hnnV0hvKPX88u3unrbl3y9yWgnQxxDTjGH4Un
+wDIZXF0xMHGUMwU0lBkVqmmxwLa0o8UwCOp/
jVTEMXBysrS2M2FTtLBjYZlv77rR98gmAMcQgY4j6QKTpU1tfFbhCTc38wD7TPgFAXRB
NYjlQafpFw0vNP0gXVYFPeEqBWBtDtDOGmAxEmj6ZfeCKBkDFEOWKIT4DkblkLA0uAj
CGysXNdXZfV3nglgDF4g5vYWTccoI9BxiB2NI4JOPksYxA7z9LLQct1NgAxM4aoZgxi52z
R/tMXxj/U7VkfVAXiZ9zgj3989+kfrvGsAjEyhmmhDGLnrLJu/
apjJ895XAEVg9zGNxg8v39Tc4vHFZAXiJ1xkxUrirW8ASoGsTMmt7zpiAb0MYSP4Uv9iUS
hRXxG7aG
+vnjJ0pLDx097mlFodj5HrGLo7mznz6oVZYHy8O1G3mDcA0TSGKGUMQbmwYcMSq/
aseUNpebFhSoCKQbkAw5SAqBtDZDKGcyODiYSmNsOUAMagYkgbw
+A7HR1NVuocDFMqW+ipBmQM
+ciOneuuXDlkpc7Jm0G98Q1QMeTjPIJDQ7E1Oid4TTSgJ8E
+EhQNQHq7n6PxBjf7SDmnr44maxdtLjUAVYg9MYQgYzh3Mjg5s0rLM0OsAKMQcZgH
8kBVoaXyBjSl4XhAgvFokLTMoDcG8OcOXO6urom/
WtSh0hWDM4jmaQEMiY7GsPixYuvXr36oMYQ9ozBPILeekOJxjcg58Zw
+fLINWvW3NYY3nvvccff/
zhx9etGjRI19+GaWKwT6SM6yAPoY7GkNw3bRp09mzZ281hrVr137xxRfB5+PHjzc2Nk
YmY3AeKRRnWKurS7zeB8hB53PKBiYmJhYsWPDNN9/cZSvp0UcfjUzFYB/
J298AxnAPYwh05MiRLVU23MUY5s6dG5mMwT4SewAYw72NldDspUsvX76caQyffPJJ
+r+dN29eNCoG55HCaw9IFQs0SAOzagyff/55UVFRpjE89thjn376afD5448/
Xr58eTQyBvtIGqQBxnc/xnB9AIJ7e/
onjzzySE9PzxNPPDF37tyysrLPPvssGhWdfSQN0gBjMCvJeaRINUgVW14obwDy0RhCWj
HYR4pM3IBaUJRALyPwT4S9MEBM9/5HKuK4ezwUEPDAktqZCivdEgJyCdj2LP3tdDd
teGPe7dtW209dUgJYAwqhhu8/saWCxf2WE
+jROvWJoeUgHwxhtBIDKcGhI54ocpKGr1DSt8zjRVQMUyNkbH+F19kDFJogDHoY/
iWY8d/
1tX1omU0ktTWV1kUABXDA7P1v1aPj3dbQ6NaNOzufDW4puZ1r15dkr5WPiVpvYA
+Bhnd7RueHVSNNjU183t7W2/9+cs/
fjbVFArofFYxaHjGdX77u7efeW6T5QOMQcag4Rk3eWrxgsPHT1tBwBhUDAbn4Qbj492Ju
```

## HL7 Format

### Laboratory Testing Workflow Examples



```
sW8AYxBxmAfCX/hDcuWmckKxqBisI+EW2ayVi4qVjqAMcQ9Y/  
jNsPNlMlytdOrcie6eNiYBxhDTimFw5D2D8zDJG6qqCoLqQfYAfQwxzRj2v5400A93MYn  
S8mKzu6HzOUYVw8n+kc2bqy1/  
MLsbjEHGYHAeHoCW1rVmd4MxxKVi6DU4D2Z3gzEwhkxat64xOA9md4MxMAA8zAVI  
n+/  
SAccGEPEMwYNz5hCB1xqcPeSpSWOsYlxRLBi0PAM45Wgj0Efg8F5MF4JOp9VDPaR  
MMObS7IHMIalAZa2kWD6HhiDiuEmZ4eHnEeC1AGMQcaQ0fA8erK1tcFahux7Q6lwmay  
tfKrcegtGELKKoatn5wcfFniMHPdCLs7Xw2uZRULDOADYwhHxpB8ueGrrw5bvzBz1NT  
M7+1tNYAPjCEcFcPwpf66uilrFwzggz4GGcMNggl/  
+KucBQuzNoBv2XLtDtD5nMcVQ1AubNiwxGqF2T/PWl4pb4BZSxmZMSgX4IU/  
YAwqhptcHOurri6wQiFxtG5tkjcg1saQbxnDmfND27c3jo11Wp4gbwBjUDFc53Tfgfb2Rm  
sTjFcCY8iZMezZ+1r  
+3l6Bi4M1NfODv69ZIZBXR3aori7hEFAx5ljTbx48uNFIBAUE9DHIGK5zamCouTlx7ZpWZ5j  
PCp3PKoZvR+Ylk7WWHPjPCsagj+EG  
+zqfHx3dZ92Bt8KBMagYbsTOJiPBthIYg4x7AzbSoCkQeyMil72Z/  
XqkoDgg7OtkDGlnQFnW6FIEDsDQgjoYxA7A1MOIU6dO9Hd01Zfx2zHCTqfxc7AdW  
+oqi47Rj59l7T8KX  
+TOflvPIPxc7jEHsDAT2kHr19N3nNd2pyLgweqy8stCGFWOITsUgdgamE2sHrtDUVPH  
73x9Mn4xKX5UajCGsGYPYGZjyzNdg9a  
+rKxof777nP7ZkaYmSgjGEo2i4NzIYfGs95IBGPMgYvjwGwXc6Opo8tMDseENJ2UJrLm  
O4L7Xv6sjV//kdO9dduXLIewvMDt09bW3tuy27+hjyN2M4NzK4efMKzyowm  
+jH1vmc1xmDfSQgt  
+eaMpsnHGFiDHmRMdhHAvKqecJ8J8aQ44rhZPmhp59e6OEE8vMgrE4I72PIQcYwNNq  
/ZUuNhxAIS8NEpkmk9qBKK0qs5iqGbHJ24F0BAxBsk2hoKA4+O+MkY8gyv+h+9de/  
fsWTBoSasooF7/cNWdNVDNnhlW3r//  
CHtz1XQKj56KP96zdudzHy073TS6cLoscxpTtIKL1L/  
umhk5voYRI580eA8lAq0tK59462e9Bp960mnO810uu0hqLR/  
3Mk50saTmPjUplb0jITfFsWI4O+xIEhARJiaOBot+5ho9naakqCqTtH+kTSVzZOzd/  
y2RGRoYx4zh7IDWNgB3a6ql+dDAOFYM27Y/  
o7UNwlyx3RsShSFylhdxnB2eMilJAC5GgRya0gemEe  
+bT3FrmKwvjwQgT1oxbrv1IA8nnWkXMez4byOSAOTp1tOkk063pt  
+zU1vEq2K4ONZXW1voKwggvP6xbFnhTBQTMel1jOHn  
+aPv2xrGxTt8tANHYkrp1hFS6yMisMzKvZZWL61c3pq8lZZWpySlv/OfWv/  
rrv8mXzudZe4Pb6b4D7e2NvliA4pBb3l1b18bAGloWlgbGELtZSQMXB2tq5k9MHPU1AhB  
nzErK4PSbBw9u9J0AwBi8j  
+E6pwaGmpsT164d9p0AwBhUDNcZGT2ZTJqaBwAhMYZyBj2d2T/  
OrrPfwlAVAw3zm8IENoXAOBuxhCXPoaLY33V1QVZnJglAFE1h1h0PmtnA4CwGsMMZQ  
za2QBAxfAX7+dbufJJ7WwAEEpyHrGELhCU1PFpPe7AgBiWjEMX  
+rPfGsrACDuGUN3T5tjSACgYri5iRSUC37xAKCPQbQAANMMyhqp1Pqfam0ULABARY5ha
```

```
xpD5rqL7fD0FACAcxvDcxo1Pr/rBmh/  
9qKamJrje5XNdf2qVauCa9PadYIE3bPPPhtUG4GvBNeWlpZJn3fvfe22P5/m5/ZdHVn/  
M4Nriuz+mcHn4H9t1v9M99a9nel76z7k5Du2s33XM  
+vXJ5PJ1E9ybAxERBQZMQYilmIMRETEGllilIDEQ3a  
+6urrmzZs3Z871xyF1JWIMRLHWQw899PXX994JBgDMQa3gOL+DHwrkDEGIj  
+fKsZMAyIzWARYyBiDESMgYgxEDEGlsZAxBiIGANRNNowhU24OMQYilmiMREREjGlibg  
DERExBiliYgxERMQYililMRAREWMgllLGQEREjGlibgDERExBiliYgxERMQYililMRARUej  
1/zbUrg7h9dl+AAAAAEIFTkSuQmCC||N||F|||20220330092736|||H2500ID<CR>  
OBR|48|2203300002|2203300002|^RBC^HALIA|||20220330113916|F<CR>  
ORC|SC|2203300002|2203300002|2203300002|A||||W00002|||||  
WARD00002^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|ED|RBC^RBC||  
^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUgAAAggAAADsCAIAAACe1MmJAAAS  
I0IEQVR42u2dbVBU5xmG/WN0ZoXFilB8rEQERDgKSGzWKIikilGxA  
+QiK0SFTUaSxABSzWkX0iQMdUEiBGRwBoTVBSNiroaEtHETq0zcZz80DH+0TjTfJTS/  
EiPHoYg4ALL7tn3nHPdc83Occu0lFeey+d9znu2x2+EEEJli/  
TgR0AIQXqEEIQQyEEIQAyGEEMRACCEEMRBCCEEMhbBBCEAMhhBDEQAghBDEQ  
QghBDIQQRQohoItW/cdKAcR7N334dqP0IV8O  
+zPctrZczUv0yLXR07dXnijK3Z61g7ra41aJvz5887WQxFxSWfX7wCVvj4sxOBwSHSsq8jf5  
NWVK5b07XM/yPd/  
Q0wy1X4e00OGfPtywpqbK1B8zhfDFu2bmMZ1F4sJCu86mJoaQUZ6Z2UsJG4ATEAY  
uha8vPzWQZVF4u2vUlrN8SGhZ28cJl1RAyAGOGydfEs6kv3z3E1PM0Kmt  
+EBi5PX8Y6lgZADMwYooVUBZqR38kv2BE2OjxORNim7cXif4vKiPBqPw8rVpDJHx1W  
+pFZP8u6Z3/5tBmzrC+0REFFhScs/  
lgZADE0pLCrWuRhavTM0dLj56KnKl6dGjBojv3Pk1VVmeukVwHbhTeedenQChlNQ0wx  
Y1/Uz7JON5nUamXbLrS0xBOjp7R6X9i1BsTAjEFRMUgmeHI2siQD  
+Z2ZSXO2Fb8r8VbuXwX/  
5jvZLshlX7y5cJdul92oaclzHzvW4q2wgBADMwYrhw7fXHFmqyYaQnyHzdtL5auw18ct7/  
iM5G/baldmOdq6KQVZJKnxOi5F2y50Kuzctblb2r3ywAQA+cYmhgyNFS+GBsx/tSFy3s  
+KB8/cZJm2gUdNg1Pq/jyQrccOeAGQAx0DO2wcevOidFT5OvQ50d  
+8Pgnfz909PmRo0RuFzo5XWf7NglPYuh5UJb9wAMwZdzxikf0JOiY3/  
sLJKfuevBYWSEkaMGINQWKKldkFXTUPbVqDtQiMGQAx0DNrhrkrqpYvTBd1OGgAQA  
zMgfYghl6PYs6/NYtDb7UkAilGOWf6ldm/  
7tuCA24N9bBbDo0nDH8azrKpYa0AMzBgERajTsJfMVVn9XLtjBYmiCePfp1DJygg  
+1oAY6BgoFkrsI8k8HB60IDWNIUUMgBiYMWihWHR/  
H0lm2aSoY2cusbilARADHYO6i8W5mrMLjH26bwWJuogxmTnslilGQAzMGFRlC6XHch
```

1N9pFDI9G0E8e9QLEAlIBjkF9xcKyNrv7A4ZmXo  
+NO2WpZ30RAyCG1inZvYdlUEuxKE9MOG3ysJcYdsTHIR38hPVFDIAY6BhUfG97wYJQ7  
wN97SWGL2Ki167fyPpyjgEQAZMGFZPT42kvK0j8PCY0deFifqoAilGOQa1cPHw0zc1gRz  
FlxEdO4gcLgBg4x6BWvszJ3eLhZl8xrl  
+OYj8dADHQMaIVvRFjr7nZV8xfDtudMaqTH62AliBGYP6OFNbl+XZz75WkEmd+TI/  
XgDEQMegPq68ty+728/OYzcJADEwY9AO1XNnV/  
i4O0IM7CYBIAY6BIXe274hYOD9IF9HIEEiljKKVeYcAyAGZgxdwOmnYS1VJxYYDQ6ywu  
PdpGh2kwRZawA6BopFp7i6YcOG/m6OE8O/  
oI5iNwKxAGJgXqCmYnEwcfm9ntEUv3Js16hYVGDIAy6BhUUyxyAwb  
+GOznUDGwm4QYADEwY1BNsbhw9GS2u5tDrcC9SYgBEAMdg5qKRf3mLQ4dMDQzb  
SLPTUIMgBj4PAY1FluDMxMdPWDgpBtiAMRAx6Cme9tzB5kcPWBgN4lzDIAYmDGoBk  
efYOC5SQClgY5BZVx5e5MyAwZ2kwAQA  
+cY1EFF3FRIBgzsjgEgBjoGdZA30FuZAcPv9ybx3CQAxMCMQeQBw7p+RiWtIPHNKVM  
+O3mOHZ4AHQMrwYChibNRE/  
l2buGHD8CMgZVgWnDEf4cHpKS8xg8fgl6BIRDx3nblBwwyf4qNY635aw/  
MGFgJazjINkzCjxhasntK9N/  
KDrLWAHQMIFaxcMqAQebOuDGli61krQGYMYBYxclpAwYeqIcYADHQMYhbLJw1YND  
5EWjEAlIBGYOgXclpJxg4Ao0YADHQMYhbLJw4YND5A/  
UQAYAGP09B0GLh3AGDnneTEAMgBjoGQe9td+6AQC  
+7SZxjAMTAjEFEnHiCgXuTABADHYOliDBg4OMZADjHwEolhAgDBj6eAYCOgZUQCBE  
GDHzYJwAzBIZCoAGDc08wsJsEgBjoGBgwsJsEgBiYMTBg4N4kAMRAx6CWe9uFgJDo  
cDeJcwYAGJgxdArFTsNaqk6kGfslZQW97SZx8hkQAx2DWMXiWsGW9e5G0cSgq3uTEA  
MgBmYMYhWLY0mzDv2F1AM  
+tINQgyAGOGYxCoWBc8PvRvol6AY9LObhBgAMTBjEKhYnKmtY/  
boK6AVZOlmt2WtAfQihskiyPzBhGJR7pfqKNtrXg9Nu6UpZ61BmDGAMoViy/  
eWFnsKW7HsCM  
+ruzgJ6w1ADMGU07e9vejJ33l7yWsGL6liV67fiNrDcCMAZRjncn7P8EmYcXw85jQ1IWL  
WSYAOgZQiDO1dQvFO9rW+tkYL01kpQCYMYBCXDJX5bkbBRdD0YT7x  
+oZLEA6BhACS7v3pPfx3QxPBwetCA1jCUCYMYASnBqzZrSAf0EF4PEskIRx85cYr0A6B  
jA4Rx89ZXTg1tu13qIsZk5vAvCQBmDOB4do578VaAt/hieHQEevwE1guAjoFzDA6/  
tz07KEDke1X1M4LmHAMgBmYmopyGXefjqQoraH4EzclnQAx0DEIU/  
PHTy9066MWMWWh7BI0YADEwYxCiWFwp2b1e+EMMOhBlwZADHQMQhSLquRkMT  
+fx9oIOMlcaw3AJAExOKpYCPv5PFaOjAjfXLiLQbQphj4PAbnFovTF74S  
+fN5rJA8JYa1BhBaDD169KioqGj1TqsLZgwCFou68kMifz6P3poGxABAe0NiYuLDhw  
+7KgZmDM69t92S+Zbln89jndixEaw1gNBiuHHjxqJf9oVw/  
HjxwcNgTsrV6+EhIQHDx4wYxCH0sgJln8+jz4nDQDaEYP0mpmZabFY2oph6dKI9+7dk6  
5ramrS09PpGMRB8M/  
n0eekAUBTYmhoaJg8eXJy6OVrSRvb29mDIJw8fDRBUaDeq1A0wCgAjFIOXLKSG5urh  
UxGAwGOGZB+Gp74Yb+bqoWw6MzDaNGnbxwmdUEEFcMUmbNmnXjxo2WYrh+/  
Xrzf2oymZgxCMLBxOmfq+Fp29b5JjRwefoyVhNAADHcvXs3KiqqpRj8/f1v3bolXV  
+7di0pKYmOQRdy/H1/  
DPZTuxgk8keHIX5kZkEBxBXDo2NrhYXN77i4uJjN5pCQEIPBEbcXd  
+fOHWYMIInDh6EmVnmBoS8MQU2xYGBtKABO5+UzH4Kx727/  
cWaSBAYPGNpQ4xwCIGRmDM0/  
DlifEa2DAoLENJU4+A2KgY3Bmscgb6K2NAUPLDaWYlyaw1gCqFwMzBqcUCw2cYGiX  
9dFRqq6qiAEQAx2D04rFlzm5WzzctCeGb8eNzliVyVoDMGNADF1mb8TYq895aU8MEjO  
GhnxaU8taA6hYDHweg/LF4kxtZZnP01aQeJ+kG9K2EiVugEXAGJgXuCcYnHlvX3Z/  
Vy1KgbZDSO91oAYADEwY3DOve3Vc2dX  
+LhrWAZqPdbAOQZADMwYnMOGgJHSv6m1LQaekwFAxwCdxVJ1QpM3qrZ7rCF  
+yBD1DqlBmDGAQlx5e5OWnoSh4UE0AB0DKERF3FSNPQmjQzdMo28AYMYAVtDeKz  
DoGwDoGKBBawbNPGqvbvGAGQMwYKBvAKBjAHvf2663AYPq  
+gbOMQBiyMag9GIYHQ4Y1NU3cPIZEMdg6LFQquP2rbdVnfGlsYAJgxlAbNPmrbBi  
ojwjcX7kIMAHQMeheDhh  
+1bQNjvt67Vq9BDACIQb9i0Pajtm2j2s9jeohYs2jEAlgBMSHXLDT/  
qG1tzKIRAYAGZgzKFQs9PGpbA/ewlgZADHQMyt3brpNHbau9b  
+AcAyAGzjEoxCVz1WvcqMozMwDoGOB3MWRkFHv2pfrzzAwAZgzQxLbggNuDfSj99A0  
AdAzQtI  
+Uxf1I9A0AzBiAfST6BgA6BmAfib4BgBkdMS5mrMLjH0o9PQNAHQMNgn0nLZgVx3I  
1VeLW7gHAMgBmYMDj8Na1mbzYBBRxtKnHwGxEDH4PBiUZ6YcFrHH9mmur4BMQBj  
Ymbg8GJRMCL0+0CehKGavGExAGKgY3B4scj28aSsq6hvQAYAGBCDY4vFJXmVn  
+Wpr4BMQBIAQwOFgNH29TWNyAGQAZMGBxLDjaprq+ATEAYqBjcOC97TWiSY19A  
+cYADFwjsGBsl/EMzMA6BiAfSSemQHAjAGeQL35ofncj0TfAEDHAL8z  
+5UPBvSjcCvjhqkvjOWvHCAGZgxCc/aMJcfL/d/BflRtZaiMCN9culu/eIAY6BjE5crbm/  
J4oqqyTA9h2AClgRmDwLw7POQfz2gWDNsAKBj4BwDY2fV36TEOQZADMwYHHAalr  
GzmvSgtJ4DYqBjsHOxYOWsghte9vHetXoNYgDEwIXBCDEwdhaEaj8P28bRIAEQAx2Dn  
YsFY2e1byshBkAMiMGexaK  
+dH8KY2eVj6MRAYAGxGDPYIEZEV7tx8c7q7tvQAYAGJgx2K1YSO3CG8+6UljV3jcgBkA  
MdAx2u7eddkEbfQPnGAAXcl7BPkjtjwxXpgs8aw  
+AjgFoF3jWHGazBmC6oGqSfG0/+wZAxwC0C5x9A2DGAf3n653vzGW6wKNYAegYQO  
ZczdmtJu9/DuKos  
+4eqQTajAHap3ZGfKfN4os20oAdAycY3iEpepEsquhgfKqrW0IzjEAYmDGYpVJZ8vKfC  
W0C5rbVulKMyAGOGbbxbBIWPDtwT5UVY1tK5V+dAgxAGJgxmCLGL7e  
+U4yD1LVZOsQoiwGMAgXgN7FUFhUzDJ0SQzCjKRhr8aYozdW+ob  
+GsPzBhYiS6lgZuRNMztwT4jeJ8zzcuTO1mBGQMr0VkkXF25YjY3l2laDKN6Py09Vvt58P  
wMYMYAHYtBssKsv3uB/ISQDUvhuaJNHoAOgZ46jkGy5LXX3UxYAVt800Q77tez/  
7w5CqjB2DGAO1ArwDoAegY4Akr0CsAegBmDECvANb0kGo0zDX5YAigY9AX9aX7KYp  
C57jSKwANBDBjQAmPITDP1cAn8EDnG4g8d  
+Mf3flr3QSQMEgwb2jOSgB7KcKPAHMGNg7AuigpaCxAdoG1e8dtXtvO

**HL7 Format**

Laboratory Testing Workflow Examples



+jnHIMywmjpdPyBGPQyY5CLb4e/Bq1eHfG70eo7kV5fM1rbO2p1Ghb0c/  
JZwlajS784aAYxiNUxdLX42taMd/hr0O6XdfU7QQyIqav7Wt3UjA1fjJm0P2NoW/  
1tLr4UC2CtdTVxsa+T7PhlM4cNU/5jwK+uXLHizSyBxLDqzTcXpC1ctHjxvHnzpFcr1/  
NTU9PS0qTXJUuXydcPkF0zc3KlniM/P196IRyjevthUUZy5fn5eVp+P8j16w11/  
L1th07k1NSWtW6dq87rJnydYf/PY9raYr0v9v8PUivzp8x/  
EYIIUSYIAZCCCGCiYEQQohmghgIIYQgBkIIIY4T507dwIDA4ODg6ULfpqEEKJwdu/  
e7eLi4uvre  
+zYMfkds9ns5eVINBoPHZ5sW6HurhJS09PLypKS0uXLl3KChFCiMLJzMxsbGy0WCze3t  
7yO5lSbt26dfPmTckWthXq7orBZDI1NDRl35Z0wQoRQohT8uuvv/bq1Uu+9vf3r6+v/  
+677yZPnmXboe6uGHR27NnqghBCiMlpKSIZvXq1fG02mw0GQ2Rk5IMHD2wr1liBEELUn  
dra2oSEBKlpkP84f/78urq6oqKi5qMMSotB7IB+  
+ukntplIIUT53Lx5Mz4+vrGxsfmd5j2l5ouuFuruiEjI6OsrGzv3r0MnwkRPIERkZKbmj5jpeX  
|  
+VxmofPXS3U3RXDvXv3Bj3O3bt3WSFCCFE4PVpEfsdsNru7u0t6qKmpsa1Qc8CNEEIIYi  
CEEIIYCCGEIAZCCCGIgrBCCGIghBCCGAghhCAGQhYzIooKk8k3wbeFDm4IYIBEP2m  
Z8+ev/  
zyS9OvBGilEifAdH778CTp0YRAyH8DhDyhAwQAYH8DhCCGAhBDIQgBkIQAYGIGRDE  
QAhiAQXEGIPMBR9wD0hilEQQghilIQQqHADIYQQXEAIQQXEEIIQYEEEEIQAyGEEMRA  
CCEEMRBCCCEEMhBBCEAMhBDEQAghBDEQQghBDIQQqHADIYQQXEAIUUF  
+T8hhoelDlnLQAAAAABJRU5ErkJggg==||N||F||20220330092736||||H2500ID<CR>  
OBR|49|2203300002|2203300002|^PLT-Ox^HALIA||20220330113916|||||||||||||F<CR>  
ORC|SC|2203300002|2203300002|2203300002|A|||||||W00002|||||||  
WARD00002^^^^^^^^^W00002<CR>  
TQ1|||||20220330113916||R<CR>  
OBX|1|ED|PLT-Ox^PLT-Ox||  
^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUgAAAQAAAAEACAIAAADTED8xAAAS  
U0IEQVR42u2dyXlJRxJE+0v6rP+3OdVJZjryj/  
oy02M0URBRFekeS65Oy4PEBkEQiJcZi0fkjx/6ev769cfP3wv/p89vvn7/9X+//vV92f+663r/  
q4Nv1Osbjr8GfU1BFPsw+zFP/9p8ZgNg/PU//WDc4OJPcg/A9dfH02o+Qltdn2/9kCdp/  
pT9APcr/  
zK44lv0PcnrbxcACwOQQk7pb0nhB6QUfMzregdDAIzcxVP2UfaVfDMItxHfPswBgOOdbA  
JgvEgB0A0AOjw6AAC+gODzBL2alEcKgLkA6HCqxK3W8Cj6v  
+fB3y4ABvhFhhsWfOAny474LdV/kQBYBgDcnrIO+qd/  
tQ39KWqselUzF4CYDqnCM9U2iw147+mud+mEec3awGwc14lzOE09/5mltAAoP8GX/  
SLBMB0/IAJPvuRVlowJSMUjxDsE6YlgH/qwwwKgaCv69nhk58ZdqazbYv/  
2ukxuOgBygbYCIlJxUqFt0+KfAoZg1LtKhCAAxgAT9BxAwBBHwv4kePFJmfy4EoAzC6  
eYeuy9n84ihJUsmhgclka7A+AL69/

q4XEzyJEUkbFMLmbejBcFgDDklh7O0ZEYSTyqSAB8eLwM6FPmCsA1nD9lwo2u5jVdM  
2RXZwFaiLzJgyOBUDX7T9iE82KIQ2AQ1kd9M3qSlejTysARoYEjiaPwWntXgZNSEAvH6z  
Kle70AmCfGJjablnXQtNBoo4O/Fes4vwlgCqL9/k/PhEEspFTMs/mSQW  
+wtzqYZ33JQASjJ51eNh6E5WeB1NDzXAZLCMYdPXJoQmAeXX/olozkjC1UzE  
+1yglwFrCWwGQ37eFhKTByg67uztWdU1jKthAAHRqXGQzKmuZBYlr7cc3QwjQ/  
aj2JHPiEgIqXUt1SzmeEKkuxePj220jNRq4mmQe5ggzUMCoBAAajvHk0ilgSjukt3enLq  
qGpAh/  
lZimgVAMN6u0AAKJfdnfVnnxapzVEA1LIAduVbAOSHemxipOlwZwW1uHuDH33G8w  
+RhVPvs1oiowBQcuLm5gfu6A4wWJUEWFK4hs7JintEOgE6YeOrFrOpetasl1kjt/  
xuqg9LAFQJH5oRszvdiWjxQQvOCOLrAOhTXBMA5UpPe4e  
+zG70dKsF2wnAbrJS3b8AmFdD6/  
jk8LgTVPz7OnfdqCAIAkctryjAFQD5fbcVbYr2ycCanS8yduT4lwD4RnoJgHnPimYyBG9bA  
f1sh0o0cj64axFU/  
6SjQ0gu0BrFL3wcOejrUxKGrBoCGDk0k6Tg4dmtkC8AcoZwgl6BXVittmMqkYq7OshYoT  
46aspTEgCFAFzKqFfQCksD34iAIGofo/rg0mPI99cpAOimb98P4rtp1p4dCQnYRkqQ/  
CvvtSsT1YAOJX3EQcJaUFcaLHSAhaodWn7ngDwqxc2Q  
+wd7FIEIcHwOW9SSnrriUB0G/7T5llwo/q5wAWdaJtzgrDTrmKGB/FnFyKE0/  
BUyKIQZXJYHqCeRdjU+9RfYpAZDfA4m4TMhhspszBA/  
p7ffQpAmBwbZiq3S4dAYMvWJ4z0+UC9QaAVcWs5dAntt1UC7HA5oqbpjABkDXRhM1  
wDyxsgc2+VEXC3QyUXIUQAL2RCA73nPMEwANcXwUNbF3XHWHLtwRa2q1HxrKzJt9  
cLVLqIBAAUSGQe7NfpRRwwYPUHQ  
+j3tuKvpu1Aag4HFN80MQu21WSoZfrwoHSGrAA6PSc  
+OyTSbb8q2aerm9rHzhVRS5QJ9WQkSxaMd0J5q8QxVT/  
0eoCIHkmCtIDuVyH2u0+sXMvRMAqwlQTJKspe50bPZstDCqCVgAFHaEBbtktHE  
+UPpZAbAVAAe6OnjBGNFK9eHh9bUJgEwXKDLcfCcwqGmkRQCAK64FQGHZK4XsY  
d3FLTy3H5iAVArkcgCYln2ANBHSin0Zt0/IABYFKBIS8CKZt2hdtZMmGbVDQRATg  
+Rb7jV3tGwl6pxn6vps7QEQKjINH0A2zaBL5Xmj7TdRX5KAISaLcBZ51fBRSwz7OL4ZBQ  
wTKJmbTi6tAVANOOJd8CwauFt3HpqyMUFXDvZ1zmpjrAEITne9XtaGIAMULmAxOqP/  
fTAYi/OWCe5/Je4ruW+o1qLkOOSnbjV0dYPwAOz/  
O4R3RRUzMuZhyvQ511IlgApnXubaxkq7k71AfD0g81QAb  
+4UgBEs0BgDuRAAQJz9elbvADorYo77TSg5kFQ5k6IU3HdrgDlxyOUKgnzA/  
F5SSsS8NO6RIAzlmf4NCbk00faQAAnXh24ip+dPsBGNvLM6TIBb8C4/  
DTIDg38ilv7YsUIQUANyBr5iuUHaLqu6Ci1hFXIAJsuUCQ2sSR3Tu8+st2QVD  
+DF4WkBW6X/22xxDc2SaN2un/yHwapUFz1G+lqvFMFVBwlKTP  
+fQNEIUWKKQdxy943LvfN3Hur1v/  
4ygFCID8NjFhtKCCAQgRpYjems6qAOBOALCHVekgfB6EbzgKeDglgNqpMnatXsEuruts  
hr94vyV1jAgArgV+9UGfwzXszekPrOw52BQWBWCDclhkBsQJDe  
+5lrmLvBAG5SoN4yMTAM7asNL/WSOj02dG4D00x7IAjniL9WK1fBUxx+QBXY/  
B0T3B8cvYZLjVrhEiimYrx8ZepiCYq6FoFd2zRCmp7MwE1VUjAJyNAYqAi8alsG0ArCBU

## HL7 Format

### Laboratory Testing Workflow Examples



g7EyC8NaFScJ3kiANABoMly5MklWnKiJsD3+eG+qBmNFq  
+tgYCAwHPoIHwDUQS0APEOy4gOhtKjGMVDSQ4mLkpvID+kJNlxJ  
+T9FoXCzzYUaN20/vwDw33ck6y/  
qpAFHLeE1YA3H9ec6Nfuts3KuORc114qOIEOvez74evNk0LnD5HA1LghAv6b4PsPdh0ig  
pQPtHADt+Qlqcl5AJv1AlAjDBQkRA4Bn/6Z+IDev68YQJXgqcdmgf3y7rGkAsA/  
JFRrBgDwSRO3p8H+AGS5ZDLWnklgxxQCX6gmACSHXqNdhU+gJ  
+mWiK5W0m04ho4qscaEQJmZ/aGOCEpnjldupEoPiwvabFI62qpwfBVAOk8UnlpvuMy/  
V1qFIORXqC22NnfBk3LXcWCDdoUBpkbG26J7gtM2T7ObSd  
+94cCoDISCWbQKuA1Adc1/AFItvdt2BUqUH3k0PjJzPWPORINTwXdl  
+bAKi6GUm9YfKa2A+OTH3Dv99jNuhshDjGUERlvt34gk  
+hBPXS7GjeHgBMQoLvZch8cxnw3bXjENID6tEeQfC0pwTYy6F0Zy4Axn7FTNAhk  
+GorTntBljsr7MBoEz/DMofvBv7Im/JngiA+WsCygsIov04AM1xD/  
aBgHu2qgNwPo9uQ6rTPyPteEbB3qfsEgChlFiHgLv/  
HTxOkdDZAYBmg3qaHmXrRQDg8W4QgMzJcBubuzqAS  
+tiIBLRThn5zmcB0E4aaDjcdVgM3mWEalFY9tQT/NGMxrt99yWGlw1iNtdjdlSU/  
wHIY4QD0FvnlwKQCQKYtPrpQ6wtQG2bsUgF8A09yhDJ0iKOGSANxcRC0fAOCdaxQA  
eC2562zQ/cjihb8pjpCjEo9LJC7XNcMCAB08L7PG4UeEnPHMhK+nTC2RUHEE/  
Gy0WH0UmNe3SwHuqYkHAeDofaG8f5GAjLi6WleUUhqhplKrmVw66IKMagAU1LJTrqgx  
hITzV7C7UjFfCQZVLUCjVp4JECN9WRjjOm0QI68VWkQLADc/g/  
eWI0tdl3N0UAQNGVysDudFCdheBxQtVYIC2F/  
lfrVmAxQE1eCJLQTBPFA50CQN5OyRTEyCHvsHsq0hAA9CAgLydNJKm9KkbHGeCJs  
PRDZDCgGrkReJUX07JPQpFoxG5fKis39fsgu/utgQIGVCAK  
+cSAKAUfyv2BsgdYt0exJlMnNCB41FSCNsYAEXG4DxDqm3FTtJtbgxYIKPoCgXBIRnf  
GViSxSNFj  
+CsMRwANqHpe1VHAIB3iGpRbLCdpb6+7dz613T9AB22fPzdFwwOSVxE7OCrXiUeBUc  
A4G4HO00Bwaob8JrXDBNA1BHmBEDJH0R9gAxnr4i+gFwSBHgcF8IbM/FS1qUcL/  
D8XI0AIjPo2DAoW+jhAzN97zaBuZ1gar//qzZrppsnpud7JlaPB0Agbej  
+7n60ykKEqoAyJJ9SxE9tqbrmBac66AXPQPirS0GQNHhGGwD32B8A+6yg7N3cifT  
+CwN+YsWS4PmAoBUOjXD0JZe1n0u1VutOsJ  
+lqNmzpllwjaU4MKHOgchfvKcOxwXPNA3ntaWpXgTALtNh6am3c  
+fpsTBuPa691YAOGfCrWX3+OsHZ24GpZ0LTUw7sSUSBwCZYtkhU  
+R4cvsvYjOeCwHQVQ26qxwoZRT4JL3q4E1B7vFhi7q4UoNCH7BjONxwHuKNEEMAGI  
WWACCm0/TRDrEOFVKiCuY095C9CABPM8AM  
+VBqYFEzYikFYDmvWJPh0IIIM1g5UrF2Z+JTyrpLAPDtLxUAPy/  
gKrWeHSfx2bTBkHcGcZsAGDMOqDnwzJ2UDM6cYjm5gKmx/  
d5ywXyzFqybZc6GdhYFk9GIUqHvRPW8fXj25fSoODAj6DKkmpCZ6ePzLD3r2Huf3/  
pBOAASJcqcAWHZtrKoSs+0Nzfd/zTAUDkLkXDky94bGBQcrzoFNKcK7d  
+VmlQvGMQGG+IbOrNXi382oimSEnmrjoAkQxFpm  
+7701x52rA8tZOFp9r7kcD4Ejz214EVZzC54m7i3d7pC8NQ//vn/9fvb/

uZADsNCg7N9M9MrY6XVMxyad6g399sAAY5gv5dDjN9PzkWvnl/  
JITAchVrVRMz0xhYDnPfkJznygGyGopKgKgmXXx/S0dLv  
+Rue8PQJ3y2TY7tnEW6bVftBtrCXNHYobjYgDcR49fCpQlwlQW/2leM+/  
uAsAvhsMdd7BYu5bFI7t7RU5m5F+6JQDB+wbdAKw1Y2dR910AVAm/  
fN2Dq9RiTzb3zSvBWZftUGKERLC+iBmZu6QQH79cQ6CQ  
+2gjF0dXACBzVwxAuEZlw2F6FSJ3p5e59z4B9IAag1XbKFquBMvgxYAOaeEr1Ali19ivSZ  
wd3CBUIQVszWRvFeaugkkBcA  
+AMQTPgMTIF8AuludbjzO4eqUQlg8TljRtxEAhwKQZbgVqX059AKga4Zn  
+CuRxQuAJbd/  
5W0EwHEtkRGcZPECYENbNwDQNN8oAlfUwrYxetUKlgRgBsySIFrbfMCYEKXaCAAT11  
/  
8m0EwAFvhLx5AfAEwE5XjMS9ebkWAmaxAJ4cG99OLwDkAs3u33besaolAGYBYLIMZR  
aQAvvclHjpEFYACADlbbQEglxeSwDgdj/  
cSZDLscEbNS8AFTu9ABAAaQBUDNP67riUhQnXYQCkDDlzb/avjwx+wK8/vqWtsH/  
UlcCEXKBIOj/  
Fw3ktbN1+YPinuD0A0jGmiAFyPftmZXfspsyQNgTasesfm7gUAOCcolS91nEu0CpG39yY  
Z9i5dXosA8Aoo//16z+/l+N3PUXMFcan2Wx7AvDN3D9t8dUuX/9jNgDch4DDKGXh  
+wBgeDi5AOA/  
aLzyxLLAcGsWRdMB0MG3MXj4+o5xZWfcaOYxOwFwtBju09xfT5V3JCYEQFa7CQA9rf9  
p73+1e9A7aka61Vu1ABAAtHNvOD+3ANxa+dd33v  
+jev8eIT4SbD0AsM0X2aF9YfthD0nAKP8fgQBcC/  
asBhAHznAw7ANzuYKpOjtTwA3YJdMAV0vdwkhxv3kNig+vAR1TsA0Ez4fPvO53/  
caqG7FXr7A6BjrR8AFdb/tMd/M/  
qngPjpcKAakMJUqxMATw5MEwBjTSJrkyEKAMviDYfet26fcyoAhMRWAPi2f7t0FbF7Ow  
g2uilXck6E0BoAGJkZ3Nxxh+eWlgoAdW1IAKwNgOGusHYPIiDeX6RsSCsBANz/  
ud3yb62cPRMQcbXaLLW6ApDi1IMngJbW1ABUMNChm0xL6xEOk8I+ve  
+E0CEal0EoNThsvtebgFQDKC1MABBDemXEEifolYaAPae2mf7FwBaOwNg1LbYo0Auk  
FY/AOqynIpotfYHILF3TEurBIBS6zckPVLdaO0JgN0JkPqXqBYDW7AAU7fcKgrXSAMjN//  
Qp3AoArUk6KNss+9H0tIaDECHtI8A0JoaANYLUnpUaxgATbukcyy  
+ubaqEmjNCMCtdTaj3qciQNY5oCBYayQAt8bqrnY5dn0BoNUVAMR2O9yMJAC0CvsB4  
rbbDHyzXHwBoJWcBUr5+mrXeu3bun2AvvQ1/Ot/  
NOynA3LG9agAAAAASUVORK5CYII=|||N|||F|||20220330092736|||H2500ID<CR>  
<FS>  
<CR>

## Output Format for Host Connection

Ref: RAA088BEN

### 2.3.5. LIS to Yumizen P8000 acknowledgment

```
<VT>
MSH|^~\&|^|^|YP8K|^|20170626173511||ACK|83990442966255|P|2.5|||||<CR>
MSA|AA|83990442966255<CR>
<FS>
<CR>
```

### 2.3.6. Automatic and Manual Check-in Examples

#### Automatic check-in

```
MSH|^~\&|^|^|YP8K||LIS||20220516143539||OUL^R22^OUL_R22|274846414128113|P|
2.5<CR>
PID|1||00000001^^^PI||PATIENT 1^TEST||19740705000000|M<CR>
SPM|1|2022051300000101^2022051300000101||BLOOD||||MAIN
LAB|||||||20220513135352<CR>
SAC|||||||20220516143527|CHECK_IN^Automatic CheckIn^YH2500 - YH2500||748265|
7<CR>
OBR|1|1151||2022051300000101<CR>
ORC|SC|1151||22051300000101||||20220513135352<CR>
```

#### Manual check-in

```
MSH|^~\&|^|^|YP8K||LIS||20220613103837||OUL^R22^OUL_R22|1641198913631860|P|
2.5<CR>
PID|1||00000011^^^PI||PATIENT 11^TEST||19831211000000|M<CR>
SPM|1|2022060700000011^2022060700000011||BLOOD||||MAIN
LAB|||||||20220607094433<CR>
SAC|||||||20220613103833|CHECK_IN^Manual CheckIn^SUPPORT<CR>
OBR|1|1381||2022060700000011<CR>
```

```
ORC|SC|1381||22060700000011||||20220607094433|||||||1<CR>
```

## 2.4. Patient Detailed Structure

MSH - Message Header Segment	Segment ID .SEQ	Field name	Example
MSH ^~\& ^ ^ YP8K  LIS  20210413112501  ORU^R01 1873659553185571 Q 2.5<CR>	MSH .0	Segment ID	MSH
	MSH .1	Field Separator	
	MSH .2	Encoding characters	^~\&
	MSH .3	Sending Application	YP8K
	MSH .5	Receiving Application	LIS
	MSH .7	Date/Time of message	20220408121010
	MSH .10	Message Control ID	1317701668710
	MSH .11	Processing ID	P

PID - Patient Identification Segment	Segment ID .SEQ	Field name	Example
PID   0002^^^LIS^PI   PATIENT 2^TEST^^^^^^ ^^^^^^  19830421 F   Main Street^^Springfield^N Y^65466^USA^ATC1       ABC123         Y<CR>	PID .0	Segment ID	PID
	PID .3	Patient Identifier List	0002^^^LIS^PI
	PID .5	Patient Name	PATIENT 2^TEST^^^^^^
	PID .7	Date/Time of Birth	19830421
	PID .8	Administrative Sex	F
	PID .11	Patient Address	Main Street^^Springfield^N Y^65466^USA^ATC1
	PID .18	Patient Account Number	ABC123
	PID .31	Identity Unknown Indicator	Y

PV1 - Patient Visit Segment	Segment ID .SEQ	Field name	Example
SPM 1 202204081014   BLOOD   MAIN LAB<CR>	PV1 .0	Segment ID	PV1
	PV1 .2	Patient Class	N
	PV1 .3	Assigned Patient Location	WARD00004^^
	PV1 .7	Attending Doctor	ATD^
	PV1 .17	Admitting Doctor	ADD^
	PV1 .19	Visit Number	ABC123
	PV1 .44	Admit Date/Time	20220408115540
	PV1 .45	Discharge Date/Time	20220408115540

SPM - Specimen Segment	Segment ID .SEQ	Field name	Example
SPM 1 202204081014   BLOOD   MAIN LAB<CR>	SPM .0	Segment ID	SPM
	SPM .1	Set ID - SPM	1
	SPM .2	Sample ID	202204081014
	SPM .4	Specimen Type	BLOOD
	SPM .8	Specimen Source Site	MAIN LAB

The next 4 segments (OBR, ORC, TQ1, OBX) is a group linked to the comment of the sample order.

Data linked to the result of parameters are given in the following groups.

Data type is differentiated by OBX Value Type (OBX .2), ST stands for String Data.

OBR - Observation Request Segment	Segment ID .SEQ	Field name	Example
OBR 1 2204081014  2204081014  ORDER_COMMENT^R equestcomment^HALI A   20220408115540               F<CR>	OBR .0	Segment ID	OBR
	OBR .1	Set ID - OBR	1
	OBR .2	Placer Order Number	2204081014
	OBR .3	Filler Order Number	2204081014
	OBR .4	Universal Service Identifier	ORDER_COMMENT^R equestcomment^HALI A
	OBR .7	Observation Date/Time	20220408115540
	OBR .25	Results Status	F

ORC - Common Order Segment	Segment ID .SEQ	Field name	Example
ORC SC 2204081014 2204081014 2204081014 A     20220408121004     W00004     WARD00004^^^^^^^^^^^W00004<CR>	ORC .0	Segment ID	ORC
	ORC .1	Order Control	SC
	ORC .2	Placer Order Number	2204081014
	ORC .3	Filler Order Number	2204081014
	ORC .4	Placer Group Number	2204081014
	ORC .5	Order Status	A
	ORC .9	Date/Time of Transaction	20220408121004
	ORC .13	Enterer's Location (Ward name)	W00004
ORC .21	Ordering Facility Name (Ward ID^^^^^^^^^^Wardname)	WARD00004^^^^^^^^^^^W00004	

TQ1 - Time/Quantity Segment	Segment ID .SEQ	Field name	Example
TQ1       20220408115540  R<CR>	TQ1 .0	Segment ID	TQ1
	TQ1 .7	Start date/time	20220408115540
	TQ1 .9	Priority	R

OBX - Observation Result Segment	Segment ID .SEQ	Field name	Example
OBX 1 ST ORDER_COMMENT^C ommentaire de la demande  Woman: >21     F  20220408121004<CR>	OBX .0	Segment ID	OBX
	OBX .1	Set ID - OBX	1
	OBX .2	Value Type	ST
	OBX .3	Observation Identifier	ORDER_COMMENT^C ommentaire de la demande
	OBX .5	Observation Value	Woman: >21
	OBX .11	Observation Result Status	F
	OBX .14	Date/Time of the observation	20220408121004

Next 4 segments (OBR, ORC, TQ1, OBX) is a group linked to the result of one parameter - 'n' groups (= "n" parameters) are sent all containing OBR, ORC, TQ1, OBX segments.

OBR set ID (OBR.1) is incremented at each parameter.

Data type is differentiated by OBX Value Type (OBX .2) - NM stands for Numeric Results.

OBR - Observation Request Segment	Segment ID .SEQ	Field name	Example
OBR 2 2204081014 2204081014 MCV^MCV^HALIA     20220408115540       F     ruleResult<CR>	OBR .0	Segment ID	OBR
	OBR .1	Set ID - OBR	2
	OBR .2	Placer Order Number	2204081014
	OBR .3	Filler Order Number	2204081014
	OBR .4	Universal Service Identifier	MCV^MCV^HALIA
	OBR .7	Observation Date/Time	20220408115540
	OBR .25	Results Status	F
	OBR .32	Principal Results Interpreter	ruleResult

ORC - Common Order Segment	Segment ID .SEQ	Field name	Example
ORC SC 2204081014 2204081014 2204081014 A     20220408121004     W00004     WARD00004^^^^^^^^^^^W00004<CR>	ORC .0	Segment ID	ORC
	ORC .1	Order Control	SC
	ORC .2	Placer Order Number	2204081014
	ORC .3	Filler Order Number	2204081014
	ORC .4	Placer Group Number	2204081014
	ORC .5	Order Status	A
	ORC .9	Date/Time of Transaction	20220408121004
	ORC .13	Enterer's Location (Ward name)	W00004
ORC .21	Ordering Facility Name (Ward ID^^^^^^^^^^Wardname)	WARD00004^^^^^^^^^^^W00004	

TQ1 - Time/Quantity Segment	Segment ID .SEQ	Field name	Example
TQ1       20220408115540  R<CR>	TQ1 .0	Segment ID	TQ1
	TQ1 .7	Start date/time	20220408115540
	TQ1 .9	Priority	R

OBX - Observation Result Segment	Segment ID .SEQ	Field name	Example
OBX 1 NM MCV^MCV  101.8 fL 75.0 - 97.0 H    F  20220408120943    YH1500 - YH2500<CR>	OBX .0	Segment ID	OBX
	OBX .1	Set ID - OBX	1
	OBX .2	Value Type	NM
	OBX .3	Observation Identifier	MCV^MCV
	OBX .5	Observation Value	101.8
	OBX .6	Units	fL
	OBX .7	References range	75.0 - 97.0
	OBX .8	Abnormal Flags	H
	OBX .11	Observation Result Status	F
	OBX .14	Date/Time of the observation	20220408120943
OBX .18	Equipment Instance Identifier	YH1500 - YH2500	

Next 4 segments (OBR, ORC, TQ1, OBX) is a group linked to the result of 1 kind of graphic (LMNE,RET, etc) - 'n' groups (=n graphics) are sent all containing OBR, ORC, TQ1, OBX segments.

OBR set ID (OBR.2) is incremented at each graph.

Data type is differentiated by OBX Value Type (OBX .2) - ED stands for Encapsulated Data

OBR - Observation Request Segment	Segment ID .SEQ	Field name	Example
OBR 52 2204081014 2204081014 ^RET%^HALIA     20220408115540       F<CR>	OBR .0	Segment ID	OBR
	OBR .1	Set ID - OBR	52
	OBR .2	Placer Order Number	2204081014
	OBR .3	Filler Order Number	2204081014
	OBR .4	Universal Service Identifier	^RET%^HALIA
	OBR .7	Observation Date/Time	20220408115540
	OBR .25	Results Status	F

ORC - Common Order Segment	Segment ID .SEQ	Field name	Example
ORC SC 2204081014 2204081014 2204081014 A     W00004     WARD00004^^^^^^^W00004<CR>	ORC .0	Segment ID	ORC
	ORC .1	Order Control	SC
	ORC .2	Placer Order Number	2204081014
	ORC .3	Filler Order Number	2204081014
	ORC .4	Placer Group Number	2204081014
	ORC .5	Order Status	A
	ORC .13	Enterer's Location (Ward name)	W00004
ORC .21	Ordering Facility Name (Ward ID^^^^^^^Wardname)	WARD00004^^^^^^^W00004	

TQ1 - Time/Quantity Segment	Segment ID .SEQ	Field name	Example
TQ1     20220408115540  R<CR>	TQ1 .0	Segment ID	TQ1
	TQ1 .7	Start date/time	20220408115540
	TQ1 .9	Priority	R

OBX - Observation Result Segment	Segment ID .SEQ	Field name	Example
OBX 1 ED RET%^RET%  ^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUgAAAQAAAEACIAAADTED8xAAAHbEIEQVR42u3dO3LcRhRAUa3EsW0tyWtQpCptxYtw5NShdzSRKdFFjQaN1w+NbnwP6gYUSY3m0wdoFgb06VPT8nh87xrLI6/...  N  F  20220408120949   YH1500 - YH2500<CR> <sup>a</sup>	OBX .0	Segment ID	OBX
	OBX .1	Set ID - OBX	1
	OBX .2	Value Type	ED
	OBX .3	Observation Identifier	RET%^RET%
	OBX .5	Observation Value	^IM^PNG^Base64^iVBORw0KGgoAAAANSUHEUgAAAQAAAEACIAAADTED8xAAAHbEIEQVR42u3dO3LcRhRAUa3EsW0tyWtQpCptxYtw5NShdzSRKdFFjQaN1w+NbnwP6gYUSY3m0wdoFgb06VPT8nh87xrLI6/...
	OBX .8	Abnormal Flags	N
	OBX .11	Observation Result Status	F
OBX .14	Date/Time of the Observation	20220408120949	
OBX .18	Equipment Instance Identifier	YH1500 - YH2500	

<sup>a</sup>: This example is not exhaustive as there are too many characters to display.

## 2.5. Minimal Lower Layer Protocol

This section describes a minimal HL7 lower level protocol to be used in a pure network environment. It is an adaptation of the hybrid lower layer protocol. It is assumed that this HL7 protocol is used only in a network environment. Most of the details of error detection and correction are handled by the lower levels of any reasonable network protocol and do not require any supplementation.

The goal of this lower level protocol (LLP) is to provide an interface between HL7 and the network that uses minimal overhead while remaining compatible with other LLPs.

Other types of links, such as RS-232 to a communication server, require another protocol to guarantee their integrity. This version of the lower LLP differs significantly from other lower level protocols in that it has only a single byte to signal the start of a message and two bytes to signal the end of a message. There is no other lower level header or trailer information. There are no other characters added to the HL7 message.

Notation conventions:

1. Single ASCII characters are enclosed in single quotes.
2. Special characters or non-printing ASCII characters are enclosed in angle brackets, <>.

Special characters are the LLP Start Block and End Block characters.

Non-printing ASCII characters may be written as their abbreviation, e.g., ESC for the Escape character. They also may be written as their hex value in the form 0xXX where X is a hexadecimal digit.

For example in Standard ASCII, <ESC> is <0x1B>.

### BLOCK FORMAT

HL7 messages are enclosed by special characters to form a block. The format is as follows:

<SB>dddd<EB><CR>

<SB>

Start Block character (1 byte)

ASCII <VT>, i.e., <0x0B>.

This should not be confused with the ASCII characters SOH or STX.

dddd

Data (variable number of bytes)

This is the HL7 data content of the block.

The data can contain any displayable ASCII characters and the carriage return character, <CR>.

<EB>

End Block character (1 byte)

ASCII <FS>, i.e., <0x1C>.

This should not be confused with the ASCII characters ETX or EOT.

<CR>

Carriage Return (1 byte)

The ASCII carriage return character, i.e., <0x0D>.

## 2.6. References

Title	Version	Date	Author
HL7 Messaging Standard	2.5	2013	HL7 Int.
IHE Laboratory (LAB) Technical Framework	5.0	2013	IHE Int.
IHE IT Infrastructure (ITI) Technical Framework	10.1	2013	IHE Int.

## 3. QC to LIS

### 3.1. Overview

The YumizenLis driver can send QC results on the same channel as patient results among which:

- QC sample ID
- QC LIS test code
- QC result value
- QC result date and time
- QC result units
- QC result device/instrument ID

Detailed information is sent to LIS, such as:

- Comments
- Target values
- Deviations
- Violated rules
- Laboratory's technician user ID validator

By default, all those information is sent.

A detailed structure is available, refer to the *QC to LIS - Detailed Structure* chapter.

To deactivate this function, refer to the technical manual.

### Message structure

The following table lists the detailed structure for QC message from Yumizen P8000 to the LIS:

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
OBR	Observation Order	R	[1..1]
OBX .1	Observation Result (first parameter)	R	[1..n]
OBX .2	Observation Result (second parameter)	R	[1..n]
OBX .n	Observation Result (n parameter)	R	[1..n]
FS	File Separator	R	[1..1]
CR	Carriage return	R	[1..1]



- QC can be resent to LIS.  
In case of resending, only controlled parameters are resent to LIS.

### 3.2. LIS to Yumizen P8000 ACK / NAK

The reply expected from the LIS in response to a QC result message is the standard HL7.



- The reply must be sent by the LIS on the same channel used by the Yumizen P8000 to send QC result messages (standard results channels).  
The LIS can accept or refuse the QC result message. This data is saved on the Yumizen P8000 so that the QC result will no longer be sent to the LIS.

### Message structure

The following table lists the ACK structure for a message from the LIS to the Yumizen P8000:

Segment	Meaning	Usage	Card.
MSH	Message Header	R	[1..1]
MSA	Message Acknowledgment	R	[1..1]
[ERR]	Error	R	[0..1]
FS	File Separator	R	[1..1]
CR	Carriage return	R	[1..1]

#### Example of a reply from the LIS:

```
MSH|^~\&|^|^|YP8K|^|20170626173511||ACK|1873659553185571|P|2.5|||||<CR>
MSA|AA|1873659553185571<CR>
<FS>
<CR>
```

### 3.3. MSH - Message Header Segment

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R	[1..1]	Segment ID	<b>MSH</b>
1	1	R	[1..1]	Field Separator	(Pipe)
2	4	R	[1..1]	Encoding characters	^~\& ^: sub field delimiter ~: repeat sub field delimiter \: ESCAPE sequence &: sub filed component delimiter

SEQ	LEN	OPT	Card.	Element name	Example
3	20	R	[1..1]	Sending Application	<b>YP8K</b> NameSpace ID: YP8K
4		R	[1..1]	Sending Facility	NameSpace ID: Empty
5		R	[1..1]	Receiving Application	<b>LIS</b> NameSpace ID: Empty
6	14	R	[1..1]	Receiving Facility	NameSpace ID: Empty
7		R	[1..1]	Date/Time of message	<b>20210413112501</b>
8		X		Security	
9		R	[1..1]	Message Type	<b>ORU^R01</b> Message code ID: ORU Trigger event ID: R01 Message structure ID: ORU_R01
10		R	[1..1]	Message Control ID	<b>1873659553185570</b>
11	1	R		Processing ID	<b>Q</b>
12		R		Version ID	<b>2.5</b>

#### Example of a QC message sent to the LIS:

```
MSH|^~\&|YP8K||LIS||20210413112501||ORU^R01|1873659553185571|Q|2.5
```

### 3.4. OBR - Observation Request Segment

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>OBR</b>
1	1	R	[0..1]	Set ID - OBR	<b>1</b>

SEQ	LEN	OPT	Card.	Element name	Example
2		RE	[0..1]	Placer Order Number	
3	9	R	[1..1]	Filler Order Number	<b>PX416H</b>
4		R	[1..1]	Universal Service Identifier	<b>QC</b>
5		X		Priority	
6		X		Requested Date/Time	
7		X		Observation Date/Time	
8	14	R		Observation End Date/Time	<b>20210413112453</b>
9		X		Collection Volume	
10		X		Collector Identifier	
11		X		Specimen Action Code	
12		X		Danger Code	
13		X		Relevant Clinical Information	
14		X		Specimen Received Date/Time	
15		R	[1..1]	Specimen Source	<b>^^^^^^Q</b>

**Example of a QC message sent to the LIS:**

OBR|1||PX416H|QC|||||20210413112453|||||^^^^^^Q

**OBR .3 field: description of control sample tube as follows:**

1. Two or three letters for the QC control type: PX (ABX Difftrol), RX (ABX Minotrol Retic), BFI (BFCTRL), SX (ESRtrol).
2. A series of numbers for the QC lot number.
3. The last characters give the QC level: L or -1 (Low), N or -2 (Medium), H or -3 (High).

Examples:

OBR .3	Control Type	Lot Number	Level	Description
PX114L	ABX Difftrol	114	Low	ABX Difftrol 114 level Low.
PX114N	ABX Difftrol	114	Normal	ABX Difftrol 114 level Normal

OBR .3	Control Type	Lot Number	Level	Description
PX114H	ABX Difftrol	114	High	ABX Difftrol 114 level High
RX115-1	ABX Minotrol Retic	115	Low	ABX Minotrol Retic 115 level Low
RX115-2	ABX Minotrol Retic	115	Normal	ABX Minotrol Retic 115 level Normal
RX115-3	ABX Minotrol Retic	115	High	ABX Minotrol Retic 115 level High
BFI0095-2	BFCTRL	0095	Normal	BFCTRL 0095 level Normal
BFI0095-3	BFCTRL	0095	High	BFCTRL 0095 level High
SX123N	ESRtrol	123	Normal	ESRtrol 123 level Normal
SX123H	ESRtrol	123	High	ESRtrol 123 level High

### 3.5. OBX - Observation Result Segment

SEQ	LEN	OPT	Card.	Element name	Example
0	3	R		Segment ID	<b>OBX</b>
1	1	R	[1..1]	Set ID - OBX	<b>2</b>
2	2	C	[0..1]	Value Type	<b>NM</b>
3		R	[1..1]	Observation Identifier	<b>HGB</b>
4		X		Observation Sub-ID	
5		C	[0..1]	Observation Value	<b>165.2001</b>
6		C	[0..1]	Units	<b>g/L</b>

SEQ	LEN	OPT	Card.	Element name	Example
7		X		References Range	
8		X		Abnormal Flags	
9		X		Probability	
10		X		Nature of Abnormal Test	
11	1	RE	[1..1]	Observation Result Status	<b>F</b>
12		X		Effective Date of Reference Range	
13		X		User Defined Access Checks	
14	14	RE	[0..1]	Date/Time of the Observation	<b>20210413104907</b>
15		X		Producer's ID	
16		X		Responsible Observer	
17		X		Observation Method	
18		C		Equipment Instance Identifier	<b>Yumizen 1</b>

**Example of a QC message sent to the LIS:**

OBX|2|NM|HGB||165.2001|g/L||||F|||20210413104907|||Yumizen 1<CR>

QC result is sent to the LIS as soon as it is considered as valid or confirmed when the user has taken a decision (accepted, rejected or ignored), see details below.

Three OBX segments follow the OBX segments linked to the n parameters.

They describe the status of the QC results as in Yumizen P8000 software:

- QC\_RESULT\_STATUS
- QC\_RESULT\_IN\_STATISTICS
- LOT\_DESCRIPTION

QC result conditions	QC_RESULT_STATUS	QC_RESULT_IN_STATISTICS
All parameters are within tolerance ranges or comply with Westgard rules	VALID	TRUE
A failed QC has been accepted by the user	ACCEPTED	TRUE
A failed QC has been rejected by the user	REJECTED	TRUE or FALSE depending on user decision
A failed QC has been ignored by the user	IGNORED	FALSE

	Values	Note
LOT_DESCRIPTION	DiffTrol	-
	Minotrol Retic	-
	BFTrol	-
	ESRTrol	-
	UNKNOWN	If the lot is not defined in the Yumizen P8000

**Example of VALID result**

OBX|42|ST|QC\_RESULT\_STATUS||VALID|||||F<CR>

OBX|43|ST|QC\_RESULT\_IN\_STATISTICS||TRUE|||||F<CR>

OBX|44|ST|LOT\_DESCRIPTION||DiffTrol|||||F<CR>

### 3.6. ORU^R01

The following tables illustrate the structure for message ORU^R01 (compliant with HL7) used to send QC results from the Yumizen P8000 to the LIS:

MSH - Message Header Segment: Date, Time / Message Control ID / Processing ID		Example
MSH ^~\& YP8K LIS  20210413112501 ORU^R01  1873659553185571 Q 2.5<CR>	MSH .7 Date/Time of message	<b>20210413112501</b>
	MSH .10 Message Control ID	<b>1873659553185571</b>
	MSH .11 Processing ID	<b>Q</b>

OBR - Observation Request Segment: QC Lot Number / Sample Type / Observ. Date, Time / Specimen Source		Example
OBR 1  PX416H QC    20210413112453    ^^^^^Q<CR>	OBR .3 Filler Order Number: QC lot number	<b>PX416H</b>
	OBR .4 Universal Service Identifier: Result type	<b>QC</b>
	OBR .8 Observation Date/time	<b>20210413112453</b>
	OBR .15 Specimen source	<b>^^^^^Q</b>

(from 1 to n parameter) OBX segments including for each: parameters / value / unit / QC run date and time / Device		Example
OBX 1 NM MPV  10.50001  fL    F  20210413104907    Yumizen 1<CR>	OBX .2 Value Type: Numeric result	<b>NM</b>
	OBX .3 Observation Identifier: QC parameter	<b>MPV</b>
	OBX .5 Observation Value: QC parameter value	<b>10.5</b>
	OBX .6 Unit	<b>fL</b>
	OBX .14 Date/Time of the Observation: Date/ time of device QC run	<b>20210413104907</b>
	OBX 18 Equipment Instance Identifier: Device name that ran QC	<b>Yumizen 1</b>

(from 1 to n parameter) OBX segments including for each: parameters / value / unit / QC run date and time / Device		Example
OBX 2 NM HGB  165.2001 g/ L    F  20210413104907    Yumizen 1<CR>		
OBX 3 NM RDW-CV  12.3  %    F  20210413104907    Yumizen 1<CR>		
OBX 4 NM PLTO  500.00001  1E09/L    F   20210413104907   Yumizen 1<CR>		
OBX 5 NM MCV  89.00001  fL    F  20210413104907    Yumizen 1<CR>		
OBX 6 NM LYM#  1.97001  1E09/L    F   20210413104907   Yumizen 1<CR>		
OBX 7 NM PCT  0.066%    F  20210413104907    Yumizen 1<CR>		
OBX 8 NM BAS#  0.81101  1E09/L    F   20210413104907   Yumizen 1<CR>		
OBX 9 NM EOS#  1.0 1E09/ L    F  20210413104907    Yumizen 1<CR>		
OBX 10 NM IMG%  2.4%    F  20210413104907    Yumizen 1<CR>		
OBX 11 NM PLT  520.00001  1E09/L    F   20210413104907   Yumizen 1<CR>		
OBX 12 NM LIC%  3.3%    F  20210413104907    Yumizen 1<CR>		

(from 1 to n parameter) OBX segments including for each: parameters / value / unit / QC run date and time / Device	Example
OBX 13 NM EOS%  5.35  %    F  20210413104907     Yumizen 1<CR>	
OBX 14 NM NRBC#  0.845  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 15 NM MON%   9.37001 %    F   20210413104907    Yumizen 1<CR>	
OBX 16 NM NEU#  11.04  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 17 NM NEU%  70.9  %    F  20210413104907     Yumizen 1<CR>	
OBX 18 NM MIC  2 %    F   20210413104907    Yumizen 1<CR>	
OBX 19 NM ALY%  2 %    F   20210413104907    Yumizen 1<CR>	
OBX 20 NM MON#  1.54001  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 21 NM BAS%  4.33  %    F  20210413104907     Yumizen 1<CR>	
OBX 22 NM LIC#  0.05 1E09/ L    F  20210413104907     Yumizen 1<CR>	
OBX 23 NM HCT  0.4725 L/ L    F  20210413104907     Yumizen 1<CR>	

(from 1 to n parameter) OBX segments including for each: parameters / value / unit / QC run date and time / Device	Example
OBX 24 NM ALY#  0.25  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 25 NM NRBC%   5.15001 %    F   20210413104907    Yumizen 1<CR>	
OBX 26 NM IMG#  0.04  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 27 NM RDW-SD  43.8  fL    F  20210413104907     Yumizen 1<CR>	
OBX 28 NM MCHC   352.0001 g/L    F   20210413104907    Yumizen 1<CR>	
OBX 29 NM MAC  2.7 %     F  20210413104907     Yumizen 1<CR>	
OBX 30 NM WBC  17.54  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 31 NM IMM#  0.00  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 32 NM IML%  0.6 %     F  20210413104907     Yumizen 1<CR>	
OBX 33 NM TNC  18.90  1E09/L    F   20210413104907    Yumizen 1<CR>	

(from 1 to n parameter) OBX segments including for each: parameters / value / unit / QC run date and time / Device	Example
OBX 34 NM IML#  0.01  1E09/L    F   20210413104907    Yumizen 1<CR>	
OBX 35 NM RBC  5.285  1E12/L    F   20210413104907    Yumizen 1<CR>	
OBX 36 NM IMM%  0.3 %     F   20210413104907     Yumizen 1<CR>	
OBX 37 NM PDW  17.6 fL     F   20210413104907     Yumizen 1<CR>	
OBX 38 NM MCH  31.25  pg    F   20210413104907     Yumizen 1<CR>	
OBX 39 NM LYM%   11.90001 %    F   20210413104907    Yumizen 1<CR>	

(n+1 ... n+3) OBX segments including for each: QC status / QC in statistics / QC material description	Example	
OBX 40 ST  QC_RESULT_STATUS   VALID    F<CR>	OBX .2 String Data	<b>ST</b>
	OBX .3 Observation Identifier	<b>QC_RESULT_STATUS</b>
	OBX .5 Observation Value	<b>VALID</b>
OBX 41 ST  QC_RESULT_IN_STATISTICS   TRUE    F<CR>	OBX .2 String Data	<b>ST</b>
	OBX .3 Observation Identifier	<b>QC_RESULT_IN_STATISTICS</b>
	OBX .5 Observation Value	<b>TRUE</b>
OBX 42 ST  LOT_DESCRIPTION  DiffTrol     F<CR>	OBX .2 String Data	<b>ST</b>
	OBX .3 Observation Identifier	LOT_DESCRIPTION
	OBX .5 Observation Value	<b>ABX Difftrol</b>

(n+1 ... n+3) OBX segments including for each: QC status / QC in statistics / QC material description	Example
<FS>	
<CR>	

ACK / NACK message from the LIS to the Yumizen P8000	Example	
MSH ^~\& ^ ^ YP8K ^  20210413112501  ACK  1873659553185571 P 2.5     <CR>	MSH.2 Encoding Characters	^~\&
	MSH.5 Receiving Application	YP8K
	MSH.7 Date/Time Of Message	20210413112501
	MSH.9 Message type	ACK
	MSH.10 Message control id	1873659553185571
	MSH.11 Processing id, always P	P
MSA AA  1873659553185571<CR>	MSA.1 Acknowledgment code (AA\AE\CA\CE)	AA
	MSA.2 Message Control id related to the ACK/NAK	1873659553185571
<FS>		
<CR>		

### 3.7. QC to LIS - Detailed Structure

#### Detailed structure

```
MSH|^~\&|^|^|YP8K||LIS||20221219174423||ORU^R01|873930713772998|Q|2.5<CR>
OBR|1||PX914H|QC||||20221219174422|||||^|Q<CR>
NTE|0|P||ACCEPTED^ACCEPTED<CR>
OBX|1|NM|ALY%||2|%||||F|||20221209111602||SUPPORT||YH500 - YH550<CR>
```

OBX|2|NM|RBC||5.485|1E12/L| **5.3; 0.08333| 1:2s,2:2s**||||F||20221209111602|| **SUPPORT**||  
YH500 - YH550<CR>

OBX|3|NM|EOS%||9.35|%|5.9;1.96667|4:1s||||F||20221209111602||SUPPORT||YH500 -  
YH550<CR>

OBX|4|NM|PDW||17.6|fL||||F||20221209111602||SUPPORT||YH500 - YH550<CR>

...

OBX|37|NM|NRBC%||6.15001|%|4.0;0.5|1:2s,1:3s||||F||20221209111602||SUPPORT||  
YH500 - YH550<CR>

OBX|38|ST|QC\_RESULT\_STATUS||ACCEPTED|||||F<CR>

OBX|39|ST|QC\_RESULT\_IN\_STATISTICS||TRUE|||||F<CR>

OBX|40|ST|LOT\_DESCRIPTION||DiffTro|||||F<CR>

<FS>

<CR

### Outline

OBX|1|NM|MCV||80|mg/dl|**expected\_value;deviation|rule\_alarm**||||F||20191028145728||  
**user**||Yumizen 1

## 4. QC Export

### 4.1. QC Results Availability

Specific settings must be done in Yumizen P8000 to setup the QC export functionality, RAS1115 - Configuring the Application procedure.

Contact your HORIBA Medical technical representative if necessary.

To share QC results with external systems (like LIS or QC expert system), the Yumizen P8000 is generating csv (or xml) files with all results send by the instruments.

The files are available through a FTP connection, on the Yumizen P8000 Server.

The name of the user for the FTP server is **qc\_p8000**. The password is **123456789**.

For each export, one or several files are available in CSV or XML format.

As soon as the file has been downloaded, the external system deletes automatically the CSV or XML file.

#### 4.1.1. CSV - File Example (Semicolon Separated) – Used by Default

Qc Sample Id;QC Yumizen P8000 Level;Lot Id;Device Id;Instrument Id;Yumizen P8000 Test Code;Result;Qualitative Result;Reagent;Units;Note;Coded Comment Text;Result Date;Expected Value;Deviation

PX412L;1;PX412;Proto03;Proto03;HGB;69;;;g/L;;;2018-07-04 11:14:37;69.0;1.33333

PX412L;1;PX412;Proto03;Proto03;WBC;2.06;;;1E09/L;;;2018-07-04 11:14:37;2.1;0.13333

PX412L;1;PX412;Proto03;Proto03;NRBC#;0.30;;;1E09/L;;;2018-07-04 11:14:37;0.27;0.04333

PX412L;1;PX412;Proto03;Proto03;NRBC%;14.8;;;%;2018-07-04 11:14:37;13.0;1.33333

PX412L;1;PX412;Proto03;Proto03;IMG#;0.01;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;IMG%;0.4;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;IMM#;0.00;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;ALY#;0.13;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;ALY%;6.3;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;BAS#;0.10;;;1E09/L;;;2018-07-04 11:14:37;0.1;0.03333

PX412L;1;PX412;Proto03;Proto03;MCHC;327;;;g/L;;;2018-07-04 11:14:37;351.0;10.0

PX412L;1;PX412;Proto03;Proto03;IMM%;0.0;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;IML#;0.00;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;IML%;0.1;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;BAS%;4.7;;;%;2018-07-04 11:14:37;4.6;1.53333

PX412L;1;PX412;Proto03;Proto03;EOS#;0.17;;;1E09/L;;;2018-07-04 11:14:37;0.13;0.04333

PX412L;1;PX412;Proto03;Proto03;EOS%;8.1;;;%;2018-07-04 11:14:37;6.2;2.06667

PX412L;1;PX412;Proto03;Proto03;MCV;84.2;;;fL;;;2018-07-04 11:14:37;82.0;1.66667

PX412L;1;PX412;Proto03;Proto03;MIC;4.2;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;MAC;2.2;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;RDW-SD;41.7;;;fL;;;2018-07-04 11:14:37;37.0;1.33333

PX412L;1;PX412;Proto03;Proto03;RDW-CV;13.1;;;%;2018-07-04 11:14:37;14.0;1.33333

PX412L;1;PX412;Proto03;Proto03;NEU#;0.99;;;1E09/L;;;2018-07-04 11:14:37;1.25;0.11667

PX412L;1;PX412;Proto03;Proto03;NEU%;48.5;;;%;2018-07-04 11:14:37;59.4;3.33333

PX412L;1;PX412;Proto03;Proto03;TNC;2.36;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;LIC#;0.01;;;1E09/L;;;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;LIC%;0.5;;;%;2018-07-04 11:14:37;;

PX412L;1;PX412;Proto03;Proto03;PDW;22.8;;;fL;;;2018-07-04 11:14:37;;  
PX412L;1;PX412;Proto03;Proto03;PCT;0.098;;;%;;;2018-07-04 11:14:37;;  
PX412L;1;PX412;Proto03;Proto03;MPV;11.9;;;fL;;;2018-07-04 11:14:37;10.5;0.66667  
PX412L;1;PX412;Proto03;Proto03;PLT;83;;;1E09/L;;;2018-07-04 11:14:37;70.0;6.66667  
PX412L;1;PX412;Proto03;Proto03;MCH;27.5;;;pg;;;2018-07-04 11:14:37;28.8;0.66667  
PX412L;1;PX412;Proto03;Proto03;HCT;0.210;;;L/L;;;2018-07-04 11:14:37;0.197;0.005  
PX412L;1;PX412;Proto03;Proto03;RBC;2.49;;;1E12/L;;;2018-07-04 11:14:37;2.4;0.05333  
PX412L;1;PX412;Proto03;Proto03;MON#;0.18;;;1E09/L;;;2018-07-04 11:14:37;0.14;0.03  
PX412L;1;PX412;Proto03;Proto03;MON%;8.6;;;%;;;2018-07-04 11:14:37;6.5;1.43333  
PX412L;1;PX412;Proto03;Proto03;LYM#;0.61;;;1E09/L;;;2018-07-04 11:14:37;0.49;0.11  
PX412L;1;PX412;Proto03;Proto03;LYM%;29.6;;;%;;;2018-07-04 11:14:37;23.3;4.0

#### 4.1.2. XML File Example

```
<?xml version="1.0" encoding="UTF-8"?>
<!DOCTYPE qcdata PUBLIC "qcdataid" "qcdata.dtd">
<qcdata xmlns="http://www.isotech.com/2007/qcdata" version="0.2"
domain="hematology" created="2018-07-06 12:07:11">
<labinfo>
<instrument manufacturerid="Horiba Medical" instrumentid="Yumizen H2500"
serialnumber="M2XH00185">Proto03</instrument>
</labinfo>
<sampleset controlproductid="" lot="PX412L" lotqcperiodidstart="20180704"
lotqcperiodidend="20180704">
<result time="2018-07-04 11:14:37">
<value units="g/L" refunits="g/L" paramid="HGB" param="HGB" decimal="0">69</
value>
<value units="10^9/L" refunits="10^9/L" paramid="WBC" param="WBC"
decimal="2">2.06</value>
```

```
<value units="10^9/L" refunits="10^9/L" paramid="NRBC#" param="NRBC#"
decimal="2">0.30</value>
<value units="" refunits="" paramid="NRBC%" param="NRBC%"
decimal="1">14.8</value>
<value units="10^9/L" refunits="10^9/L" paramid="IMG#" param="IMG#"
decimal="2">0.01</value>
<value units="" refunits="" paramid="IMG%" param="IMG%" decimal="1">0.4</
value>
<value units="10^9/L" refunits="10^9/L" paramid="IMM#" param="IMM#"
decimal="2">0.00</value>
<value units="10^9/L" refunits="10^9/L" paramid="ALY#" param="ALY#"
decimal="2">0.13</value>
<value units="" refunits="" paramid="ALY%" param="ALY%" decimal="1">6.3</
value>
<value units="10^9/L" refunits="10^9/L" paramid="BAS#" param="BAS#"
decimal="2">0.10</value>
<value units="g/L" refunits="g/L" paramid="MCHC" param="MCHC"
decimal="0">327</value>
<value units="" refunits="" paramid="IMM%" param="IMM%" decimal="1">0.0</
value>
<value units="10^9/L" refunits="10^9/L" paramid="IML#" param="IML#"
decimal="2">0.00</value>
<value units="" refunits="" paramid="IML%" param="IML%" decimal="1">0.1</
value>
<value units="" refunits="" paramid="BAS%" param="BAS%" decimal="1">4.7</
value>
<value units="10^9/L" refunits="10^9/L" paramid="EOS#" param="EOS#"
decimal="2">0.17</value>
<value units="" refunits="" paramid="EOS%" param="EOS%" decimal="1">8.1</
value>
<value units="fL" refunits="fL" paramid="MCV" param="MCV" decimal="1">84.2</
value>
<value units="" refunits="" paramid="MIC" param="MIC" decimal="1">4.2</value>
```

```
<value units="%" refunits="%" paramid="MAC" param="MAC" decimal="1">2.2</value>
<value units="fL" refunits="fL" paramid="RDW-SD" param="RDW-SD" decimal="1">41.7</value>
<value units="%" refunits="%" paramid="RDW-CV" param="RDW-CV" decimal="1">13.1</value>
<value units="10^9/L" refunits="10^9/L" paramid="NEU#" param="NEU#" decimal="2">0.99</value>
<value units="%" refunits="%" paramid="NEU%" param="NEU%" decimal="1">48.5</value>
<value units="10^9/L" refunits="10^9/L" paramid="TNC" param="TNC" decimal="2">2.36</value>
<value units="10^9/L" refunits="10^9/L" paramid="LIC#" param="LIC#" decimal="2">0.01</value>
<value units="%" refunits="%" paramid="LIC%" param="LIC%" decimal="1">0.5</value>
<value units="fL" refunits="fL" paramid="PDW" param="PDW" decimal="1">22.8</value>
<value units="%" refunits="%" paramid="PCT" param="PCT" decimal="3">0.098</value>
<value units="fL" refunits="fL" paramid="MPV" param="MPV" decimal="1">11.9</value>
<value units="10^9/L" refunits="10^9/L" paramid="PLT" param="PLT" decimal="0">83</value>
<value units="pg" refunits="pg" paramid="MCH" param="MCH" decimal="1">27.5</value>
<value units="L/L" refunits="L/L" paramid="HCT" param="HCT" decimal="3">0.210</value>
<value units="10^12/L" refunits="10^12/L" paramid="RBC" param="RBC" decimal="2">2.49</value>
<value units="10^9/L" refunits="10^9/L" paramid="MON#" param="MON#" decimal="2">0.18</value>
```

```
<value units="%" refunits="%" paramid="MON%" param="MON%" decimal="1">8.6</value>
<value units="10^9/L" refunits="10^9/L" paramid="LYM#" param="LYM#" decimal="2">0.61</value>
<value units="%" refunits="%" paramid="LYM%" param="LYM%" decimal="1">29.6</value>
</result>
</sampleset>
<sampleset controlproductid="" lot="PX413L" lotqcperiodidstart="20180704" lotqcperiodidend="20180704">
<result time="2018-07-04 11:33:38">
<value units="%" refunits="%" paramid="IML%" param="IML%" decimal="1">0.2</value>
<value units="10^9/L" refunits="10^9/L" paramid="IMG#" param="IMG#" decimal="2">0.01</value>
<value units="%" refunits="%" paramid="IMG%" param="IMG%" decimal="1">0.4</value>
<value units="%" refunits="%" paramid="PCT" param="PCT" decimal="3">0.100</value>
<value units="fL" refunits="fL" paramid="PDW" param="PDW" decimal="1">21.1</value>
<value units="%" refunits="%" paramid="MAC" param="MAC" decimal="1">2.5</value>
<value units="10^9/L" refunits="10^9/L" paramid="IML#" param="IML#" decimal="2">0.00</value>
<value units="10^9/L" refunits="10^9/L" paramid="WBC" param="WBC" decimal="2">2.02</value>
<value units="%" refunits="%" paramid="BAS%" param="BAS%" decimal="1">4.7</value>
<value units="%" refunits="%" paramid="EOS%" param="EOS%" decimal="1">9.7</value>
<value units="%" refunits="%" paramid="NRBC%" param="NRBC%" decimal="1">16.5</value>
```

```
<value units="fL" refunits="fL" paramid="RDW-SD" param="RDW-SD"
decimal="1">40.5</value>
<value units="L/L" refunits="L/L" paramid="HCT" param="HCT" decimal="3">0.210</
value>
<value units="fL" refunits="fL" paramid="MCV" param="MCV" decimal="1">84.1</
value>
<value units="10^9/L" refunits="10^9/L" paramid="EOS#" param="EOS#"
decimal="2">0.20</value>
<value units="10^9/L" refunits="10^9/L" paramid="BAS#" param="BAS#"
decimal="2">0.09</value>
<value units="10^12/L" refunits="10^12/L" paramid="RBC" param="RBC"
decimal="2">2.49</value>
<value units="10^9/L" refunits="10^9/L" paramid="NRBC#" param="NRBC#"
decimal="2">0.33</value>
<value units="10^9/L" refunits="10^9/L" paramid="NEU#" param="NEU#"
decimal="2">1.01</value>
<value units="g/L" refunits="g/L" paramid="MCHC" param="MCHC"
decimal="0">328</value>
<value units="10^9/L" refunits="10^9/L" paramid="PLT" param="PLT"
decimal="0">78</value>
<value units="g/L" refunits="g/L" paramid="HGB" param="HGB" decimal="0">69</
value>
<value units="fL" refunits="fL" paramid="MPV" param="MPV" decimal="1">12.9</
value>
<value units="%" refunits="%" paramid="RDW-CV" param="RDW-CV"
decimal="1">12.7</value>
<value units="pg" refunits="pg" paramid="MCH" param="MCH" decimal="1">27.6</
value>
<value units="%" refunits="%" paramid="LYM%" param="LYM%" decimal="1">26.5</
value>
<value units="10^9/L" refunits="10^9/L" paramid="TNC" param="TNC"
decimal="2">2.35</value>
```

```
<value units="%" refunits="%" paramid="MIC%" param="MIC%" decimal="1">4.6</
value>
<value units="10^9/L" refunits="10^9/L" paramid="LIC#" param="LIC#"
decimal="2">0.01</value>
<value units="10^9/L" refunits="10^9/L" paramid="IMM#" param="IMM#"
decimal="2">0.00</value>
<value units="%" refunits="%" paramid="MON%" param="MON%" decimal="1">8.5</
value>
<value units="%" refunits="%" paramid="ALY%" param="ALY%" decimal="1">5.1</
value>
<value units="%" refunits="%" paramid="LIC%" param="LIC%" decimal="1">0.6</
value>
<value units="10^9/L" refunits="10^9/L" paramid="ALY#" param="ALY#"
decimal="2">0.10</value>
<value units="10^9/L" refunits="10^9/L" paramid="LYM#" param="LYM#"
decimal="2">0.54</value>
<value units="10^9/L" refunits="10^9/L" paramid="MON#" param="MON#"
decimal="2">0.17</value>
<value units="%" refunits="%" paramid="IMM%" param="IMM%" decimal="1">0.0</
value>
</result>
<result time="2018-07-04 11:48:08">
<value units="10^9/L" refunits="10^9/L" paramid="WBC" param="WBC"
decimal="2">1.93</value>
<value units="%" refunits="%" paramid="NRBC%" param="NRBC%"
decimal="1">22.8</value>
<value units="%" refunits="%" paramid="ALY%" param="ALY%" decimal="1">4.9</
value>
<value units="%" refunits="%" paramid="MON%" param="MON%" decimal="1">7.5</
value>
<value units="10^9/L" refunits="10^9/L" paramid="MON#" param="MON#"
decimal="2">0.14</value>
```

```
<value units="10^9/L" refunits="10^9/L" paramid="LYM#" param="LYM#"
decimal="2">0.65</value>
<value units="10^9/L" refunits="10^9/L" paramid="PLT" param="PLT"
decimal="0">77</value>
<value units="10^9/L" refunits="10^9/L" paramid="BAS#" param="BAS#"
decimal="2">0.09</value>
<value units="10^9/L" refunits="10^9/L" paramid="EOS#" param="EOS#"
decimal="2">0.13</value>
<value units="%" refunits="%" paramid="BAS%" param="BAS%" decimal="1">4.6</
value>
<value units="%" refunits="%" paramid="EOS%" param="EOS%" decimal="1">6.8</
value>
<value units="fL" refunits="fL" paramid="MCV" param="MCV" decimal="1">84.8</
value>
<value units="%" refunits="%" paramid="RDW-CV" param="RDW-CV"
decimal="1">12.3</value>
<value units="10^9/L" refunits="10^9/L" paramid="NEU#" param="NEU#"
decimal="2">0.91</value>
<value units="L/L" refunits="L/L" paramid="HCT" param="HCT" decimal="3">0.207</
value>
<value units="%" refunits="%" paramid="LIC%" param="LIC%" decimal="1">0.3</
value>
<value units="10^9/L" refunits="10^9/L" paramid="IMG#" param="IMG#"
decimal="2">0.00</value>
<value units="10^9/L" refunits="10^9/L" paramid="TNC" param="TNC"
decimal="2">2.37</value>
<value units="%" refunits="%" paramid="IML%" param="IML%" decimal="1">0.3</
value>
<value units="fL" refunits="fL" paramid="PDW" param="PDW" decimal="1">22.1</
value>
<value units="%" refunits="%" paramid="MIC" param="MIC" decimal="1">4.3</value>
<value units="%" refunits="%" paramid="IMG%" param="IMG%" decimal="1">0.0</
value>
```

```
<value units="%" refunits="%" paramid="MAC" param="MAC" decimal="1">2.9</
value>
<value units="%" refunits="%" paramid="PCT" param="PCT" decimal="3">0.100</
value>
<value units="%" refunits="%" paramid="IMM%" param="IMM%" decimal="1">0.0</
value>
<value units="10^9/L" refunits="10^9/L" paramid="IML#" param="IML#"
decimal="2">0.01</value>
<value units="10^9/L" refunits="10^9/L" paramid="LIC#" param="LIC#"
decimal="2">0.01</value>
<value units="10^9/L" refunits="10^9/L" paramid="IMM#" param="IMM#"
decimal="2">0.00</value>
<value units="10^9/L" refunits="10^9/L" paramid="ALY#" param="ALY#"
decimal="2">0.09</value>
<value units="g/L" refunits="g/L" paramid="HGB" param="HGB" decimal="0">68</
value>
<value units="pg" refunits="pg" paramid="MCH" param="MCH" decimal="1">27.9</
value>
<value units="%" refunits="%" paramid="LYM%" param="LYM%" decimal="1">33.5</
value>
<value units="fL" refunits="fL" paramid="RDW-SD" param="RDW-SD"
decimal="1">38.1</value>
<value units="g/L" refunits="g/L" paramid="MCHC" param="MCHC"
decimal="0">329</value>
<value units="fL" refunits="fL" paramid="MPV" param="MPV" decimal="1">13.0</
value>
<value units="10^9/L" refunits="10^9/L" paramid="NRBC#" param="NRBC#"
decimal="2">0.44</value>
<value units="10^12/L" refunits="10^12/L" paramid="RBC" param="RBC"
decimal="2">2.44</value>
</result>
<result time="2018-07-04 14:15:08">
```

<value units="%" refunits="%" paramid="IMM%" param="IMM%" decimal="1">0.0</value>  
<value units="10^9/L" refunits="10^9/L" paramid="IMM#" param="IMM#" decimal="2">0.00</value>  
<value units="10^9/L" refunits="10^9/L" paramid="NRBC#" param="NRBC#" decimal="2">0.34</value>  
<value units="10^9/L" refunits="10^9/L" paramid="NEU#" param="NEU#" decimal="2">0.98</value>  
<value units="g/L" refunits="g/L" paramid="MCHC" param="MCHC" decimal="0">334</value>  
<value units="fL" refunits="fL" paramid="RDW-SD" param="RDW-SD" decimal="1">39.3</value>  
<value units="%" refunits="%" paramid="MIC" param="MIC" decimal="1">5.4</value>  
<value units="%" refunits="%" paramid="LYM%" param="LYM%" decimal="1">30.1</value>  
<value units="%" refunits="%" paramid="MON%" param="MON%" decimal="1">8.2</value>  
<value units="10^9/L" refunits="10^9/L" paramid="ALY#" param="ALY#" decimal="2">0.12</value>  
<value units="%" refunits="%" paramid="LIC%" param="LIC%" decimal="1">0.2</value>  
<value units="%" refunits="%" paramid="NRBC%" param="NRBC%" decimal="1">16.4</value>  
<value units="10^9/L" refunits="10^9/L" paramid="LIC#" param="LIC#" decimal="2">0.00</value>  
<value units="pg" refunits="pg" paramid="MCH" param="MCH" decimal="1">28.3</value>  
<value units="10^9/L" refunits="10^9/L" paramid="LYM#" param="LYM#" decimal="2">0.61</value>  
<value units="10^9/L" refunits="10^9/L" paramid="MON#" param="MON#" decimal="2">0.17</value>  
<value units="10^9/L" refunits="10^9/L" paramid="PLT" param="PLT" decimal="0">76</value>

<value units="%" refunits="%" paramid="MAC" param="MAC" decimal="1">3.3</value>  
<value units="%" refunits="%" paramid="PCT" param="PCT" decimal="3">0.092</value>  
<value units="g/L" refunits="g/L" paramid="HGB" param="HGB" decimal="0">69</value>  
<value units="L/L" refunits="L/L" paramid="HCT" param="HCT" decimal="3">0.205</value>  
<value units="fL" refunits="fL" paramid="MPV" param="MPV" decimal="1">12.2</value>  
<value units="10^9/L" refunits="10^9/L" paramid="IML#" param="IML#" decimal="2">0.00</value>  
<value units="10^9/L" refunits="10^9/L" paramid="IMG#" param="IMG#" decimal="2">0.00</value>  
<value units="%" refunits="%" paramid="RDW-CV" param="RDW-CV" decimal="1">12.6</value>  
<value units="%" refunits="%" paramid="IML%" param="IML%" decimal="1">0.2</value>  
<value units="%" refunits="%" paramid="IMG%" param="IMG%" decimal="1">0.0</value>  
<value units="fL" refunits="fL" paramid="MCV" param="MCV" decimal="1">84.7</value>  
<value units="fL" refunits="fL" paramid="PDW" param="PDW" decimal="1">16.9</value>  
<value units="%" refunits="%" paramid="EOS%" param="EOS%" decimal="1">8.8</value>  
<value units="%" refunits="%" paramid="BAS%" param="BAS%" decimal="1">4.7</value>  
<value units="10^9/L" refunits="10^9/L" paramid="EOS#" param="EOS#" decimal="2">0.18</value>  
<value units="10^9/L" refunits="10^9/L" paramid="BAS#" param="BAS#" decimal="2">0.10</value>

```
<value units="%" refunits="%" paramid="NEU%" param="NEU%" decimal="1">48.0</value>
<value units="%" refunits="%" paramid="ALY%" param="ALY%" decimal="1">5.8</value>
<value units="10^9/L" refunits="10^9/L" paramid="TNC" param="TNC"
decimal="2">2.38</value>
<value units="10^9/L" refunits="10^9/L" paramid="WBC" param="WBC"
decimal="2">2.04</value>
<value units="10^12/L" refunits="10^12/L" paramid="RBC" param="RBC"
decimal="2">2.42</value>
</result>
</sampleset>
</qcdata>
```