

REF 1300141427

REAGENT 2 x 20 mL

IVD  2797

HORIBA ABX SAS
Parc Euromédecine
Rue du Caducée
BP 7290
34184 Montpellier Cedex 4
FRANCE

Yumizen C Calcium

- Yumizen C230
- Yumizen C240

Diagnostic reagent for quantitative *in vitro* determination of Calcium in serum, plasma and urine by colorimetry.

Intended Use

Yumizen C Calcium reagent is intended for the quantitative *in vitro* diagnostic determination of calcium in human serum, plasma and urine based on a colorimetric method.

Clinical laboratories use.

Measurement of calcium is used in the diagnosis and treatment of parathyroid disease, a variety of bone diseases and chronic renal disease and tetany (intermittent muscular contractions or spasms).

Measurement of physiological and pathological variations of Calcium in human serum, plasma and urine is useful for screening or follow-up of these diseases and also in the assessment of electrolyte homeostasis and the acid-base balance of the body.


Clinical Interest (1, 2, 3)

Calcium plays an essential role in many cell functions: intracellularly in muscle contraction and glycogen metabolism, extracellularly, in bone mineralization, in blood coagulation and in transmission of nerve impulses. Calcium is present in plasma in three forms: free, bound to proteins or complexed with anions as phosphate, citrate and bicarbonate. Under physiological conditions, calcium balance is determined by the relationship between calcium intake and calcium absorption and excretion. Urinary excretion is an important determinant of calcium retention in the body. Decreased total calcium levels can be associated with diseases of the bone apparatus (especially osteoporosis), kidney diseases (especially under dialysis), defective intestinal absorption and hypoparathyroidism. Increased total calcium can be measured in hyperparathyroidism, malignant diseases with metastases and sarcoidosis. Calcium measurements also help in monitoring of calcium supplementation mainly in the prevention of osteoporosis.

Method (4, 5, 6, 7)

Many colourimetric methods for determining calcium have been used in the past. Connerty and Briggs described methods using alizarin 3-sulphonate (4) and cresolphthalein complexone (5) whilst Gindler and King have described a method using thymol blue (6). There have been many subsequent modifications to these methods. The method used here is based on the metallochromogen Arsenazo III. Calcium ions (Ca^{2+}) react with Arsenazo III (2,2'-[1,8-Dihydroxy-3,6-disulphonaphthylene-2,7-bisazo]-bisbenzenearsonic acid) at pH 6.75 to form an intense purple coloured chromophore. The absorbance of the Ca-Arsenazo III complex is measured bichromatically at 660/700 nm.

The resulting increase in absorbance of the reaction mixture is directly proportional to the calcium concentration in the sample. Arsenazo III has a high affinity ($K^{\circ} = 1 \times 10^{-7}$) for calcium ions (7) and shows no interference from other cations normally present in serum, plasma or urine.

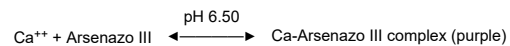

Reagents

Yumizen C Calcium is ready-to-use.

Reagent:

MES pH6.50	100 mmol/L
Arsenazo III	200 µmol/L

Yumizen C Calcium should be used according to this notice. The manufacturer cannot guarantee its performance if used otherwise.

Yumizen C Calcium

Handling

1. Remove the cap of the cassette.
2. If present, remove foam by using a plastic pipette.
3. Place the cassette into the refrigerated reagent compartment.

Calibrator

For calibration, use:

ABX Pentra Multical (A11A01652) (not included)
10 x 3 mL (lyophilisate)

Control

For internal quality control, use:

- **ABX Pentra N MultiControl** (1300054414) (not included)
10 x 5 mL (lyophilisate)
- **ABX Pentra P MultiControl** (1300054415) (not included)
10 x 5 mL (lyophilisate)
- **Yumizen C Urine Level 1 Control** (1300023946) (not included)
6 x 5 mL
- **Yumizen C Urine Level 2 Control** (1300023947) (not included)
6 x 5 mL

Each control should be assayed daily and/or after a calibration.

The frequency of controls and the confidence intervals should correspond to laboratory guidelines and country-specific directives. You should follow federal, state and local guidelines for testing quality control materials. The results must be within the range of the defined confidence limits. Each laboratory should establish a procedure to follow if the results exceed these confidence limits.

Materials Required but not Provided

- Automated clinical chemistry analyzer:
Yumizen C230/C240
- Calibrator: **ABX Pentra Multical** (A11A01652)
- Controls:
ABX Pentra N MultiControl (1300054414)
ABX Pentra P MultiControl (1300054415)
Yumizen C Urine Level 1 Control (1300023946)
Yumizen C Urine Level 2 Control (1300023947)
- Standard laboratory equipment.

Specimen (8, 9)

This device intended testing population is general population.

Specimen types

- Serum.
- Plasma in lithium heparin.
- Urine.

Do not use EDTA plasma: EDTA anticoagulant is unsuitable for analysis because this compound chelates calcium, making it unavailable for reaction with the reagent.

Anticoagulants other than those listed have not been tested by HORIBA and are therefore not recommended for use with this assay.

24 hours urine specimens have to be collected with HCl 6N (9). Non acidified urines which have been refrigerated should be acidified and/or heated at 56°C for 15 minutes to redissolve any precipitate.

Stability (8)

Serum, plasma (8)

- At 20-25°C: 7 days
- At 4-8°C: 3 weeks
- At -20°C: 8 months

Urine (9)

- At 20-25°C: 2 days
- At 4-8°C: 4 days
- At -20°C: 3 weeks

Reference Range (2)

Each laboratory should establish its own reference ranges. The values given here are used as guidelines only.

Serum, plasma

2.15 - 2.57 mmol/L (8.6 - 10.3 mg/dL)

Urine (10)

Women: < 6.24 mmol/24h (< 250 mg/24h)

Men: < 7.49 mmol/24h (< 300 mg/24h)

Clinical sensitivity and specificity, positive predictive value and negative predictive value are not commonly reported for this analyte. This is largely attributed to the fact that this analyte is not sole indicator for the intended purpose

Yumizen C Calcium

and patient treatment decision making. To arrive at a diagnosis and a course of treatment, results from others routine clinical chemistry tests should be used in conjunction with other diagnostic information and the attending health-care professional's evaluation of the patient's condition.

Storage and Stability

Stability before opening:

Stable up to the expiry date on the label if stored at 2-8°C. Store protected from light.

Stability after opening:

Refer to the paragraph "Performance on Yumizen C230/C240".

Waste Management

Please refer to local legal requirements.

General Precautions

- This reagent is for professional *in vitro* diagnostic use only.
For laboratory use.
- For prescription use only.
- This reagent is classified as non-hazardous in compliance with regulation (EC) N°.1272/2008.
- As calcium is an ubiquitous ion, essential precaution must be taken against accidental contamination. Only use disposable materials.
- Observe the standard laboratory precautions for use.
- The reagent cassettes are disposable and should be disposed of in accordance with the local legal requirements.
- Please refer to the SDS associated with the reagent.
- Do not use the product if there is visible evidence of biological, chemical or physical deterioration.
- Do not use the product if the recommended storage conditions, including temperature, are not followed.
- User must be trained by a HORIBA representative before attempting to operate the device.
- It is the user's responsibility to verify that this document is applicable to the reagent used.
- For technical assistance, you can call +33 (0)4 67 14 15 16.
- Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the competent authority of the country in which the user and/or the patient is established.

Performance on Yumizen C230/C240

Serum, plasma

The performance data listed below have been obtained on the Yumizen C230/C240 analyzer.

Number of tests: approximately 2 x 140 tests

On Board Reagent Stability

Once opened, the reagent cassette placed in the refrigerated Yumizen C230/C240 compartment is stable for 70 days.

Sample volume: 2 µL/test

Lowest Detectable Level

The lowest detectable level represents the lowest measurable level of analyte that can be distinguished from zero. It is calculated as the absolute mean plus three standard deviations of 20 replicates of an analyte free sample. The lowest detectable level is estimated at 0.011 mmol/L (0.044 mg/dL).

Limit of Quantitation

The limit of quantitation is determined according to CLSI (NCCLS), EP17-A2 protocol (11) and equals 0.025 mmol/L (0.10 mg/dL).

Accuracy and Precision

Repeatability (*within-run precision*)

Repeatability according to the recommendations found in the CLSI (NCCLS), EP05-A3 protocol (12) with samples tested 20 times:

- 2 controls
- 3 specimens (low / medium / high levels)

	Mean value mmol/L	Mean value mg/dL	CV %
Control specimen 1	2.15	8.63	1.11
Control specimen 2	3.25	13.03	1.16
Specimen 1	1.85	7.40	1.12
Specimen 2	2.47	9.91	1.11
Specimen 3	3.39	13.59	1.13

Yumizen C Calcium

Reproducibility (total precision)

Reproducibility according to the recommendations found in the CLSI (NCCLS), EP05-A3 protocol (12) with samples tested in duplicate for 20 days (2 series per day):

- 2 controls
- 3 specimens (low / medium / high levels)

	Mean value mmol/L	Mean value mg/dL	CV %
Control specimen 1	2.14	8.58	1.8
Control specimen 2	3.23	12.95	1.1
Specimen 1	1.94	7.78	2.3
Specimen 2	2.37	9.50	1.5
Specimen 3	3.51	14.08	1.4

Measuring Range

The assay confirmed a measuring range from 0.025 mmol/L (0.10 mg/dL) to 4.50 mmol/L (18.05 mg/dL). The measuring range is extended up to 18.00 mmol/L (72.18 mg/dL) with the automatic post-dilution. The reagent linearity has been assessed up to 4.50 mmol/L (18.05 mg/dL) according to the recommendations found in the CLSI (NCCLS), EP06-Ed2 protocol (13).

Correlation

Patient samples: Serum

Number of patient samples: 98

Specimens are correlated with a commercial reagent taken as reference according to the recommendations found in the CLSI (NCCLS), EP09c protocol (14).

Values ranged from 1.73 mmol/L (6.94 mg/dL) to 3.55 mmol/L (14.24 mg/dL).

The equation for the allometric line obtained using Passing-Bablok regression procedure (15) is:

$$Y = 0.9472 X + 0.029 \text{ (mmol/L)}$$

$$Y = 0.9472 X + 0.116 \text{ (mg/dL)}$$

with a correlation coefficient $r^2 = 0.986$.

Interferences

Haemoglobin: No significant influence is observed up to 434 µmol/L (750 mg/dL).

Triglycerides: Do not use lipemic samples.

Total Bilirubin: No significant influence is observed up to 1000 µmol/L (58.5 mg/dL).

Direct Bilirubin: No significant influence is observed up to 320 µmol/L (18.72 mg/dL).

Other limitations are given by Young as a list of drugs and preanalytical variables known to affect this methodology (16, 17).

Calibration Stability

The reagent is calibrated on Day 0. The calibration stability is checked by testing 2 control specimens. The calibration stability is 20 days.

Note: A recalibration is recommended when reagent lots change, and when quality control results fall outside the range established.

Conversion Factor

$$\text{mmol/L} \times 4.01 = \text{mg/dL}$$

Urine

The performance data listed below have been obtained on the Yumizen C230/C240 analyzer.

Number of tests: approximately 2 x 140 tests

On Board Reagent Stability

Once opened, the reagent cassette placed in the refrigerated Yumizen C230/C240 compartment is stable for 70 days.

Sample volume: 2 µL/test

Lowest Detectable Level

The lowest detectable level represents the lowest measurable level of analyte that can be distinguished from zero. It is calculated as the absolute mean plus three standard deviations of 20 replicates of an analyte free sample. The lowest detectable level is estimated at 0.011 mmol/L (0.044 mg/dL).

Limit of Quantitation

The limit of quantitation is determined according to CLSI (NCCLS), EP17-A2 protocol (11) and equals 0.025 mmol/L (0.10 mg/dL).

Accuracy and Precision

Repeatability (within-run precision)

Repeatability according to the recommendations found in the CLSI (NCCLS), EP05-A3 protocol (12) with samples tested 20 times:

- 2 controls
- 3 specimens (low / medium / high levels)

Yumizen C Calcium

	Mean value mmol/L	Mean value mg/dL	CV %
Control specimen 1	1.89	7.59	1.02
Control specimen 2	2.85	11.41	1.50
Specimen 1	1.00	4.00	1.41
Specimen 2	2.83	11.36	0.94
Specimen 3	4.46	17.88	1.01

Reproducibility (total precision)

Reproducibility according to the recommendations found in the CLSI (NCCLS), EP05-A3 protocol (12) with samples tested in duplicate for 20 days (2 series per day):

- 2 controls
- 2 specimens (medium / high levels)

	Mean value mmol/L	Mean value mg/dL	CV %
Control specimen 1	1.91	7.66	2.4
Control specimen 2	2.86	11.47	1.3
Specimen 1	2.58	10.35	1.6
Specimen 2	3.88	15.56	1.4

Measuring Range

The assay confirmed a measuring range from 0.025 mmol/L (0.10 mg/dL) to 4.50 mmol/L (18.05 mg/dL). The measuring range is extended up to 18.00 mmol/L (72.18 mg/dL) with the automatic post-dilution.

The reagent linearity has been assessed up to 4.50 mmol/L (18.05 mg/dL) according to the recommendations found in the CLSI (NCCLS), EP06-Ed2 protocol (13).

Correlation

Patient samples: urine

Number of patient samples: 101

Specimens are correlated with a commercial reagent taken as reference according to the recommendations found in the CLSI (NCCLS), EP09c protocol (14).

Values ranged from 0.07 mmol/L (0.28 mg/dL) to 4.26 mmol/L (17.08 mg/dL).

The equation for the allometric line obtained using Passing-Bablok regression procedure (15) is:

$$Y = 0.9505 X + 0.0195 \text{ (mmol/L)}$$

$$Y = 0.9505 X + 0.078 \text{ (mg/dL)}$$

with a correlation coefficient $r^2 = 0.998$.

Interferences

Haemoglobin: No significant influence is observed up to 579 $\mu\text{mol/L}$ (1000 mg/dL).

Direct Bilirubin: No significant influence is observed up to 640 $\mu\text{mol/L}$ (37.44 mg/dL).

Other limitations are given by Young as a list of drugs and preanalytical variables known to affect this methodology (16, 17).

Calibration Stability

The reagent is calibrated on Day 0. The calibration stability is checked by testing 2 control specimens.

The calibration stability is 20 days.

Note: A recalibration is recommended when reagent lots change, and when quality control results fall outside the range established.

Conversion Factor:

$$\text{mmol/L} \times 4.01 = \text{mg/dL}$$

Reference

1. Thomas L. Clinical Laboratory Diagnostics. 1st ed. Frankfurt: TH-Books Verlagsgesellschaft (1998): 192-202.
2. Endres DB, Rude RK. Mineral and bone metabolism. In: Burtis C.A., Ashwood E.R., editors. Tietz Textbook of Clinical Chemistry. 3rd ed. Philadelphia: WB Saunders Company (1999): 1395-1457.
3. Matkovic V, Llich JZ, Andon MB, Hsieh LC, Tzagournis MA, Laggar BJ, Goel PK, Am. J. Clin. Nutr. (1995) **62** (2): 417-25.
4. Connerty HV, Briggs AR. Clin. Chem. (1965) **11**: 716-28.
5. Connerty HV, Briggs AR. Am. J. Clin. Path. (1966) **45**: 290-6.
6. Gindler EM, Kin JD, Am. J. Clin. Path. (1972) **58**: 376-82.
7. Bauer PJ. Anal. Biochem. (1981) **110**: 61-72.
8. Ehret W, Heil W, Schmitt Y, Töpfer G, Wisser H, Zawta B et al. Use of anticoagulants in diagnostic laboratory investigations and stability of blood, plasma and serum samples. WHO publication WHO/DIL/LAB/99.1 Rev.2: 25 (2002).
9. NCCLS. Urinalysis and collection, transportation and preservation of urine specimen; Approved guideline - 2nd Edition, NCCLS document GP16-A2, **21** (19).
10. Thomas L. Clinical Laboratory Diagnostics. 1st ed. Frankfurt: THBooks Verlagsgesellschaft (1998): 231-241.

Yumizen C Calcium

11. Evaluation of detection capability for clinical laboratory measurement procedures. Approved Guideline, 2nd ed., CLSI (NCCLS) document EP17-A2 (2012) **32** (8).
12. Evaluation of Precision of Quantitative Measurement Procedures. Approved Guideline, CLSI (NCCLS) document EP05-A3 (2014) **24** (25).
13. Evaluation of Linearity of Quantitative Measurement Procedures. 2nd Edition, CLSI (NCCLS) guideline EP06-Ed2 (2020) **40** (16).
14. Measurement Procedure Comparison and Bias Estimation Using Patient Samples. Approved Guideline, 3rd ed., CLSI (NCCLS) document EP09c (2018) **38** (12).
15. Passing H, Bablok W. A new biometrical procedure for testing the equality of measurements from two different analytical methods. J. Clin. Chem. Clin. Biochem. (1983) **21**: 709-720.
16. Young DS. Effects of Drugs on Clinical Laboratory Tests. 5th Edition, Washington, DC, AACC Press (2000).
17. Young DS. Effects of Preanalytical Variables on Clinical Laboratory Tests. 2nd Edition, Washington, DC, AACC Press (1997) **3**: 120-132.