


Yumizen G DDi 2 for YUMIZEN G1500/G1550



D 4910301363910322D



D 5400246530018952D

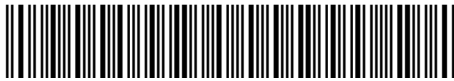

D 6001225100055905D



D 7002650250135178D

Reference No.:		1300036391	
Lot No.:		910301 (rev.1)	
Expiry date:		2022-03-31	
Coagulometers		MASTER CURVE	
Optical	Yumizen G1500/G1550	µg(FEU)/mL	ΔOD
		4,00	0,2465
		3,00	0,1895
		2,00	0,1225
		1,00	0,0559
		0,50	0,0265
		0,25	0,0135
Slope value		0,06	

Yumizen G DDi 2 for YUMIZEN G800



D 4910301363910322D



D 5400230030018002D



D 6001150100048005D



D 7002500250130131D

Reference No.:		1300036391	
Lot No.:		910301 (rev.1)	
Expiry date:		2022-03-31	
Coagulometers		MASTER CURVE	
Optical	Yumizen G800	µg(FEU)/mL	ΔOD
		4,00	0,2300
		3,00	0,1800
		2,00	0,1150
		1,00	0,0480
		0,50	0,0250
		0,25	0,0130
Slope value		0,06	



Yumizen G DDi 2 for YUMIZEN G400 DDi/G200

Reference No.:		1300036391	
Lot No.:		910301 (rev.1)	
Expiry date:		2022-03-31	
Coagulometers		MASTER CURVE	
		µg(FEU)/mL	ΔOD
Optical	Yumizen G400 DDi/G200	4,00	0,1015
		3,00	0,0757
		2,00	0,0500
		1,00	0,0247
		0,50	0,0105
		0,25	0,0068
		Slope value	0,068
Important notice			
<p>For Yumizen G400 DDi and Yumizen G200, make sure that the end time is set to 180 seconds in the D-Dimer setup and the curve type Lin-Lin p-p is selected in Calibration > D-Dimer > ug FEU/ml. For more information, please contact your local HORIBA Medical representative.</p>			